

## Laboratory Filtration Products



turning science **into solutions**

## ■ Sartorius Laboratory Filtration Products

Filtration and ultrafiltration are unavoidable process steps in nearly all environmental, chemistry and bio scientific laboratory applications.

Sartorius supplies a wide range of individual filter papers, microporous membranes, filtration devices, ultrafiltration devices and protein purification devices to suit these applications. This catalogue provides a condense overview of the Sartorius Lab Filtration product range, please contact us directly for specialty catalogues, available for in depth technical information.

## ■ Table of Contents

### □ Ultrafiltration

- 6 Protein Concentration
- 34 DNA Concentration
- 40 Protein Purification
- 46 Virus Purification and Concentration

### □ Filtration Devices

- 54 Minisart® Syringe Filters
- 65 Sartolab® P20 and Sartolab® P20 Plus
- 68 Sartolab® 150v
- 69 Sartolab® RF | BT
- 72 Chemical Compatibility

### □ Basic Filtration

- 76 Filter Papers
- 80 Glass & Quartz Microfiber Filters
- 84 Extraction Thimbles
- 85 Membrane Filters
- 91 Blotting | Chromatography Papers & Blotting Membranes
- 92 Filtration Equipment
- 102 Combisart® Manifolds
- 105 CONFIDENCE® – Validation Services







## ■ Ultrafiltration

### Table of Contents

- 6 Protein Concentration
- 34 DNA Concentration
- 40 Protein Purification
- 46 Virus Purification and Concentration



 Protein Concentration

- |                       |  |
|-----------------------|--|
| 7 Vivaspin® 500       | 18 Vivaclear Centrifugal Filters             |
| 8 Vivaspin® 2         | 19 Vivacell 70                               |
| 10 Centrisart I       | 22 Vivacell 100                              |
| 11 Vivaspin® 4        | 24 Vivacell 250                              |
| 12 Vivaspin® 6        | 26 Vivaflow 50                               |
| 13 Vivaspin® 15R      | 28 Vivaflow 200                              |
| 14 Vivaspin® Turbo 15 | 30 Vivapore Solvent Absorption Concentrators |
| 16 Vivaspin® 20       | 32 Ultrafiltration Membrane Filters          |



# ■ Vivaspin® 500



## 100 µl to 500 µl Samples

Vivaspin® 500 µl centrifugal filter units offer a simple, one step procedure for sample preparation. They can effectively be used in a fixed angle rotors accepting 2.2 ml centrifuge tubes.

The patented vertical membrane design and thin channel filtration chamber (US 5,647,990), minimises membrane fouling and provides high speed concentrations, even with particle laden solutions.



## □ Technical Specifications

### Vivaspin® 500

Concentrator capacity	Swing bucket rotor	do not use
	Fixed angle rotor	500 µl
Dimensions	Total length	50 mm
	Width	11 mm
	Active membrane area	0.5 cm <sup>2</sup>
	Hold-up volume, membrane and support	< 5 µl
	Dead stop volume	5 µl
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polypropylene
	Concentrator cap	Polycarbonate
	Membrane	Polyethersulfone


## □ Ordering Information

Vivaspin® 500 Polyethersulfone	Qty/Pk	Prod. No.
3,000 MWCO	25	VS0191
3,000 MWCO	100	VS0192
5,000 MWCO	25	VS0111
5,000 MWCO	100	VS0112
10,000 MWCO	25	VS0101
10,000 MWCO	100	VS0102
30,000 MWCO	25	VS0121
30,000 MWCO	100	VS0122
50,000 MWCO	25	VS0131
50,000 MWCO	100	VS0132
100,000 MWCO	25	VS0141
100,000 MWCO	100	VS0142
300,000 MWCO	25	VS0151
300,000 MWCO	100	VS0152
1,000,000 MWCO	25	VS0161
1,000,000 MWCO	100	VS0162
0.2 µm	25	VS0171
0.2 µm	100	VS0172
Starter pack (5 of each 5 k, 10 k, 30 k, 50 k, 100 k)	25	VS01S1

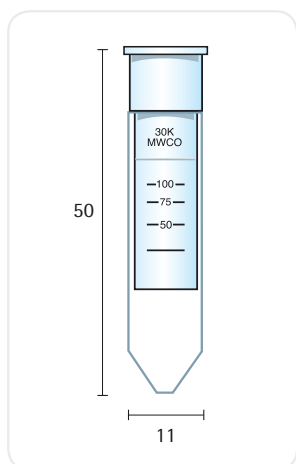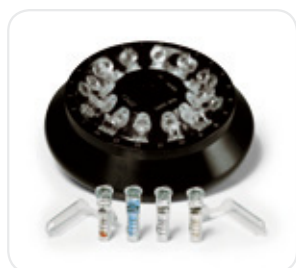

# Vivaspin® 2

Choice of Membranes

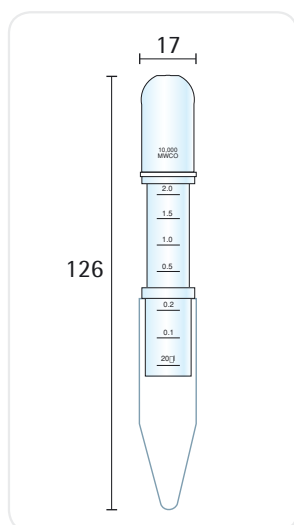

## 0.4–2 ml Samples

The Vivaspin® 2 bridges the gap between the 500 µl and 4 ml centrifugal concentrators. This device combines the speed of the classic Vivaspin® products with low internal surface and membrane area for superior recoveries from very dilute solutions.

Available with a choice of PES, Cellulose Triacetate and Hydrosart® membranes, Vivaspin® 2 offers the highest flexibility for process optimisation.

Also unique to the Vivaspin® 2, is the choice of directly pipetting the concentrate from the dead stop pocket built into the bottom of the concentrator, or alternatively reverse spinning into the concentrate recovery cap which can then be sealed for storage. Both methods result in near total concentrate recoveries.




## Technical Specifications

### Vivaspin® 2

Concentrator capacity	Swing bucket rotor	3 ml
	Fixed angle rotor	2 ml
Dimensions	Total length	126 mm
	Width	17 mm
	Active membrane area	1.2 cm <sup>2</sup>
	Hold-up volume of membrane	< 10 µl
	Dead stop volume	8 µl
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polycarbonate
	Concentrator cap	Polycarbonate
	Membrane	PES, CTA, HY

## Ordering Information

Vivaspin® 2 Polyethersulfone	Qty/Pk	Prod. No.
3,000 MWCO	25	VS0291
3,000 MWCO	100	VS0292
5,000 MWCO	25	VS0211
5,000 MWCO	100	VS0212
10,000 MWCO	25	VS0201
10,000 MWCO	100	VS0202
30,000 MWCO	25	VS0221
30,000 MWCO	100	VS0222
50,000 MWCO	25	VS0231
50,000 MWCO	100	VS0232
100,000 MWCO	25	VS0241
100,000 MWCO	100	VS0242
300,000 MWCO	25	VS0251
300,000 MWCO	100	VS0252
1,000,000 MWCO	25	VS0261
1,000,000 MWCO	100	VS0262
0.2 µm	25	VS0271
0.2 µm	100	VS0272
Starter pack (5 of each 5 k, 10 k, 30 k, 50 k, 100 k)	25	VS02S1

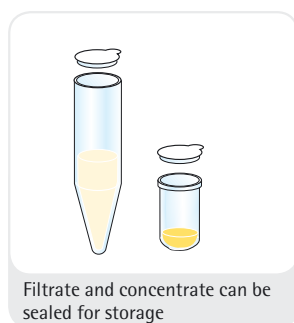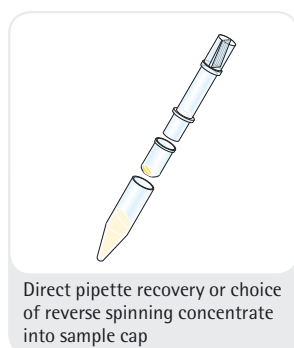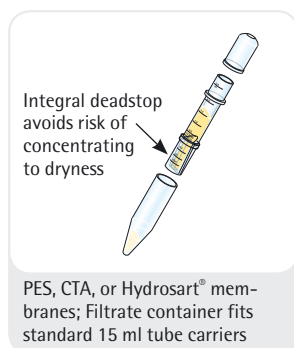


Vivaspin® 2 Cellulose Triacetate	Qty/Pk	Prod. No.
5,000 MWCO	25	VS02U1
5,000 MWCO	100	VS02U2
10,000 MWCO	25	VS02V1
10,000 MWCO	100	VS02V2
20,000 MWCO	25	VS02X1
20,000 MWCO	100	VS02X2

Vivaspin® 2 Hydrosart	Qty/Pk	Prod. No.
2,000 MWCO	25	VS02H91
2,000 MWCO	100	VS02H92
5,000 MWCO	25	VS02H11
5,000 MWCO	100	VS02H12
10,000 MWCO	25	VS02H01
10,000 MWCO	100	VS02H02
30,000 MWCO	25	VS02H21
30,000 MWCO	100	VS02H22

### Ordering Tips

- Choose a membrane pore size at least 50% smaller than the size of the molecule to be retained.
- Usually choose Polyethersulfone membranes for fastest concentrations.
- Usually choose Cellulose Triacetate for Protein Removal | Ultrafiltrate recovery.
- Usually choose Hydrosart® membranes for highest recovery with Ig fractions.

## Centrisart I

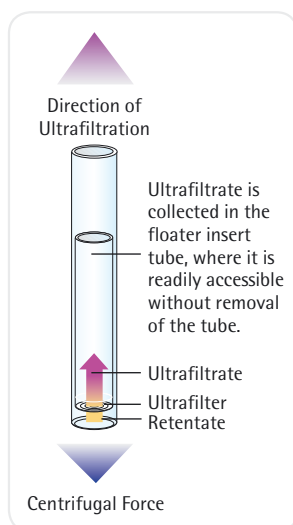

### 0.5–2.5 ml Samples

Centrisart I is a ready to use unit for small volume centrifugal ultrafiltration to separate proteins from low molecular weight substances in biological samples.

Centrisart I features a unique design, ultrafiltration in the opposite direction to the centrifugal force. This is so effective in preventing premature blockage of the filter that even whole blood samples can be deproteinized.

The ultrafiltrate is collected in the float insert tube, where it is readily accessible without removing the tube.

### Typical Applications Include:

- drug binding studies
- determination of metabolites in serum
- protein removal from blood samples
- cleaning of liposomes
- virus removal

## Technical Specifications

### Centrisart I

Concentrator capacity	Swing bucket rotor	2.5 ml
	Fixed angle rotor	2.5 ml
Dimensions	Total length	93 mm
	Width	14 mm
	Active membrane area	0.79 cm <sup>2</sup>
	Hold-up volume of membrane	< 5 µl
	Dead stop volume	100 µl
Materials of construction	Centrifuge tube	Polystyrene
	Float tube	Cellulose propionate
	Cap	Polyethylene
	Membrane	CTA, PES

## Ordering Information

	Qty/Pk	Prod. No.
5,000 MWCO CTA	12	13229-E
10,000 MWCO CTA	12	13239-E
20,000 MWCO CTA	12	13249-E
100,000 MWCO PES	12	13269-E
300,000 MWCO PES	12	13279-E
Starter pack (3 units each of 5 k, 10 k, 20 k, 100 k)	12	13209-E

### References

P. Nebinger and M. Koel  
Determination of acyclovir by ultrafiltration and high-performance liquid chromatography.  
J. Chromatography 619, 342–344 (1993)

F. da Fonseca-Wollheim, K.-G. Heinze, K. Lomsky and H. Schreiner  
Serum ultrafiltration for the elimination of endogenous interfering substances in creatinine determination.  
J. Clin. Chem. Clin. Biochem. 26, 523–525 (1988)

R. H. Christenson, S. D. Studenberg, S. Beck-Davis and F. A. Sedor  
Digoxin-like immunoreactivity eliminated from serum by centrifugal ultrafiltration before fluorescence polarization immunoassay of digoxin.  
Clinical Chemistry 33, 606–608 (1987)

## ■ Vivaspin® 4




### 1–4 ml Samples

Vivaspin® 4 ml concentrators are disposable ultrafiltration devices for the concentration of biological samples. Maximum initial sample volumes range from 1 ml to 4 ml. They can be effectively used in either swing bucket or fixed angle rotors accepting 15 ml centrifuge tubes.

The patented vertical membrane design and thin channel filtration chamber (US 5,647,990) minimises membrane fouling and provides high speed concentrations, even with particle laden solutions.

Vivaspin® 4 is available with the high flux polyethersulfone membrane range which is recommended for most solutions.

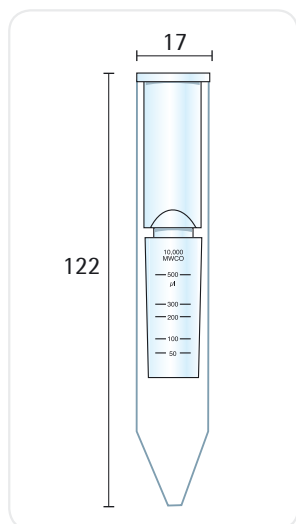

### □ Technical Specifications

#### Vivaspin® 4

Concentrator capacity	Swing bucket rotor	4 ml
	Fixed angle rotor	4 ml
Dimensions	Total length	122 mm
	Width	17 mm
	Active membrane area	2.0 cm <sup>2</sup>
	Hold-up volume of membrane	< 10 µl
	Dead stop volume	20 µl
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polypropylene
	Concentrator cap	Polycarbonate
	Membrane	Polyethersulfone

### □ Ordering Information

Vivaspin® 4 Polyethersulfone	Qty/Pk	Prod. No.
5,000 MWCO	25	VS0413
5,000 MWCO	100	VS0414
10,000 MWCO	25	VS0403
10,000 MWCO	100	VS0404
30,000 MWCO	25	VS0423
30,000 MWCO	100	VS0424
50,000 MWCO	25	VS0433
50,000 MWCO	100	VS0434
100,000 MWCO	25	VS0443
100,000 MWCO	100	VS0444
0.2 µm	25	VS0473
0.2 µm	100	VS0474
Starter pack (5 of each 5 k, 10 k, 30 k, 50 k, 100 k)	25	VS04S3

## ■ Vivaspin® 6




### 2–6 ml Samples

Vivaspin® 6 ml concentrators have been developed to offer increased volume flexibility and performance.

Vivaspin® 6 can process an impressive 6 ml in either swing bucket or fixed angle rotors accepting standard 15 ml conical bottom test tubes.

The Vivaspin® 6 features twin vertical membranes for unparalleled filtration speeds and 100× plus concentrations. Remaining volume is easy to read off the printed scale on the side of the concentrator and the modified dead stop pocket further simplifies direct pipette recovery of the final concentrate.

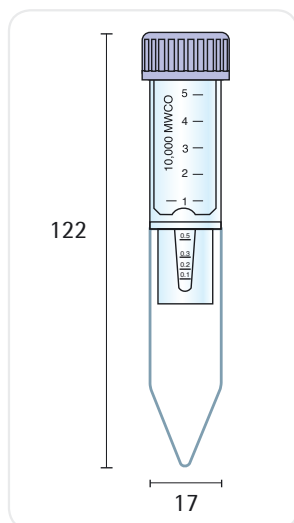

### □ Technical Specifications

#### Vivaspin® 6

Concentrator capacity	Swing bucket rotor	6 ml
	Fixed angle rotor	6 ml
Dimensions	Total length	122 mm
	Width	17 mm
	Active membrane area	2.5 cm <sup>2</sup>
	Hold-up volume of membrane	< 10 µl
	Dead stop volume	30 µl
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polycarbonate
	Concentrator cap	Polypropylene
	Membrane	Polyethersulfone

### □ Ordering Information

Vivaspin® 6 Polyethersulfone	Qty/Pk	Prod. No.
3,000 MWCO	25	VS0691
3,000 MWCO	100	VS0692
5,000 MWCO	25	VS0611
5,000 MWCO	100	VS0612
10,000 MWCO	25	VS0601
10,000 MWCO	100	VS0602
30,000 MWCO	25	VS0621
30,000 MWCO	100	VS0622
50,000 MWCO	25	VS0631
50,000 MWCO	100	VS0632
100,000 MWCO	25	VS0641
100,000 MWCO	100	VS0642
300,000 MWCO	25	VS0651
300,000 MWCO	100	VS0652
1,000,000 MWCO	25	VS0661
1,000,000 MWCO	100	VS0662
0.2 µm	25	VS0671
0.2 µm	100	VS0672
Starter pack (5 of each 5 k, 10 k, 30 k, 50 k, 100 k)	25	VS06S1


## ■ Vivaspin® 15R



### 2–15 ml Samples

Vivaspin® 15R is targeting the volume segment 2 to 15 ml with a modified regenerated cellulose membrane; Hydrosart®. This membrane is ideal where extremely high recovery with very low adsorption is needed, for example in applications such as desalting and concentration of Ig fractions.

- Ultimate recovery at low adsorption (95–98%)
- Extremely short concentration time (30 × in 15 min.)
- Convenient application protocol with easy handling
- Easy scale-up to Vivaflow 200 with Hydrosart® membrane for volumes up to 5 litres
- Very small hold up volume (< 20 µl)

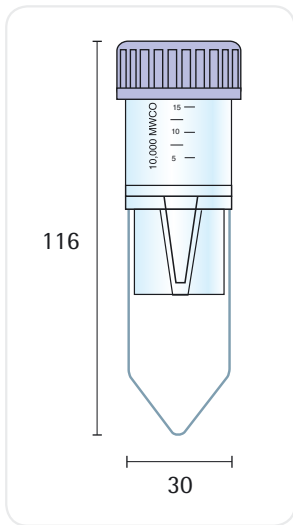


### □ Technical Specifications

#### Vivaspin® 15R


Concentrator capacity	Swing bucket rotor	15 ml
	Fixed angle rotor	12.5 ml
Dimensions	Total length	116 mm
	Width	30 mm
	Active membrane area	3.9 cm <sup>2</sup>
	Hold-up volume of membrane	< 20 µl
	Dead stop volume	30 µl
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polypropylene
	Concentrator cap	Polycarbonate
	Membrane	Hydrosart

### □ Ordering Information

Vivaspin® 15R Hydrosart	Qty/Pk	Prod. No.
2,000 MWCO	12	VS15RH91
2,000 MWCO	48	VS15RH92
5,000 MWCO	12	VS15RH11
5,000 MWCO	48	VS15RH12
10,000 MWCO	12	VS15RH01
10,000 MWCO	48	VS15RH02
30,000 MWCO	12	VS15RH21
30,000 MWCO	48	VS15RH22

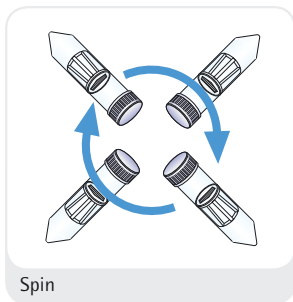

# Vivaspin® Turbo 15

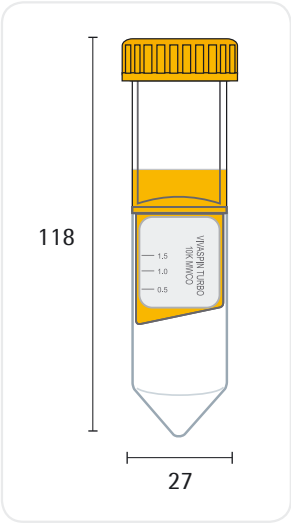

### 4 – 15 ml

Vivaspin® Turbo 15 is the newest member of the ultrafiltration family and allows fastest sample concentration with highest recoveries. This device can handle up to 15 ml sample volume in swing bucket rotors and 11 ml in fixed angle rotors accepting 50 ml centrifuge tubes.

The UV joining technology allows for a smooth joint transition between membrane and plastic housing, allowing the collection of the complete concentrated sample into the unique pipette friendly dead stop pocket.

The Vivaspin® Turbo 15 optimized design and sleek internal profile ensure maximum process speeds right the way down to the last few micro litres leading to > 100 fold concentration.

## Technical Specifications

### Vivaspin® Turbo 15

Materials	Body	Styrene butadiene copolymere
	Filtrate vessel	Polypropylene
	Concentrator cap	Polypropylene
	Membrane	Polyethersulfone (PES)
Dimensions	Total length (concentrator insert)	77 mm
	Total length (in tube with cap)	118 mm
	Diameter (concentrator insert)	27 mm
	Active membrane area	7.2 cm <sup>2</sup>
	Hold up volume of membrane	<10 µl
	Dead stop volume in swing out	100 µl
	Dead stop volume in fixed angle	60 µl
Concentrator capacity	Swing bucket rotor	15 ml
	Fixed angle rotor (25°)	11 ml
Maximum speed	4000 × g	4000 × g
Sterilization	ETO or 70% EtOH	
Removal of endotoxins [Depyrogenization]	Flushing with 1N NaOH	

 Ordering Information

<b>Vivaspin® Turbo 15 Polyethersulfone</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
3,000	12	VS15T91
3,000	48	VS15T92
5,000	12	VS15T11
5,000	48	VS15T12
10,000	12	VS15T01
10,000	48	VS15T02
30,000	12	VS15T21
30,000	48	VS15T22
50,000	12	VS15T31
50,000	48	VS15T32

## ■ Vivaspin® 20



### 5–20 ml Samples

Vivaspin® 20 ml centrifugal concentrators have been developed to offer increased volume flexibility and performance.

Vivaspin® 20 handles up to 20 ml in swing bucket centrifuges and 14 ml in 25° fixed angle rotors accepting 50 ml centrifuge tubes.



Featuring twin vertical membranes for unparalleled filtration speeds the Vivaspin® 20 can achieve 100x plus concentrations.


Remaining volume is easy to read off the printed scale on the side of the concentrator and the modified dead stop pocket further simplifies direct pipette recovery of the final concentrate.

### More Process Flexibility

Vivaspin® 20 is available with unique accessories and operating methods that are designed to provide more process flexibility and further time saving.

### Gas Pressure Filtration

When an appropriate centrifuge is unavailable, or for single sample processing, Vivaspin® 20 can be filled with up to 15 ml and then pressurised for bench top concentration. For even faster processing, gas pressure can be combined with centrifugal force. "Pressure-fugation" is particularly suitable for difficult or viscous samples such as serum, or when using a low process temperature which reduces filtration speed, and generally when minimum process time is essential.

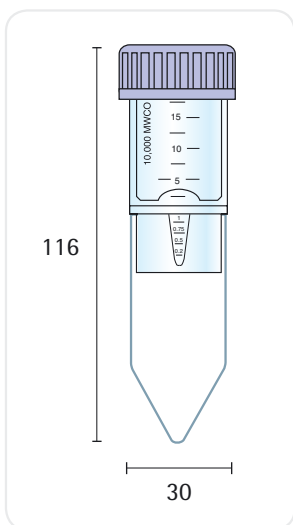

## □ Technical Specifications

### Vivaspin® 20


Concentrator capacity	Swing bucket rotor	20 ml
	Fixed angle rotor	14 ml
	With pressure head	15 ml
Dimensions	Total length	116 mm
		125 mm with pressure head
	Width	30 mm
	Active membrane area	6.0 cm <sup>2</sup>
	Hold-up volume of membrane	< 20 µl
Materials of construction	Dead stop volume	50 µl
	Body	Polycarbonate
	Filtrate vessel	Polycarbonate
	Concentrator cap	Polypropylene
	Pressure head	Acetal   aluminium
	Membrane	Polyethersulfone

## □ Ordering Information

Vivaspin® 20 Polyethersulfone	Qty/Pk	Prod. No.
3,000 MWCO	12	VS2091
3,000 MWCO	48	VS2092
5,000 MWCO	12	VS2011
5,000 MWCO	48	VS2012
10,000 MWCO	12	VS2001
10,000 MWCO	48	VS2002

<b>Vivaspin® 20 Polyethersulfone</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
30,000 MWCO	12	VS2021
30,000 MWCO	48	VS2022
50,000 MWCO	12	VS2031
50,000 MWCO	48	VS2032
100,000 MWCO	12	VS2041
100,000 MWCO	48	VS2042
300,000 MWCO	12	VS2051
300,000 MWCO	48	VS2052
1,000,000 MWCO	12	VS2061
1,000,000 MWCO	48	VS2062
0.2 µm	12	VS2071
0.2 µm	48	VS2072
Starter pack (2 of each 5 k, 10 k, 30 k, 50 k, 100 k, 0.2 µm)	12	VS20S1

<b>Vivaspin® 20 Accessories</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
Air pressure controller (APC)	1	VCA002
Charge valve for pressure head	1	VCA005
Diafiltration cups	12	VSA005
Female connector	1	VCA010
Male connector	1	VCA011
4 mm OD pneumatic tube (3 m)	1	VCA012
Vivaspin® 20 pressure head	1	VCA200

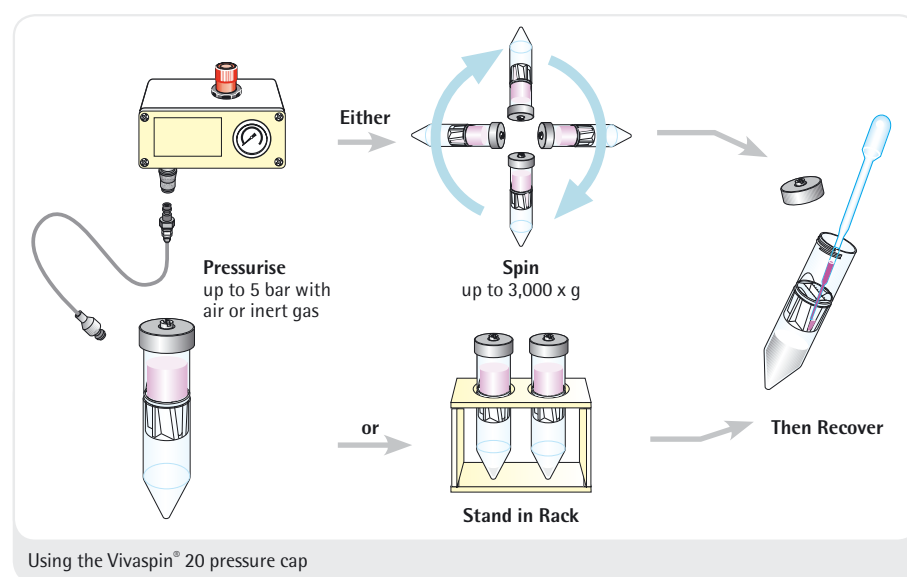

## ■ Vivaclear Centrifugal Filters



Vivaclear centrifugal filters are disposable microfiltration devices for the fast and reliable clarification | filtration of biological samples in the range 100 µl to 500 µl. They can be used in fixed angle rotors accepting 2.2 ml centrifuge tubes.

### Product Features

- High-flux Polyethersulphone membrane
- 0.8 µm pore size

- Low hold up volume (<5 µl)
- Fast and reproducible performance

### Applications

- Clarification of samples before loading onto Vivapure protein purification spin columns
- Removal of particles and particulates
- Filtration of plasma and serum
- Filtration of cells or cell debris

## □ Technical Specifications

### Vivaclear Centrifugal Filters

Rotor	40–45° Fixed angle rotor 500 µl	
Pore size	0.8 µm	
Dimensions	Total length	43 mm
	Filtrate collection tube diameter	11 mm
	Active membrane area	0.34 cm <sup>2</sup>
	Hold-up volume, membrane plus support	< 5 µl
	Maximum RCF	2,000 × g
Materials of construction	Body	Polypropylene
	Membrane	Polyethersulphone
	Filtrate collection tube	Polypropylene

## □ Ordering Information

	Qty/Pk	Prod. No.
Vivaclear Mini 0.8 µm PES	100	VK01P042

## ■ Vivacell 70



### 10–70 ml Samples

Vivacell 70 combines the ease of use of centrifugal devices with the flexibility and control provided by pressurised ultrafiltration cells. Vivacell 70 is inexpensive, quick and easy to assemble, requires no tubing connections or stirring mechanisms and can be adapted to equipment availability or to specific user preferences.

For convenience, simply spin in a large capacity centrifuge (rotors accepting 250 ml bottles). For highest speeds particularly with difficult samples, pressurise the device with air or inert gas before centrifuging.

For more process control or for single samples, combine gas pressure with a gentle orbital shake, or you can even pressurise and then leave standing on a bench top or in a refrigerator for highest simplicity with minimum equipment requirements.

The longitudinal membrane inhibits fouling, whilst the built-in dead stop will hinder further concentration when residual volume drops below 150  $\mu$ l.

### □ Technical Specifications

#### Vivacell 70

Concentrator capacity	Swing bucket rotor	70 ml
	Fixed angle rotor	50 ml
	With pressure head	70 ml
	With pressure-fuge head	50 ml
Di mensions	Total length	119 mm standard centrifugal 185 mm with pressure head 125 mm with pressure fuge head
	Width	62 mm
	Active membrane area	20 cm <sup>2</sup>
	Hold up volume of membrane	< 200 $\mu$ l
	Dead stop volume	150 $\mu$ l
Operating requirements	Rotor type	Swing bucket or fixed angle
	Minimum rotor angle	25°
	Rotor cavity	To fit 250 ml (62 mm) centrifuge bottles
	Maximum speed	1,000 g
Materials of construction	Maximum pressure	5 bar (75 psi)
	Body	Polycarbonate
	Filtrate vessel	Polycarbonate
	Concentrator cap	Santoprene
	Pressure head   pressure fuge head	Acetal
	Membrane	Polyethersulfone

### Ordering Information


#### **Vivacell 70 Polyethersulfone – Concentrator Bodies With Polycarbonate Filtrate Bottles**

	Qty/Pk	Prod. No.
5,000 MWCO	2	VS6011
10,000 MWCO	2	VS6001
30,000 MWCO	2	VS6021
50,000 MWCO	2	VS6031
100,000 MWCO	2	VS6041
0.2 µm	2	VS6071

#### **Vivacell 70 Polyethersulfone – Concentrator Body Only**


5,000 MWCO	10	VS6012
10,000 MWCO	10	VS6002
30,000 MWCO	10	VS6022
50,000 MWCO	10	VS6032
100,000 MWCO	10	VS6042
0.2 µm	10	VS6072

<b>Vivacell 70 Accessories</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
Air pressure controller (APC) complete with pressure gauge, regulator, over-pressure safety valve, female connector to Sartorius Stedim Biotech pressure products and 1 m extension line (4 mm pneumatic tubing) with male and female connectors and 1 m of 6 mm inlet tubing	1	VCA002
250 ml centrifuge bottle – standard caps	4	VSA003
Modified caps for use in fixed angle rotors with 250 ml centrifuge bottles	2	VCA004
Charge valve for pressure-fuge head	1	VCA005
Replacement seals for pressure-fuge head (VCA701)	10	VCA007
Female connector	1	VCA010
Male connector	1	VCA011
4 mm pneumatic tubing (3 m)	1	VCA012
Vivacell 70 pressure head with reservoir and filtrate bottle (bench top use)	1	VCA700
Vivacell 70 pressure-fuge head (for use in centrifuge)	2	VCA701

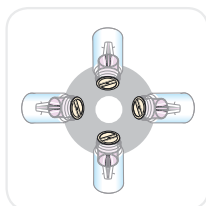


#### **Centrifuge**

- Process convenience
- Low shear, no foaming
- Less visual control

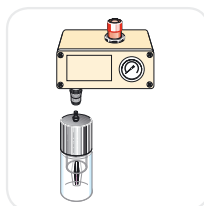


#### **Pressure**

- Simplicity and highest process control
- Ideal for refrigerated use
- Slower concentrations

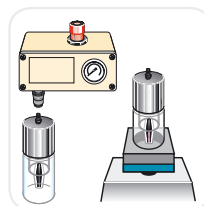

#### **Pressure-Shake**

- Speed and process control
- Ideal for single samples
- If left unattended can concentrate to dryness

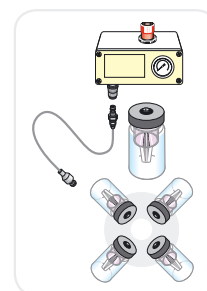

#### **Pressure-Fuge**

- Fastest processing
- Ideal with low MWCO or with difficult solutions
- Less visual control

# Vivacell 100





**20–100 ml Samples**  
Vivacell 100 is the latest member of the Vivacell family and bridges the volume range between the Vivacell 70 and the Vivacell 250.

The patented vertical membrane design allows highest performance and unmatched flexibility.



Vivacell 100 is a unique and innovative concentrator for volumes from 20 ml to 100 ml, which utilizes pressure, centrifuge or pressure-shake to rapidly concentrate even samples with very high particle loading.

Vivacell 100 is designed for centrifugal concentration of samples up to 100 ml which makes it the largest centrifugal unit available. At the same time, the new construction design allows for maximum centrifugal force of 4,000× g to be used for even faster concentration.

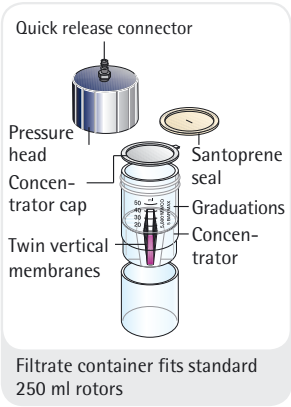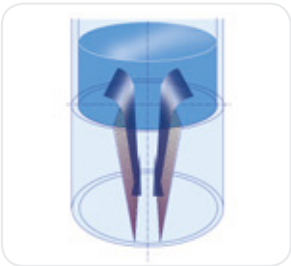

- Vivacell 100 Utilizes:**
- Pressure
  - Centrifuge
  - Pressure-shake

Like the smaller Vivacell 70 unit, Vivacell 100, when used as a centrifugal device, fits only into swing bucket rotors accepting 250 ml bottles.

Vivacell 100 units can also be used for single or extremely sensitive samples in the pressurized mode only and left on the bench or placed on a laboratory shaker for faster concentration. It can also be kept in a pressurized mode in the refrigerator. Handling is made easy by use of quick connectors. In whichever mode Vivacell 100 is used, the vertical membrane design inhibits membrane fouling while the built-in dead stop impedes concentration to dryness and loss of sample.


## Technical Specifications

### Vivacell 100

Concentrator capacity	Swing bucket rotor	90 ml
	With pressure head	98 ml
Dimensions	Total length	123 mm centrifugal 197 mm with pressure head
	Width	62 mm
	Active membrane area	23.5 cm <sup>2</sup>
	Hold-up volume of membrane	< 250 µl
	Dead stop volume	350 µl
Operating requirements	Rotor type	Swing bucket
	Rotor cavity	To fit 250 ml (62 mm) centrifuge bottles (maximum cavity depth 105 mm)
	Maximum speed	2,000 g
	Maximum pressure	5 bar (75 psi)
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polycarbonate
	Concentrator cap	Santoprene
	Pressure head	Acetal
	Membrane	Polyethersulfone


## Ordering Information

<b>Vivacell 100 Polyethersulfone with Polypropylene Concentrator Cap</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
5,000 MWCO	2	VC1011
5,000 MWCO	10	VC1012
10,000 MWCO	2	VC1001
10,000 MWCO	10	VC1002
30,000 MWCO	2	VC1021
30,000 MWCO	10	VC1022
50,000 MWCO	2	VC1031
50,000 MWCO	10	VC1032
100,000 MWCO	2	VC1041
100,000 MWCO	10	VC1042
300,000 MWCO	2	VC1051
300,000 MWCO	10	VC1052
1,000,000 MWCO	2	VC1061
1,000,000 MWCO	10	VC1062
0.2 µm	2	VC1071
0.2 µm	10	VC1072
<b>Vivacell 100 Accessories</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
Air pressure controller (APC) complete with pressure gauge, regulator, over-pressure safety valve, female connector, 1 m extension line (4 mm pressure tubing) with male and female connectors and 1 m of 6 mm inlet tubing	1	VCA002
Female connector	1	VCA010
Male connector	1	VCA011
4 mm pressure tubing (3 m)	1	VCA012
Santoprene replacement seals	10	VCA014
Vivacell 100 pressure head with replacement seals (5)	1	VCA800

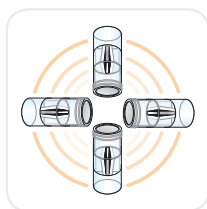


### Centrifuge

- Process convenience
- Low shear, no foaming
- Less visual control

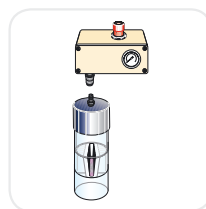

### Pressure

- Simplicity and highest process control
- Ideal for refrigerated use
- Slower concentrations

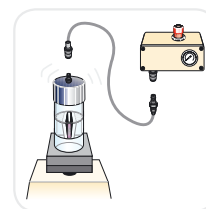

### Pressure-Shake

- Speed and process control
- Ideal for single samples

# Vivacell 250



### 50–250 ml Samples

The Vivacell 250 is a totally new concept for the concentration of larger biological samples. This product offers numerous advantages when compared to stirred cells.

- One size handles a volume range from under 50 ml to 250 ml.
- Use free standing on a bench top or in a refrigerator for maximum simplicity, or use on laboratory shaker for fastest concentrations.
- The unique conical dead stop built into the bottom of the membrane insert allows concentrations to under 1 ml.


- The gentle vortex action controls membrane polarisation whilst greatly reducing the shear effects typical of stirring mechanisms.
- Set up or membrane replacement takes just a few seconds. Quick connect fittings and simple screw closure further enhance ease of use.

Unique membrane module takes seconds to replace. Concentrate can be easily monitored through the graduated inspection window.

## Technical Specifications

### Vivacell 250

Concentrator capacity	250 ml	
Max pressure	4 bar (60 psi)	
Dimensions	Width	116 mm
	Height (incl. pressure indicator)	235 mm
	Active membrane area	40 cm <sup>2</sup>
	Hold-up vol. memb. & support	< 200 µl
	Dead stop volume	600 µl
Materials of construction	Screw closure	Acetal
	Pressure head	Acetal
	Quick release connector	Acetal
	Concentrator body   sleeve	Polycarbonate
	Filtrate container	Polycarbonate

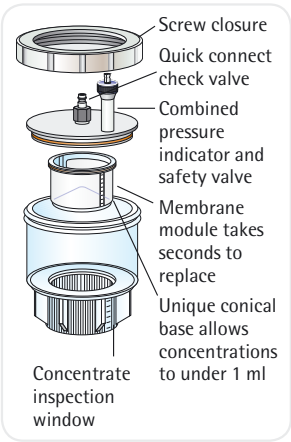

## □ Ordering Information

<b>Vivacell 250</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
Vivacell 250 complete with pressure head, pressure indicator   over-pressure release valve, quick release connection to APC, 2 sample reservoirs, filtrate container & insert tool	1	VCA250
<b>Vivacell 250 Polyethersulfone Inserts</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
5,000 MWCO	5	VC2511
10,000 MWCO	5	VC2501
30,000 MWCO	5	VC2521
50,000 MWCO	5	VC2531
100,000 MWCO	5	VC2541
0.2 µm	5	VC2571
Starter kit (1 of each 5 k, 10 k, 30 k, 50 k, 100 k)	5	VC25S1
<b>Vivacell 250 Accessories</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
Air pressure controller (APC) complete with pressure gauge, regulator, over-pressure safety valve, female connector to Sartorius Stedim Biotech pressure products and 1 m extension line (4 mm pneumatic tubing) with male and female connector and 1 m of 6 mm inlet tubing	1	VCA002
Replacement pressure indicator   over pressure relief valve	1	VCA008
Vivacell 250 maintenance kit (includes one sample reservoir and filtrate container, and "O" ring seals for pressure head)	1	VCA009
Female connector	1	VCA010
Male connector	1	VCA011
4 mm OD pressure tubing (3 m)	1	VCA012
Replacement pressure head & screw closure	1	VCA015

# Vivaflow 50



### 100 ml to 5 Litres

The unique and patent protected Vivaflow 50 system provides a standard of ease of use, performance, flexibility and economy which is unrivalled by any laboratory or pilot scale filtration system on the market.


### Unique Features

- Thin channel flip-flow recirculation path provides high cross flow velocities with minimum pump requirements.
- No need for pressure holders.
- Crystal clear for simple control of remaining hold up and membrane status.
- Unique Interlocking modules with series connectors for easy scale up.
- Disposable.

### Unique Performance

- A single 50 cm<sup>2</sup> module will typically reduce 500 ml to less than 15 ml in under 50 minutes.
- Less than 10 ml minimum system recirculation for highest concentrations.
- Less than 500 µl non recoverable hold up volume.
- Near total recoveries achievable with a single 10 ml rinse.


Unique “flip-flow” thin channel flow path results in high turbulence and linear velocity for exceptional flux even at high concentrations

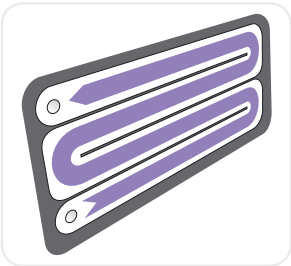


## Technical Specifications

### Vivaflow 50

Dimensions	Overall L   H   W	107   84   25 mm
	Channel W   H	15 mm   0.3 mm
	Active membrane area	50 cm <sup>2</sup>
	Hold up volume (module)	1.5 ml
	Minimum recirculation volume	< 10 ml
Operating conditions	Non recoverable hold-up	< 0.5 ml
	Pump flow	200-400 ml/min
	Maximum pressure	3 bar (45 psi)
Materials of construction	Maximum temperature	60°C
	Main housing	Polycarbonate
	Flow channel	TPX (PMP)
	Membrane support	TPX (PMP)
	Seals and O rings	Silicone
	Pressure indicator	Polypropylene, SS spring
	Flow restrictor	Polypropylene
Single module	Fittings	Nylon
	Tubing	PVC (medical grade)

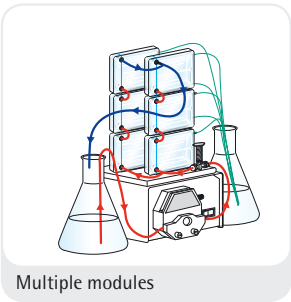

Multiple modules

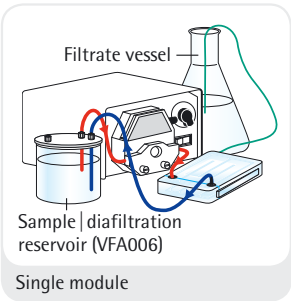

Single module

## □ Ordering Information

<b>Vivaflow 50*</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
3,000 MWCO PES	2	VF05P9
5,000 MWCO PES	2	VF05P1
10,000 MWCO PES	2	VF05P0
30,000 MWCO PES	2	VF05P2
50,000 MWCO PES	2	VF05P3
100,000 MWCO PES	2	VF05P4
1,000,000 MWCO PES	2	VF05P6
0.2 µm PES	2	VF05P7
100,000 MWCO RC	2	VF05C4

\* Vivaflow 50 modules include filtrate tube, size 16 peristaltic tubing, flow restrictor and fittings.

### **Vivaflow 50 Complete System Comprises**

Pump (240 V), Easy load pump head (size 16), tubing, 500 ml sample   diafiltration reservoir, module stand, pressure indicator, T connectors, series interconnectors	1	VFS502
Pump (115 V), Easy load pump head (size 16), tubing, 500 ml sample   diafiltration reservoir, module stand, pressure indicator, T connectors, series interconnectors	1	VFS504

### **PVC Tubing and Fittings**

Size 16 PVC pump tubing (3 metres, 3.2 × 1.6 mm)	VFA004
Flow restrictor set (2 × 0.4, 0.6, 0.8 mm)	VFA009
T connectors for running 2 stacks (2 pieces)	VFA030
Series interconnectors (6 pieces)	VFA031
Female luer fittings (10 pieces)	VFA032
VF50 tubing Kit (2 × 1 m size 16 PVC tubing with inlet fittings, 2 × 50 cm size 16 PVC tubing with 0.6 mm flow restrictors, 1 × series interconnector)	VFA034
Flow restrictor 0.6 mm (6 pieces)	VFA035

### **Accessories**

Masterflex economy drive variable speed peristaltic pump (240 V)	VFP001
Masterflex economy drive variable speed peristaltic pump (115 V)	VFP002
500 ml sample and   or diafiltration reservoir	VFA006
Masterflex easy load pump head – size 16	VFA012
Vivaflow 50 stand	VFA016
Pressure indicator (1-3 bar)	VFA020

# Vivaflow 200




### 0.5 to 5 Litres

Concentrate 250 ml to under 20 ml in just a few minutes or concentrate one litre 50 times in less than 30 minutes. Alternatively, use two Vivaflow 200's in parallel and concentrate 5 litres in under 75 minutes.

Near total sample recoveries can be expected with most solutions.

The economical standard package comes complete with tubing, pressure indicator, flow restrictor and high pressure pump tubing. All you need is a peristaltic pump capable of handling 6.4 mm OD (size 16) tubing. Should your pump head require larger tubing, link your own peristaltic tube up to the standard product, using the interconnector provided.

Two modules in parallel will concentrate 5 litres in under 75 minutes

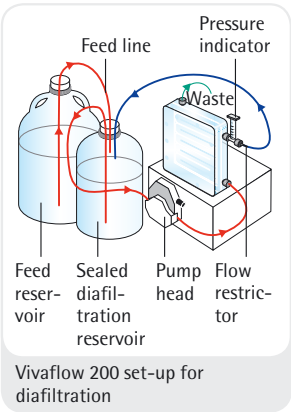

## Technical Specifications

### Vivaflow 200

Dimensions	Overall L   H   W	126   138   38 mm
	Channel W   H	10 mm   0.4 mm
	Active membrane area	200 cm <sup>2</sup>
	Hold up volume (module)	5.3 ml
	Min. recirculation volume	< 20 ml
	Non recoverable hold-up	< 2 ml
Materials of construction	Main housing	Acrylic
	Flow channel	Acrylic
	Membrane support	Polypropylene
	Seals and O rings	Silicone
	Pressure indicator	Polypropylene, SS spring
	Flow restrictor	Polypropylene
	Fittings	Nylon
	Tubing	PVC (medical grade)
Operating conditions	Pump flow	200–400 ml/min
	Maximum pressure	4 bar (60 psi)
	Maximum temperature	60°C

## □ Ordering Information

<b>Vivaflow 200*</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
2,000 MWCO Hydrosart®	1	VF20H9
3,000 MWCO PES	1	VF20P9
5,000 MWCO PES	1	VF20P1
10,000 MWCO PES	1	VF20P0
30,000 MWCO PES	1	VF20P2
50,000 MWCO PES	1	VF20P3
100,000 MWCO PES	1	VF20P4
0.2 µm PES	1	VF20P7
100,000 MWCO RC	1	VF20C4
2,000 MWCO Hydrosart®	1	VF20H9
5,000 MWCO Hydrosart®	1	VF20H1
10,000 MWCO Hydrosart®	1	VF20H0
30,000 MWCO Hydrosart®	1	VF20H2

\* Vivaflow 200 modules include pressure indicator, flow restrictor and size 16 pvc peristaltic tubing and fittings.

### **Vivaflow 200 Complete System Comprises**

Pump (240 V), Easy load pump head (size 16), tubing, 500 ml sample   diafiltration reservoir	1	VFS202
Pump (115 V), Easy load pump head (size 16), tubing, 500 ml sample   diafiltration reservoir	1	VFS204

### **Accessories**

Masterflex economy drive variable speed peristaltic pump (240 V)	VFP001
Masterflex economy drive variable speed peristaltic pump (115 V)	VFP002
500 ml sample and   or diafiltration reservoir	VFA006
Masterflex easy load pump head – size 16	VFA012
Masterflex easy load pump head – size 15	VFA013

### **Tubing and Fittings**

Size 15 pvc pump tubing and Luer fittings (3 m, 4.8 × 2.6 mm)	VFA003
Size 16 pvc pump tubing and Luer fittings (3 m, 3.2 × 1.6 mm)	VFA004
Y connector (size 15 to 2 × size 16)	VFA005
Flow restrictor set (2 × 0.4, 0.6, 0.8 mm)	VFA009
Female luer fittings size 16 (10 pieces)	VFA032
Flow restrictors 0.6 mm (6 pieces)	VFA035
Female luer fittings size 15 (10 pieces)	VFA036

# Vivapore Solvent Absorption Concentrators




### 3 ml–20 ml Samples

With no need for additional equipment, pressure or vacuum, solvent absorption is the most economic and user friendly concentration technique available to the clinician and research scientist.

Just fill the unit with the solution to be concentrated, wait for the desired concentration level to be achieved and then pipette the concentrated sample from the bottom of the reservoir.

Vivapore is ideal for general purpose laboratory concentration or purification prior to further analysis. It is particularly suited for labile solutions that can denature with alternative shear or pressure inducing methods or that require processing in a cold room environment.

Vivapore concentrators extend the solvent absorption technique to a totally new level of performance, application potential and ease of use.

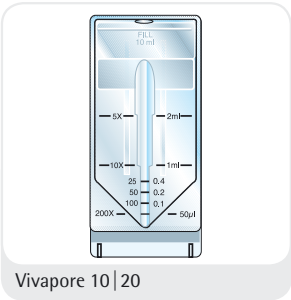

### Technical Specifications

	Vivapore 5	Vivapore 10   20
Membrane material	PES	PES
Membrane MWCO	7,500	7,500
Membrane surface area	20 cm <sup>2</sup>	28 cm <sup>2</sup>
Reservoir material	SAN	SAN
Volume range	1–5 ml	2–10 ml   20 ml*
Minimum concentrate volume	50 µl	50 µl
Vivapore overall dimensions		
Width (mm)	42	46
Height (mm)	82	100

## □ Ordering Information

<b>Vivapore 5*</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
7,500 MWCO PES	4	VP0503
7,500 MWCO PES	30	VP0501

\* includes stand and recovery pipettes

### **Requires Stand**

7,500 MWCO PES	100	VP0502
----------------	-----	--------

### **Vivapore 10 | 20\***

7,500 MWCO PES	4	VP2003
7,500 MWCO PES	30	VP2001

\* includes stand and recovery pipettes

### **Requires Stand**


7,500 MWCO PES	100	VP2002
----------------	-----	--------

### **Accessories**

Disposable stands for 4 units	6	VPA002
Plastic recovery pipettes (Vivapore 10   20)	100	VPA005
10 ml expansion reservoir (Vivapore 10   20)	10	VPA006
Plastic recovery pipettes (Vivapore 5)	100	VPA007

## Ultrafiltration Membrane Filters

PES 146, CTA 145 and Hydrosart 144

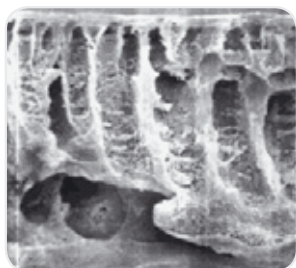

### Polyethersulfone (PES)

This is a general purpose membrane that provides excellent performance with most solutions when retentate recovery is of primary importance. Polyethersulfone membranes exhibit no hydrophobic or hydrophilic interactions and are usually preferred for their low fouling characteristics, exceptional flux and broad pH range.

microsolutes, these membranes are to be preferred for sample cleaning and protein removal and when high recovery of the filtrate solution is of primary importance.

### Cellulose Triacetate (CTA)

High hydrophilicity and very low non-specific binding characterize this membrane. Cast without any membrane support that could trap or bind passing

### Hydrosart

These membranes are also highly hydrophilic and are often preferred for their higher protein recovery when processing very dilute solutions. Resistance to autoclaving, ease of cleaning and extended chemical resistance also characterize this type of membrane.

## Technical Specifications

### Typical Performance for Polyethersulfone, Type 146

Thickness	120 µm	
pH range	1–14	
Waterflux	MWCO 10,000	0.2 ml/min/cm <sup>2</sup>
Protein retention	Cytochrome C	95%

### Specifications for Cellulose Triacetate, Type 145

Thickness	120 µm	
pH range	4–8	
Waterflux	MWCO 10,000	0.11 ml/min/cm <sup>2</sup>
Protein retention	Cytochrome C	90%

### Specifications for Hydrosart, Type 144

Thickness	180 µm	
pH range	1–13	
Waterflux	MWCO 10,000	0.08 ml/min/cm <sup>2</sup>
Protein retention	Cytochrome C	99%

## Ordering Information

### Polyethersulfone Membrane Filters, Type 146

Diameter in mm	MWCO	Qty/Pk	Prod. No.
47	1,000 Dalton	10	14609--47-----D
63	1,000 Dalton	10	14609--63-----D
76	1,000 Dalton	10	14609--76-----D
25	5,000 Dalton	10	14629--25-----D
47	5,000 Dalton	10	14429--47-----D
63	5,000 Dalton	10	14629--63-----D
76	5,000 Dalton	10	14629--76-----D
25	10,000 Dalton	10	14639--25-----D

### □ Ordering Information

#### Polyethersulfone Membrane Filters, Type 146

Diameter in mm	MWCO	Qty/Pk	Prod. No.
63	10,000 Dalton	10	14639--63-----D
76	10,000 Dalton	10	14639--76-----D
150	10,000 Dalton	10	14639-150-----D
25	30,000 Dalton	10	14659--25-----D
63	30,000 Dalton	10	14659--63-----D
76	30,000 Dalton	10	14659--76-----D
25	50,000 Dalton	10	14650--25-----D
47	50,000 Dalton	10	14650--47-----D
76	50,000 Dalton	10	14650--76-----D
25	300,000 Dalton	10	14679--25-----D
47	300,000 Dalton	10	14679--47-----D
76	300,000 Dalton	10	14679--76-----D

#### Cellulose Triacetate Membrane Filters, Type 145

Diameter in mm	MWCO	Qty/Pk	Prod. No.
25	5,000 Dalton	10	14529--25-----D
47	5,000 Dalton	10	14529--47-----D
25	10,000 Dalton	10	14539--25-----D
47	10,000 Dalton	10	14539--47-----D
50	10,000 Dalton	10	14539--50-----D
25	20,000 Dalton	10	14549--25-----D
43	20,000 Dalton	10	14549--43-----D
47	20,000 Dalton	10	14549--47-----D
47	20,000 Dalton	100	14549--47-----N
63	20,000 Dalton	10	14549--63-----D
25	30,000 Dalton	10	14459--25-----D
76	30,000 Dalton	10	14459--76-----D

#### Hydrosart Membrane Filters, Type 144

Diameter in mm	MWCO	Qty/Pk	Prod. No.
25	5,000 Dalton	10	14429--25-----D
44	5,000 Dalton	10	14429--44-----D
63	5,000 Dalton	10	14429--63-----D
76	5,000 Dalton	10	14429--76-----D
25	10,000 Dalton	10	14439--25-----D
47	10,000 Dalton	10	14439--47-----D
63	10,000 Dalton	10	14439--63-----D
76	10,000 Dalton	10	14439--76-----D
25	30,000 Dalton	100	14459--25-----D
47	30,000 Dalton	10	14459--47-----D
63	30,000 Dalton	10	14459--63-----D
76	30,000 Dalton	10	14459--76-----D



☐ DNA Concentration

35 Vivacon® 500

38 Vivacon® 2

## ■ Vivacon® 500

For DNA Sample Desalting and Concentration



### Reproducible DNA and Protein Sample Desalting and Concentration

Vivacon® 500 centrifugal concentrators offer the optimal solution for DNA and protein concentration and buffer exchange applications. For optimal performance with very dilute samples, e.g. forensic samples, Vivacon® 500 is equipped with the patented regenerated cellulose membrane Hydrosart®.

High recoveries and excellent reproducibilities are paired with convenience offered by molecular weight cut-off printed on individual devices.

The possibility of a re-spin after sample processing assures complete concentrate recovery which is especially important when working with low sample concentrations.

### New: Vivacon® 500-PCR Grade

When using DNA amplification technologies, any traces of DNA originating from the equipment have to be eliminated.

Vivacon® 500-PCR Grade units are treated with ethylene oxide (ETO) in a validated process in order to deactivate all traces of DNA that might interfere with subsequent amplification procedures.

References: K. Shaw et al., Int. J. Legal Med. (2008) 122: 29–33

### □ Technical Specifications

#### Vivacon® 500

Concentrator capacity	Fixed angle rotor	0.5 ml
Dimensions	Total length (Concentration)	45 mm
	Total length (back spin)	47.5 mm
	Width	12.4 mm
	Active membrane area	0.32 cm <sup>2</sup>
	Hold up volume of membrane and support	< 5 µl
	Dead stop volume (40° rotor)	5 µl
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polypropylene
	Membrane	Hydrosart®

### □ Conversion Table for Hydrosart® MWCO to Nucleotide Cut-Off

Membrane	MWCO	Double-Stranded Nucleotide Cut-Off (bp)
Hydrosart®	2 kDa	> 10
Hydrosart®	10 kDa	> 30
Hydrosart®	30 kDa	> 50
Hydrosart®	50 kDa	> 300
Hydrosart®	100 kDa	> 600
Cellulose Acetate	125 kDa	> 650

### Ordering Information

<b>Vivacon® 500</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
2,000 MWCO	25	VN01H91
2,000 MWCO	100	VN01H92
10,000 MWCO	25	VN01H01
10,000 MWCO	100	VN01H02
30,000 MWCO	25	VN01H21
30,000 MWCO	100	VN01H22
50,000 MWCO	25	VN01H31
50,000 MWCO	100	VN01H32
100,000 MWCO	25	VN01H41
100,000 MWCO	100	VN01H42
125,000 MWCO	25	VN01H81
125,000 MWCO	100	VN01H82
125,000 MWCO	500	VN01H83
<b>Vivacon® 500</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
Sample Kit L (4 units each of 2, 10, 30 K)	12	VN01HL12
Sample Kit H (4 units each of 30, 50, 100 K)	12	VN01HH12
<b>Vivacon® 500-PCR Grade Sample Pack</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
30,000 MWCO	4	VN01H2SETO
50,000 MWCO	4	VN01H3SETO
100,000 MWCO	4	VN01H4SETO

<b>Vivacon® 500-PCR Grade</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
30,000 MWCO	25	VN01H21ETO
30,000 MWCO	100	VN01H22ETO
30,000 MWCO	500	VN01H23ETO
50,000 MWCO	25	VN01H31ETO
50,000 MWCO	100	VN01H32ETO
50,000 MWCO	500	VN01H33ETO
100,000 MWCO	25	VN01H41ETO
100,000 MWCO	100	VN01H42ETO
100,000 MWCO	500	VN01H43ETO
125,000 MWCO	25	VN01H81ETO
125,000 MWCO	100	VN01H82ETO
125,000 MWCO	500	VN01H83ETO
<b>Accessories</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
Tubes	100	VNCT01

## Vivacon® 2

For DNA Sample Desalting and Concentration



### Reproducible DNA Sample Desalting and Concentration

Vivacon® 2 centrifugal concentrators offer the optimal solution for DNA and protein concentration and buffer exchange applications. For optimal performance with very dilute samples, e.g. forensic samples, Vivacon® 2 is equipped with the patented regenerated cellulose membrane Hydrosart®.

High recoveries and excellent reproducibilities are paired with convenience offered by volume graduation and molecular weight cut-off printed on individual devices.

The possibility of a re-spin after sample processing assures complete concentrate recovery which is especially important when working with low sample concentrations.

### New: Vivacon® 2-PCR Grade

Vivacon® 2-PCR Grade units are treated with ethylene oxide (ETO) in a validated process in order to deactivate all traces of DNA that might interfere with subsequent amplification procedures.

### Technical Specifications

Concentrator capacity	Fixed angle rotor	2 ml
Dimensions	Total length (Concentration)	125 mm
	Total length (Back-spin)	115 mm
	Width	16 mm
	Active membrane area	0.95 cm <sup>2</sup>
	Hold-up volume membrane and support	10 µl
	Dead stop volume (25° rotor)	55 µl
Materials of construction	Body	Polycarbonate
	Filtrate vessel	Polypropylene
	Back spin vial	Polypropylene
	Concentrator cap	Polypropylene
	Membrane	Hydrosart®

### Conversion Table for Hydrosart® MWCO to Nucleotide Cut-Off

Membrane	MWCO	Double-Stranded Nucleotide Cut-Off (bp)
Hydrosart®	2 kDa	> 10
Hydrosart®	10 kDa	> 30
Hydrosart®	30 kDa	> 50
Hydrosart®	50 kDa	> 300
Hydrosart®	100 kDa	> 600
Cellulose Acetate	125 kDa	> 650

## □ Ordering Information

<b>Vivacon® 2</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
2,000 MWCO	25	VN02H91
2,000 MWCO	100	VN02H92
2,000 MWCO	500	VN02H93
10,000 MWCO	25	VN02H01
10,000 MWCO	100	VN02H02
10,000 MWCO	500	VN02H03
30,000 MWCO	25	VN02H21
30,000 MWCO	100	VN02H22
30,000 MWCO	500	VN02H23
50,000 MWCO	25	VN02H31
50,000 MWCO	100	VN02H32
50,000 MWCO	500	VN02H33
100,000 MWCO	25	VN02H41
100,000 MWCO	100	VN02H42
100,000 MWCO	500	VN02H43
125,000 MWCO	25	VN02H81
125,000 MWCO	100	VN02H82
125,000 MWCO	500	VN02H83
<b>Vivacon® 2-PCR Grade</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
30,000 MWCO	25	VN02H21ETO
30,000 MWCO	100	VN02H22ETO
30,000 MWCO	500	VN02H23ETO
50,000 MWCO	25	VN02H31ETO
50,000 MWCO	100	VN02H32ETO
50,000 MWCO	500	VN02H33ETO
100,000 MWCO	25	VN02H41ETO
100,000 MWCO	100	VN02H42ETO
100,000 MWCO	500	VN02H43ETO
125,000 MWCO	25	VN02H81ETO
125,000 MWCO	100	VN02H82ETO
125,000 MWCO	500	VN02H83ETO
<b>Vivacon® 2-PCR Grade Sample Pack</b>	<b>Qty/Pk</b>	<b>Prod. No.</b>
30,000 MWCO	4	VN02H2SETO
50,000 MWCO	4	VN02H3SETO
100,000 MWCO	4	VN02H4SETO




## ☐ Protein Purification

41 Vivapure® Ion Exchange Protein Purification Products

43 Vivapure® Mini & Maxiprep Purification Kits



## Ion Exchange Protein Purification Products

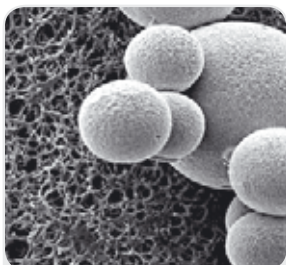

Chromatography gel beads (right) are shown on top of a membrane adsorber in this SEM picture. The membrane adsorber pores are over 50 × larger than bead pores.

### Fast and Easy-to-Use Spin Columns

Vivapure Ion Exchange (IEX) spin columns are centrifugal devices, incorporating Sartobind Membrane Adsorber technology as their chromatography matrix. Vivapure IEX spin columns make protein purification as easy as filtration. The devices are ready-to-use and do not bear the risk of running dry. For many protein purification applications, they can replace time-consuming and tedious column chromatography.

The rapid 1–2–3 bind–wash–elute protocol especially lends itself to screening applications, where many different samples are processed in parallel.

### The Sartobind Membrane Adsorber Matrix

Sartobind IEX membrane adsorbers are based on stabilized regenerated cellulose and display a microporous structure with a pore size of > 3 µm, which is orders of magnitude larger than conventional chromatographic gel materials. This allows molecules to be transported to the ligands immobilized on the membrane adsorber by convective flow, leading to very high flow rates.

In contrast to that, gel chromatography is slowed down due to diffusion limitations, as the molecules need to enter the small bead pores in order to be bound by the ligands. The porous membrane adsorber enables fast, reproducible and scalable protein purification.

### Fast and Simple to Use Spin Columns

- Devices are ready to use
- Make protein purification as simple as filtration

### Reproducible Results

- No column packing necessary devices are ready to use
- Membrane adsorber spin columns cannot crack or run dry

### Centrifugal Devices


- Offer the possibility of working in parallel

### Low Bed Volume

- Small membrane adsorber bed volumes allow working with lower buffer amounts, leading in concentrated elution fractions

### Up-Scalable Product Range

- Process scale modules are available with the same Sartobind IEX membrane adsorber matrix

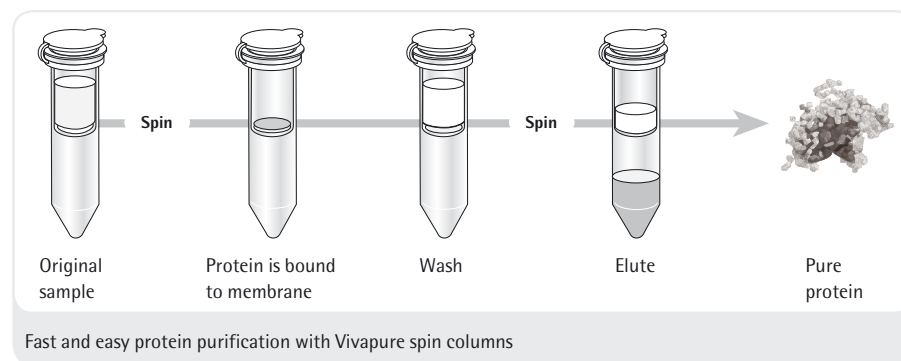



Vivapure Mini-400 | 500 µl  
Binding capacity: 1–4 mg



Vivapure Maxi-19 | 20 ml  
Binding capacity: 15–80 mg

## Technical Specifications

### Available Formats

Vivapure® IEX Products	Application
Vivapure Mini Spin Columns	<ul style="list-style-type: none"> <li>– Sample fractionation</li> <li>– Purification condition scouting</li> <li>– Small scale purification</li> </ul>
Vivapure Maxi Spin Columns	<ul style="list-style-type: none"> <li>– Large scale sample fractionation</li> <li>– One step protein purification   concentration</li> <li>– Polishing of his-tagged protein</li> </ul>

### Membrane Availability

Functional Groups	Ion Exchanger Type
Sulphonic acid (S)	Strong acidic cation exchanger: $\text{R-CH}_2\text{-SO}_3^-\text{Na}^+$
Quaternary ammonium (Q)	Strong basic anion exchanger: $\text{R-CH}_2\text{-N}^+\text{-(CH}_3)_3\text{Cl}^-$
Diethylamine (D)	Weak basic anion exchanger: $\text{R-CH}_2\text{-NH}^+\text{-(CH}_2\text{H}_5)_2$

### Performance Characteristics

Vivapure Spin Columns	Protein Binding Capacity* (mg)	Max. Volume per Centrifuge Run Using a Swing-Out Rotor (ml)	Max. Volume per Centrifuge Using a Fixed Angle Rotor Run (ml)
Vivapure Mini H	4	0.4	
Vivapure Maxi H	60–80	19	10.5

\* Actual yields depend on specific protein sample and selected pH and salt conditions. Yields established using 1 mg/ml BSA in 25 mM Tris/HCL pH 8.0 with Vivapure Q & D spin columns and 1 mg/ml cytochrome c in 25 mM sodium acetate buffer pH 5.5 with Vivapure S spin columns.

## Ordering Information

Description	Spin Columns	Centrifuge Tubes	Prod. No.
<b>Vivapure Mini Ion Exchange Spin Columns (up to 0.5 ml)</b>			
Vivapure Mini S&Q H starter kit	16	32	VS-IX01SQ16
Vivapure D Mini H	24	48	VS-IX01DH24
Vivapure Q Mini H	24	48	VS-IX01QH24
Vivapure S Mini H	24	48	VS-IX01SH24
<b>Vivapure Maxi Ion Exchange Spin Columns (up to 20 ml)</b>			
Vivapure D Maxi H	8	16	VS-IX20DH08
Vivapure Q Maxi H	8	16	VS-IX20QH08
Vivapure S Maxi H	8	16	VS-IX20SH08

# Vivapure® Mini & Maxiprep Purification Kits

For a Fast Antibody and His-Tagged Protein Purification



## Rapid Purification with High Yields

Vivapure® miniprep and maxiprep kits are spin column based kits for fast and effective purification of His-tagged proteins and antibodies.

Spin columns have the advantage of speed over gravity drip columns and batch protocols.


With the patented FlowGo regulator the sample residence time is extended to assure adequate sample binding to the resin. Due to this, Vivapure® miniprep and maxiprep spin column kits combine the merits of spin columns and gravity drip columns resulting in rapid purification with up to 95% protein recovery and purity.

All spin columns can conveniently be used in a centrifuge. For processing larger sample volumes, e. g. from diluted cell culture supernatants, the Vivapure® maxiprep spin columns can additionally be run with a peristaltic pump collar (VS-PPCSC).


The Vivapure® miniprep and maxiprep columns come in a variety of different configurations for your convenience. They can be purchased as completely ready-to-use kits with buffers and ultrafiltration devices but just as well as stand alone spin columns in small packs or large packs for frequent users.



## Vivapure® Miniprep

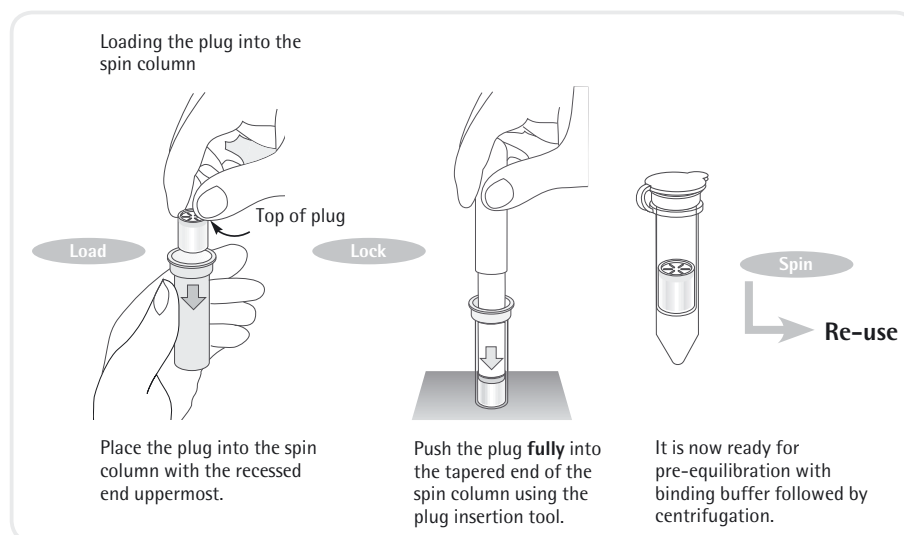

## Vivapure® Maxiprep

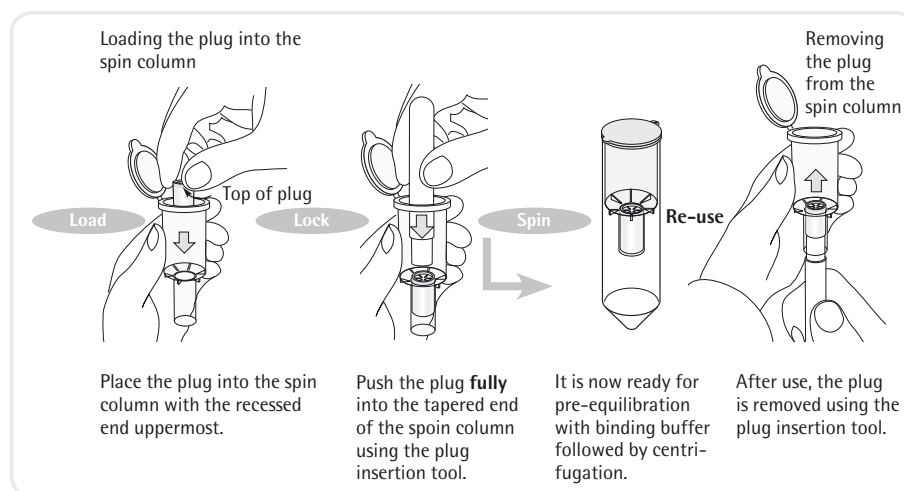

## Technical Specifications

### Protein A & G for Antibody Purification

Protein A & G Miniprep	Centrifuge
Sample size	0.65 ml
Typical Binding Capacity	1 mg IgG/column
Number of re-uses	3

Protein A & G Maxiprep	Centrifuge <sup>1</sup>
Sample size	20 ml
Typical Binding Capacity	20 mg IgG/ column
Number of re-uses	5

<sup>1</sup> Use the peristaltic pump accessory (VS-PPCSC) for larger volumes

## Ordering Information

Product Name	Qty/Pk	Prod. No.
Vivapure A Starter Pack*	2 miniprepA	VS-ARSTPKA2
Vivapure miniprepA Kit*	16 miniprepA	VS-ARAMINIK
Vivapure miniprepA Bulk Pack	48 miniprepA	VS-ARAMINIB
Vivapure maxiprepA Kit*	4 maxiprepA	VS-ARAMAXIK
Vivapure maxiprepA Bulk Pack	12 maxiprepA	VS-ARAMAXIB
Vivapure A Buffer Pack		VS-ARABUFPK
Vivapure G Starter Pack*	2 miniprepG	VS-ARSTPKG2
Vivapure miniprepG Kit*	16 miniprepG	VS-ARGMINIK
Vivapure miniprepG Bulk Pack	48 miniprepG	VS-ARGMINIB
Vivapure maxiprepG Kit*	4 maxiprepG	VS-ARGMAXIK
Vivapure maxiprepG Bulk Pack	12 maxiprepG	VS-ARGMAXIB
Vivapure G Buffer Pack		VS-ARGBUFPK
Sealing Cap & Peristaltic Pump Collor	1	VS-PPCSC

\* including UF-concentrators and buffers

## □ Technical Specifications

### Protein MC Miniprep

Kits	Centrifuge
Sample size	0.65 ml
Typical Binding	1 mg
Capacity	His-tagged protein
Number of re-uses	2

### Protein MC Maxiprep


Kits	Centrifuge <sup>1</sup>
Sample size	20 ml
Typical Binding	10 mg
Capacity	His-tagged protein
Number of re-uses	2

<sup>1</sup> Use the peristaltic pump accessory (VS-PPCSC) for larger volumes

## □ Ordering Information

Product Name	Qty/Pk	Prod. No.
Vivapure metal cheleate Starter Pack*	4	VS-MCST04
Vivapure miniprepMC Kit*	24	VS-MCMINI24
Vivapure miniprepMC Bulk Pack	72	VS-MCMINIB
Vivapure maxiprepMC Kit*	8	VS-MCMAXIK
Vivapure maxiprepMC Bulk Pack	24	VS-MCMAXIB
Vivapure metal chelate Buffer Pack		VS-MCBUFPK

\* including UF-concentrators and buffers

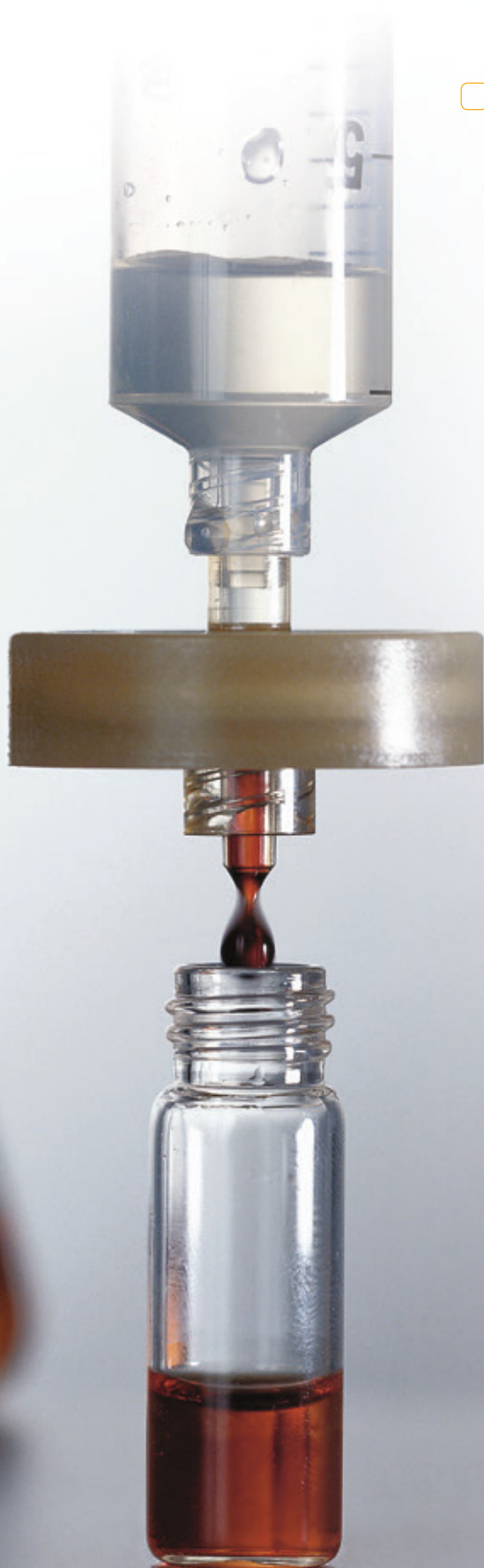

## ☐ Virus Purification and Concentration

47 Vivapure® Virus Purification and Concentration Kits

48 Adenovirus Purification with AdenoPACK Kits

50 Lentivirus Purification with LentiSELECT Kit

## ■ Vivapure<sup>®</sup> Virus Purification and Concentration Kits



Recombinant virus vectors are the preferred method for a wide range of gene delivery applications. Especially adenovirus type 5 and VSV-G pseudotyped lentivirus are two frequently utilized viral vectors for in vitro and in vivo applications.

### **Recombinant Adenovirus Vectors**

Recombinant adenovirus vectors are versatile tools in research and therapeutic applications for gene transfer and protein expression in cell lines that have low transfection efficiency with liposomes. After entering cells, the virus remains epichromosomal (i.e. does not integrate into the host chromosome, leaving the host genome unaffected). The delivery of RNAi into cells is becoming a major application for adenovirus vectors.

### **Lentivirus Vectors**

Lentivirus vectors are frequently used in gene transfer studies, due to their ability of gene transfer and integration into dividing and non-dividing cells. The pseudotyped envelope with vesicular stomatitis virus envelope G (VSV-G) protein broadens their target cell range. Lentiviral vectors have been shown to deliver genes into cell types (e.g. neurons, lymphocytes and macrophages) which other retrovirus vectors could not be used for. The lentivirus vector is increasingly used to integrate siRNA efficiently in a wide variety of cell lines and primary cells, both in vitro and in vivo.

### **Rapid Virus Purification by Membrane Chromatography**

The Sartobind<sup>®</sup> ion exchange membrane adsorber technology used in AdenoPACK and LentiSELECT is unique in its capability to efficiently and rapidly capture and recover large virus particles. When compared to chromatography media, membrane adsorbers provide large 3000 nm pores allowing unrestricted access and recovery of virus from the charged adsorber surface. Convective flow through the syringe filter devices provides high-speed separations not possible with traditional chromatography, cesium chloride density gradients and ultracentrifugation methods. Our membrane adsorbers with porous matrices, high capacities, low differential pressures, high flow rates and low unspecific adsorption show an excellent performance in small scale virus purification. Additionally, they are also scalable and confirm to cGMP facilities to large volume, high performance separation, reducing the processing time by a factor of 10 in the final process.

## ■ Adenovirus Purification with Vivapure AdenoPACK Kits

### AdenoPACK 20|100|500

The AdenoPACK adenovirus purification and concentration kits offer researchers who need to recover up to  $3 \times 10^{13}$  purified recombinant adenovirus particles for in-vitro transfection a fast, safe and easy to use solution. The kits include all reagents and devices necessary for clarification, purification and concentration of adenovirus type 5 from HEK293 cell cultures in only two hours. These straight forward kits replace time-consuming and labor-intensive 48 hour CsCl density gradients.

AdenoPACK kits are offered as AdenoPACK 20, AdenoPACK 100 and AdenoPACK 500, for the purification and concentration of adenovirus type 5 from 20 ml to 500 ml cell culture, leading to  $1 \times 10^{11}$ – $3 \times 10^{13}$  purified viral particles. For each sample volume, the most convenient handling method is offered for ultimate convenience.

To this end, preparations using AdenoPACK 20 are pursued in spin column format in a centrifuge, AdenoPACK 100 is a manually operated kit in syringe filter format\*, and AdenoPACK 500 is a pump driven kit.

\* Vivapure® AdenoPACK 100 can optionally be operated with a laboratory pump and an infusion pump, for which protocols are provided on our web page [www.sartorius-stedim.com](http://www.sartorius-stedim.com). Additionally, the tubes and adaptors needed for these operation modes can be ordered.

### AdenoPACK Advantages

#### Fast and Easy Virus Purification

- Purification completed in 2 hours
- Convenient, over 10× faster alternative to CsCl density gradient

#### Quantitative Yields

- In contrast to CsCl density gradient, the complete cell culture is used for virus purification and not only the viral pellet

#### Flexible Product Range

- Applicable from initial construct screening to large scale virus production

#### Complete Kit

- Including filtration devices, AdenoPACK units for virus purification, Vivaspin® and all buffers

#### Low Endotoxin Levels

- High cell viability and infection rates due to endotoxin levels of < 0.025 EU/ml

## □ Technical Specifications

### Adenovirus Purification Kit Specifications

Product	AdenoPACK 20	AdenoPACK 100	AdenoPACK 500
Sample Size	20 ml cell culture	20–200 ml of cell culture	500 ml of cell culture
Number of Purifications	6 × 20 ml	2 × 20–60 ml 1 × 200 ml	1 × 500 ml
Virus Particle (VP) per ml	Typically up to $1 \times 10^{11}$ – $10^{12}$	Typically up to $1 \times 10^{13}$	Typically up to $3 \times 10^{13}$
VP/IU	50–100	20–50	20–50
Processing time	Typically one hour	Typically two hours	
Endotoxin level	<0.025 EU/ml	<0.025 EU/ml	<0.025 EU/ml



□ Ordering Information

**Vivapure® AdenoPACK 20**

Vivapure® AdenoPACK 20	VS-AVPQ020
Vivapure® AdenoPACK 20 RT*	VS-AVPQ022

\* AdenoPACK 20 RT does not contain Benzonase®



**Vivapure® AdenoPACK 100**

Vivapure® AdenoPACK 100	VS-AVPQ101
Vivapure® AdenoPACK 100 RT*	VS-AVPQ102

**AdenoPACK 100 Accessories**

VS-AVPA001	Pump tubing set for Vivapure AdenoPACK 100
------------	---

\* AdenoPACK 100 RT does not contain Benzonase®



**Vivapure® AdenoPACK 500**

Vivapure® AdenoPACK 500	VS-AVPQ501
Vivapure® AdenoPACK 500 RT*	VS-AVPQ502

\* AdenoPACK 500 RT does not contain Benzonase®

## ■ Lentivirus Purification with Vivapure® LentiSELECT Kit

### LentiSELECT 40|500|1000

The LentiSELECT lentivirus purification and concentration kits offer researchers who need to recover up to  $5 \times 10^9$  infective lentivirus particles per ml for invitro transfection or animal studies a fast and easy to use solution.

These straight forward kits replace time-consuming ultracentrifugation protocols, which typically take approximately one day for large sample volumes, thus reducing the purification time to only a few hours.

LentiSELECT kits are offered as LentiSELECT 40, LentiSELECT 500 and LentiSELECT 1000 for the purification and concentration of VSV-G pseudotyped lentivirus from 40 ml to 1000 ml cell culture, leading to  $8 \times 10^8$  –  $1 \times 10^{10}$  purified infective particles. For each sample volume, the most convenient handling method is offered. To this end, 40 ml sample volumes are processed manually with LentiSELECT 40, while LentiSELECT 500 and 1000 are pump driven kits.

### LentiSELECT Advantages

#### Fast and Easy Virus Purification

- Purification completed in under one to six hours, depending on sample volume
- Kit as easy to use as filtration

#### No Need for Expensive Instruments

- Lentivirus purification with LentiSELECT is independent of equipment such as ultracentrifuges

#### High Virus Purity

- Achieve pure virus due to a chromatography purification for your experiments instead of a crude and variable cell culture supernatant pellet

#### Optimal for Multiple Virus Construct Screening

With LentiSELECT 40, four purification runs can be conducted in parallel with one kit

#### Complete Kits

- Including LentiSELECT units for virus purification, Vivaspins for concentration| buffer exchange and all buffers and syringes necessary

#### Low Endotoxin Levels

- High cell viability and infection rates due to endotoxin levels of  $< 0.025$  EU/ml

### □ Technical Specifications

Product	LentiSELECT 40	LentiSELECT 500	LentiSELECT 1000
Sample Size	40 ml cell culture	500 ml of cell culture	1000 ml of cell culture
Number of Purifications	4 × 40 ml	1 × 500 ml	1 × 1000 ml
Virus Particle (VP) per ml	Typically up to $3 \times 10^9$	Typically up to $2-5 \times 10^9$ *	Typically up to $4-6 \times 10^{13}$
VP/IU	5–15	5–15	20–50
Processing time	Typically up to 45 min	Typically up to 3 hours	Typically up to 6 hours
Endotoxin level	$<0.025$ EU/ml	$<0.025$ EU/ml	$<0.025$ EU/ml



### □ Ordering Information

#### **Vivapure® LentiSELECT 40**

Vivapure LentiSELECT 40

VS-LVPQ040



#### **Vivapure® LentiSELECT 500**

Vivapure LentiSELECT 500

VS-LVPQ500



#### **Vivapure® LentiSELECT 1000**

Vivapure LentiSELECT 1000

VS-LVPQ1000





## ■ Filtration Devices

### Table of Contents

54 Minisart® Syringe Filters

65 Sartolab® P20 and Sartolab® P20 Plus

68 Sartolab® 150v

69 Sartolab® RF | BT

72 Chemical Compatibility

## ■ Minisart® Syringe Filters

Removal of Particles and Microorganisms from Liquids and Gases



### **Sample Preparation** **HPLC | UHPLC | Analytics**

Elimination of particles from your samples prior to HPLC or other chromatographic analysis is essential in order to maintain the integrity of your chromatography column and to maximize its operating life.


Minisart® syringe filters consist of a PP housing and membrane components featuring maximum chemical compatibility and minimum extractables to ensure excellent results. See page 56.



### **Filtration of Aqueous Liquids** **Clarification | Sterile Filtration**

For clarification and sterilization of liquids, filtration is the optimal method. It removes all microorganisms and particles reliably, without any effects on the ingredients, due to adsorption or decomposition.

For optimal results Minisart® filters made of MBS housing provide high flow rates and lowest adsorption characteristics. The MBS housing is color-coded for an easy pore size identification. See page 59.

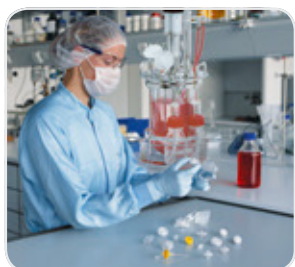

### **Medical Use and Sterile Venting** **Special Applications**


Minisart® syringe filters are available in a broad range of pore sizes. They are ideal for clarification of liquids laden with particles, e.g. for preparation of pharmaceuticals or infusion solutions. For sterilization and removal of particles from air and other

gases, syringe filters are optimal for sterile venting of containers, bioreactors, fermenters and tubing systems in medical devices. Many Minisart® syringe filters have the CE mark of conformity (European directive) and are available in a wide choice of membranes, connectors and housing materials. See page 62.

## Selection Guide

1

### Sample Composition

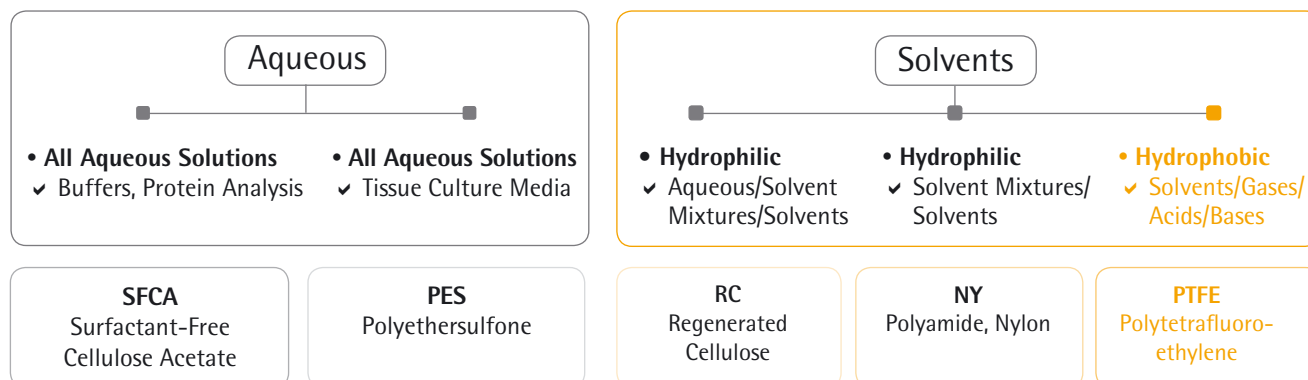

2

### Sample Volumes



3

### Pore Sizes

Sterilization		Sample Preparation · Clarification · Particle Removal					Prefiltration
<b>0.1 µm</b>	<b>0.2 µm</b>	<b>0.45 µm</b>	<b>0.65 µm</b>	<b>0.8 µm</b>	<b>1.2 µm</b>	<b>5 µm</b>	<b>GF</b> (Glass Fibre)
Small Bacteria	UHPLC, etc.	HPLC, etc.	Particles	Particles	Particles	Large	Glass Prefilter
Mycoplasma	(Columns	(Columns	Yeast Cells	Yeast Cells	Yeast Cells	Particles	Glass+Membrane
Colloids >0.1	< 3µm	> 3µm			Platelets	Rubber Grit	Highly Particle-laden Samples
	Particles)	Particles				Cells	
	Bacteria						

## Applications

Type of Filtration	1st Choice	Rec. Alternatives
HPLC, UHPLC, LCMS, IC, GC – Sample Preparation	RC	PTFE   NY
Undiluted Organic Solvents	RC   PTFE	NY
Protein Analysis, Samples with Biomolecules, Buffers	SFCA	RC   PES
Tissue Culture Media	PES	SFCA   RC
Highly Particle-laden Samples, Organic Solvents	GF   NY	–
Highly Particle-laden Samples, Aqueous Solutions	GF   SFCA	GF   NY

pH 1-14: PTFE | GF, pH 3-14: RC | NY | PTFE | GF, pH 3-12: RC | PES | NY | PTFE | GF,  
pH 4-8: SFCA | RC | PES | NY | PTFE | GF\*

\* Compatibility tested with a contact time of 24 hours at 20°C

## Sample Preparation for Analytics

Reliable Removal of Particles and Microorganisms from Liquids to Protect Your Columns and Instruments



### Minisart® Features

- 15 mm and 25 mm with pore size and type of membrane identified
- 4 mm with color-coded package
- Low adsorption of analytes
- Maximum chemical compatibility
- Minimum extractables or leachables
- Superior flow rate
- Low hold-up volumes
- 100% camera system control
- Bidirectional use
- Certified quality


## Key Specifications

Diameter	Pore Size	Membranes	Housing	Connector Outlet	Hold-Up Volume	Filter Area
25 mm	0.2   0.45 µm	RC   PTFE   NY	PP*	Male Luer Slip	100-200 µl**	4.8 cm <sup>2</sup>
15 mm	0.2   0.45 µm	RC   PTFE   NY	PP*	Male Luer Slip   Male Spike	30-100 µl**	1.7 cm <sup>2</sup>
4 mm	0.2   0.45 µm	RC   PTFE	PP*	Male Luer Slip	5-10 µl**	0.07 cm <sup>2</sup>
GF Prefilter	1.2+0.2   0.45	GF+NY	PP*	Male Luer Slip	250-500 µl**	4.8 cm <sup>2</sup>


\* PP = Polypropylene \*\*Hold-up volume after air purge. Volumes can vary depending on membrane and liquid used.

### HPLC Certification\*

#### Methanol

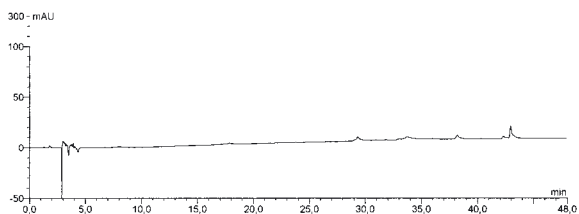


#### Filtrate Methanol

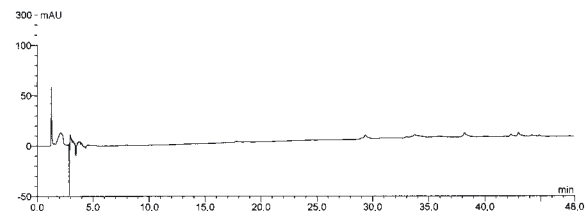


#### Methanol | Filtrate Methanol

Channel: UV\_VIS\_2      Run Time (min.): 58.00  
Wavelength: 214      Injection Volume: 100.0  
Bandwidth: 1

#### Acetonitrile

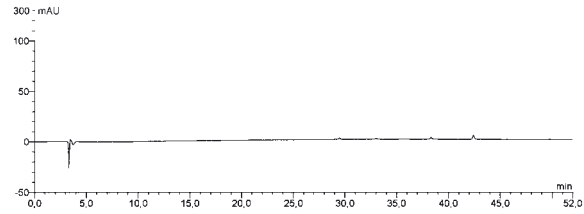

#### Filtrate Acetonitrile

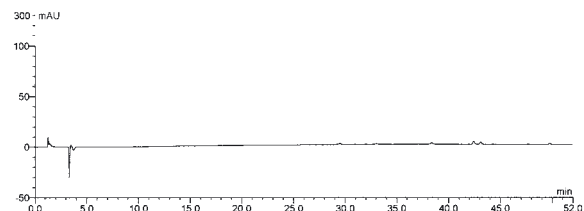

#### Acetonitrile | Filtrate Acetonitrile

Channel: UV\_VIS\_4      Run Time (min.): 58.00  
Wavelength: 280      Injection Volume: 100.0  
Bandwidth: 1

\* Minisart® used with methanol | water and acetonitrile | water did not show any artifacts or interference peaks in the range of 200–300 nm

## Minisart® Syringe Filters – Sample Preparation Chromatography

### Ordering Information

Type	Ø mm	Membrane	Housing	Pore Size	Connector Outlet	Color   Printing	Sterile*	Qty   Pk	Order No.
Minisart® RC (Regenerated Cellulose)									
25 mm	RC	PP	0.2 µm	Male Luer Slip	White, Printed	Yes	50	17764-----ACK	
25 mm	RC	PP	0.2 µm	Male Luer Slip	White, Printed	No	50	17764-----K	
25 mm	RC	PP	0.2 µm	Male Luer Slip	White, Printed	No	200	17764-----S	
25 mm	RC	PP	0.2 µm	Male Luer Slip	White, Printed	No	500	17764-----Q	
25 mm	RC	PP	0.45 µm	Male Luer Slip	White, Printed	No	50	17765-----K	
25 mm	RC	PP	0.45 µm	Male Luer Slip	White, Printed	No	200	17765-----S	
25 mm	RC	PP	0.45 µm	Male Luer Slip	White, Printed	No	500	17765-----Q	
15 mm	RC	PP	0.2 µm	Male Luer Slip	White, Printed	Yes	50	17761-----ACK	
15 mm	RC	PP	0.2 µm	Male Luer Slip	White, Printed	No	50	17761-----K	
15 mm	RC	PP	0.2 µm	Male Luer Slip	White, Printed	No	500	17761-----Q	
15 mm	RC	PP	0.45 µm	Male Luer Slip	White, Printed	No	50	17762-----K	
15 mm	RC	PP	0.45 µm	Male Luer Slip	White, Printed	No	500	17762-----Q	
4 mm	RC	PP	0.2 µm	Male Luer Slip	Blue Tray	No	50	17821-----K	
4 mm	RC	PP	0.2 µm	Male Luer Slip	Blue Tray	No	500	17821-----Q	
4 mm	RC	PP	0.45 µm	Male Luer Slip	Yellow Tray	No	50	17822-----K	
4 mm	RC	PP	0.45 µm	Male Luer Slip	Yellow Tray	No	500	17822-----Q	
Minisart® SRP (Hydrophobic PTFE)									
25 mm	PTFE	PP	0.2 µm	Male Luer Slip	White, Printed	Yes	50	17575-----ACK	
25 mm	PTFE	PP	0.2 µm	Male Luer Slip	White, Printed	No	50	17575-----K	
25 mm	PTFE	PP	0.2 µm	Male Luer Slip	White, Printed	No	200	17575-----S	
25 mm	PTFE	PP	0.2 µm	Male Luer Slip	White, Printed	No	500	17575-----Q	
25 mm	PTFE	PP	0.45 µm	Male Luer Slip	White, Printed	No	50	17576-----K	
25 mm	PTFE	PP	0.45 µm	Male Luer Slip	White, Printed	No	200	17576-----S	
25 mm	PTFE	PP	0.45 µm	Male Luer Slip	White, Printed	No	500	17576-----Q	
15 mm	PTFE	PP	0.2 µm	Male Spike	White, Printed	No	50	17558-----K	
15 mm	PTFE	PP	0.2 µm	Male Spike	White, Printed	No	500	17558-----Q	
15 mm	PTFE	PP	0.2 µm	Male Luer Slip	White, Printed	Yes	50	17573-----ACK	
15 mm	PTFE	PP	0.2 µm	Male Luer Slip	White, Printed	No	50	17573-----K	
15 mm	PTFE	PP	0.2 µm	Male Luer Slip	White, Printed	No	500	17573-----Q	
15 mm	PTFE	PP	0.45 µm	Male Spike	White, Printed	No	50	17559-----K	
15 mm	PTFE	PP	0.45 µm	Male Spike	White, Printed	No	500	17559-----Q	
15 mm	PTFE	PP	0.45 µm	Male Luer Slip	White, Printed	No	50	17574-----K	
15 mm	PTFE	PP	0.45 µm	Male Luer Slip	White, Printed	No	500	17574-----Q	
4 mm	PTFE	PP	0.2 µm	Male Luer Slip	Blue Tray	No	500	17844-----Q	
4 mm	PTFE	PP	0.45 µm	Male Luer Slip	Yellow Tray	No	50	17820-----K	
4 mm	PTFE	PP	0.45 µm	Male Luer Slip	Yellow Tray	No	500	17820-----Q	

Type	Ø mm	Membrane	Housing	Pore Size	Connector	Outlet	Color   Printing	Sterile*	Qty   Pk	Order No.
Minisart® NY (Nylon) & NY25 Plus (Glass Fiber 1.2 µm + Nylon)										
25 mm	Nylon	PP	0.2 µm	Male Luer Slip	White, Printed	Yes	50	17845-----ACK		
25 mm	Nylon	PP	0.2 µm	Male Luer Slip	White, Printed	No	500	17845-----Q		
25 mm	Nylon	PP	0.45 µm	Male Luer Slip	White, Printed	Yes	50	17846-----ACK		
25 mm	Nylon	PP	0.45 µm	Male Luer Slip	White, Printed	No	500	17846-----Q		
15 mm	Nylon	PP	0.2 µm	Male Luer Slip	White, Printed	No	50	1776B-----K		
15 mm	Nylon	PP	0.2 µm	Male Luer Slip	White, Printed	No	500	1776B-----Q		
15 mm	Nylon	PP	0.45 µm	Male Luer Slip	White, Printed	No	50	1776C-----K		
15 mm	Nylon	PP	0.45 µm	Male Luer Slip	White, Printed	No	500	1776C-----Q		
25 mm	Nylon+GF	PP	0.2 µm	Male Luer Slip	White, Printed	No	50	1784B-----K		
25 mm	Nylon+GF	PP	0.2 µm	Male Luer Slip	white, printed	No	500	1784B-----Q		
25 mm	Nylon+GF	PP	0.45 µm	Male Luer Slip	white, printed	No	50	1784C-----K		
25 mm	Nylon+GF	PP	0.45 µm	Male Luer Slip	white, printed	No	500	1784C-----Q		

\* Sterile Minisarts are individually packaged. If not stated otherwise, Minisarts have been sterilized by ethylene oxide. Minisarts NOT presterilized: RC, PTFE and Nylon can be sterilized by autoclaving at 121°C for 30 min/or by using ethylene oxide (EO).

## Technical Specifications

### Specifications for Minisart® RC, Minisart® SRP, Minisart® NY, 4 mm, 15 mm and 25 mm Membrane Diameter

Application limits	Max. recommended operating pressure 4.5 bar   65 psi
Housing burst pressure	6 bar   87 psi and higher
Max. temperature	121°C, 30 min (autoclave)
Bubble point	Min value with water (RC): $\geq 3.2$ bar   $\geq 47$ psi (0.2 µm), $\geq 2.0$ bar   $\geq 29$ psi (0.45 µm) Min value with isopropanol (SRP): $\geq 1.4$ bar   $\geq 20$ psi (0.2 µm), $\geq 0.9$ bar   $\geq 13$ psi (0.45 µm) Min value with water (NY): $\geq 3.0$ bar   $\geq 44$ psi (0.2 µm), $\geq 2.0$ bar   $\geq 29$ psi (0.45 µm)
Flow rate, 4 mm	Typical values for water at 3 bar – RC: 0.5 ml/min (0.2 µm), 1.5 ml/min (0.45 µm) Typical values for methanol at 1 bar – RC: 1.5 ml/min (0.2 µm), 3.0 ml/min (0.45 µm); SRP: 4.5 ml/min (0.45 µm) Typical values for air at 1 bar – SRP: 60 ml/min (0.45 µm)
Flow rate, 15 mm	Typical values for water at 1 bar – RC: 10 ml/min (0.2 µm), 30 ml/min (0.45 µm); NY: 20 ml/min (0.2 µm), 40 ml/min (0.45 µm) Typical values for methanol at 1 bar – RC: 55 ml/min (0.2 µm), 105 ml/min (0.45 µm); SRP: 55 ml/min (0.2 µm), 150 ml/min (0.45 µm) Typical values for air at 1 bar – SRP: 0.5 l/min (0.2 µm), 1.1 l/min (0.45 µm)
Flow rate, 25 mm	Typical values for water at 1 bar – RC: 60 ml/min (0.2 µm), 100 ml/min (0.45 µm); NY: 75 ml/min (0.2 µm), 130 ml/min (0.45 µm) Typical values for methanol at 1 bar – RC: 160 ml/min (0.2 µm), 325 ml/min (0.45 µm); SRP: 160 ml/min (0.2 µm), 260 ml/min (0.45 µm) Typical values for air at 1 bar – SRP: 1.2 l/min (0.2 µm), 1.8 l/min (0.45 µm)
Water penetration point	Minisart® SRP (hydrophobic PTFE) 4.0 bar (0.2 µm) or 3.0 bar (0.45 µm)

## Filtration of Aqueous Liquids

Filtration is the Optimal Method for Clarification and Sterilization of Liquids



### Minisart® Features

- Pore size + membrane color-coded
- Lowest adsorption
- Superior flow rate
- High total throughput
- Low hold-up volume
- Particulate-free
- Minimum extractables
- 100% camera system control
- PVC-free
- Choice of presterilized or non-sterile units
- Gamma irradiated or EO sterilized
- Fully validated, certified quality
- Bidirectional use


## Key Specifications

Diameter	Pore Size	Membranes	Housing	Connector Outlets	Hold-Up Volume	Filter Area
28 mm	0.1 µm	PES	MBS*	Male Luer Lock	100-150 µl**	6.2 cm <sup>2</sup>
28 mm	0.2 µm	SFCA   PES	MBS*	Male Luer Lock   Male Luer Slip	100-150 µl**	6.2 cm <sup>2</sup>
28 mm	0.45 µm	SFCA   PES	MBS*	Male Luer Lock   Male Luer Slip	100-150 µl**	6.2 cm <sup>2</sup>
28 mm	0.65 µm	SFCA	MBS*	Male Luer Slip	100-150 µl**	6.2 cm <sup>2</sup>
28 mm	0.8 µm	SFCA	MBS*	Male Luer Slip	100-150 µl**	6.2 cm <sup>2</sup>
28 mm	1.2 µm	SFCA	MBS*	Male Luer Slip	100-150 µl**	6.2 cm <sup>2</sup>
28 mm	5 µm	SFCA	MBS*	Male Luer Slip	100-150 µl**	6.2 cm <sup>2</sup>
GF Prefilter	1.2+0.2   0.45	GF+SFCA	MBS*	Male Luer Slip	230 µl**	6.2 cm <sup>2</sup>
GF Prefilter	1.2 µm	GF	MBS*	Male Luer Lock   Male Luer Slip	200 µl**	6.2 cm <sup>2</sup>


\* MBS = Methacrylate butadiene styrene polymerizate

\*\* Hold-up volume after air purge. Volumes can vary depending on membrane and liquid used.

### Sartorius Advantage: Water Flow at 1 bar (15.4 psi) ml/min

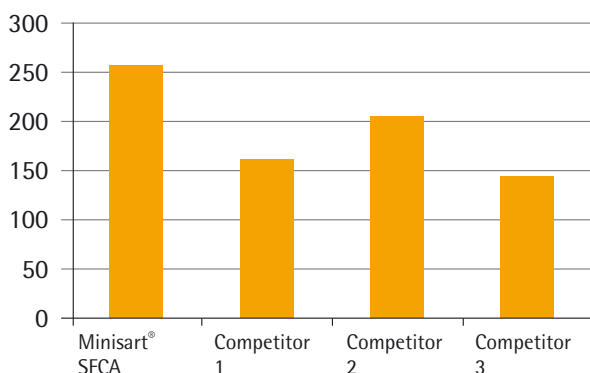

0.45 µm SFCA membrane compared with CA (Cellulose Acetate) or similar hydrophilic membranes

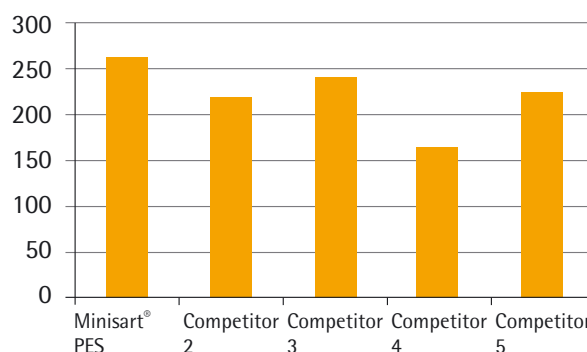

0.45 µm PES membrane compared with other PES membranes

## Minisart® Syringe Filters – Preparation of Aqueous Liquids

## □ Ordering Information

Type	Ø mm	Membrane	Housing	Pore Size	Connector Outlet	Color	Sterile*	Qty   Pk	Order No.
<b>Minisart® High Flow (PES – Polyethersulfone)</b>									
	28 mm	PES	MBS	0.1 µm	Male Luer Lock	Dark Red	Yes	50	16553-----K
	28 mm	PES	MBS	0.2 µm	Male Luer Lock	Royal Blue	Yes#	50	16532-----GUK
	28 mm	PES	MBS	0.2 µm	Male Luer Lock	Royal Blue	Yes	50	16532-----K
	28 mm	PES	MBS	0.2 µm	Male Luer Slip	Royal Blue	Yes	50	16541-----K
	28 mm	PES	MBS	0.2 µm	Male Luer Lock	Royal Blue	No	500	16532-----Q
	28 mm	PES	MBS	0.2 µm	Male Luer Slip	Royal Blue	No	500	16541-----Q
	28 mm	PES	MBS	0.45 µm	Male Luer Lock	Amber	Yes	50	16537-----K
	28 mm	PES	MBS	0.45 µm	Male Luer Lock	Amber	No	500	16537-----Q
	28 mm	PES	MBS	0.45 µm	Male Luer Slip	Amber	Yes#	50	16533-----GUK
	28 mm	PES	MBS	0.45 µm	Male Luer Slip	Amber	Yes	50	16533-----K
	28 mm	PES	MBS	0.45 µm	Male Luer Slip	Amber	No	500	16533-----Q
<b>Minisart® NML (SFCA – Surfactant-free Cellulose Acetate)</b>									
	28 mm	SFCA	MBS	0.2 µm	Male Luer Lock	Blue	Yes	50	16534-----K
	28 mm	SFCA	MBS	0.2 µm	Male Luer Lock	Blue	Yes#	50	16534-----GUK
	28 mm	SFCA	MBS	0.2 µm	Male Luer Lock	Blue	No	500	16534-----Q
	28 mm	SFCA	MBS	0.2 µm	Male Luer Slip	Blue	Yes	50	17597-----K
	28 mm	SFCA	MBS	0.2 µm	Male Luer Slip	Blue	No	500	17597-----Q
	28 mm	SFCA	MBS	0.45 µm	Male Luer Lock	Yellow	Yes	50	16555-----K
	28 mm	SFCA	MBS	0.45 µm	Male Luer Lock	Yellow	Yes#	50	16555-----GUK
	28 mm	SFCA	MBS	0.45 µm	Male Luer Lock	Yellow	No	500	16555-----Q
	28 mm	SFCA	MBS	0.45 µm	Male Luer Slip	Yellow	Yes	50	17598-----K
	28 mm	SFCA	MBS	0.45 µm	Male Luer Slip	Yellow	No	500	17598-----Q
	28 mm	SFCA	MBS	0.65 µm	Male Luer Slip	Pink	Yes	50	16569-----K
	28 mm	SFCA	MBS	0.8 µm	Male Luer Lock	Green	Yes	50	16592-----K
	28 mm	SFCA	MBS	0.8 µm	Male Luer Lock	Green	Yes#	50	16592-----GUK
	28 mm	SFCA	MBS	0.8 µm	Male Luer Lock	Green	No	500	16592-----Q
	28 mm	SFCA	MBS	1.2 µm	Male Luer Lock	Red	Yes	50	17593-----K
	28 mm	SFCA	MBS	1.2 µm	Male Luer Lock	Red	No	500	17593-----Q
	28 mm	SFCA	MBS	5 µm	Male Luer Lock	Brown	Yes	50	17594-----K
	28 mm	SFCA	MBS	5 µm	Male Luer Lock	Brown	No	500	17594-----Q

Type	Ø mm	Membrane	Housing	Pore Size	Connector Outlet	Color	Sterile*	Qty   Pk	Order No.
<b>Minisart® NML Plus (Glass Fiber 1.2 µm + SFCA)</b>									
28 mm		GF+SFCA	MBS	0.2 µm	Male Luer Lock	Blue	Yes	50	17823-----K
28 mm		GF+SFCA	MBS	0.2 µm	Male Luer Lock	Blue	No	500	17823-----Q
28 mm		GF+SFCA	MBS	0.45 µm	Male Luer Lock	Yellow	Yes	50	17829-----K
28 mm		GF+SFCA	MBS	0.45 µm	Male Luer Lock	Yellow	No	500	17829-----Q
28 mm		GF+SFCA	MBS	1.2 µm	Male Luer Lock	Red	No	500	17825-----Q
28 mm		GF	MBS	0.7 µm	Male Luer Lock	White	No	50	17824-----K
28 mm		GF	MBS	0.7 µm	Male Luer Lock	White	No	500	17824-----Q

\* Sterilized Minisarts are individually packaged. If not stated otherwise, Minisarts are sterilized by ethylene oxide.

#-Mark indicates sterilization by gamma irradiation.

Minisarts not presterilized: PES, SFCA, GF+SFCA and GF can be sterilized by ethylene oxide or gamma irradiation.

Would you like to filter solvents, acids or bases?

Do you need to filter small volumes?

Would you like to use PP housings and other membranes?

– Please see Minisart® RC, Minisart® NY or Minisart® SRP  
for the highest chemical compatibility; also available in 4 mm  
or 15 mm.

## Technical Specifications

### Specifications for Minisart® High Flow, Minisart® NML, Minisart® NML Plus

Application limits	Max. recommended operating pressure 4.5 bar   65 psi
Housing burst pressure	6 bar   87 psi and higher
Max. temperature	50°C
Bubble point	Min value with water ≥ 5 bar   ≥ 73 psi (0.1 µm), ≥ 3.2 bar   ≥ 46 psi (0.2 µm), ≥ 2.0 bar   ≥ 29 psi (0.45 µm), ≥ 1.3 bar   ≥ 19 psi (0.65 µm), ≥ 0.8 bar   ≥ 12 psi (0.8 µm), ≥ 0.7 bar   ≥ 10 psi (1.2 µm), ≥ 0.4 bar   ≥ 6 psi (5 µm), ≥ 0.5 bar   ≥ 7 psi (Minisart® GF)
Flow rate	Typical values for water at 1 bar   14.5 psi, 30 ml/min (0.1 µm), 60 ml/min (0.2 µm), 160 ml/min (0.45 µm), 250 ml (0.65 µm), 350 ml/min (0.8 µm), 400 ml/min (1.2 µm), 500 ml/min (5 µm), 450 ml/min (Minisart® GF)
Cytotoxicity	Detectably no inhibition with MRC-5 (human lung cells)



## Medical Use & Venting – Special Applications

Make Your Choice from a Broad Range of Pore Sizes, Materials and Formats

Minisart® in many configurations



### Minisart® Features

- Certified quality
- Fully validated
- 100% camera system control
- Integrity testable before or after use
- Low adsorption
- Minimum extractables
- Particulate-free
- PVC-free
- Presterilized or non-sterile
- Gamma irradiated or EO sterilized
- CE-marked types
- Various inlet connectors
- Various outlet connectors
- Bidirectional

### Key Specifications

Diameter	Pore Size	Membranes	Housing	Connector Outlets	Hold-Up Volume	Filter Area
28 mm	0.2   0.45 µm	SFCA	MBS*	Male Luer Lock   Male Luer Slip	100-150 µl***	6.2 cm <sup>2</sup>
28 mm	5 µm	SFCA	MBS*	Male Luer Lock	100-150 µl***	6.2 cm <sup>2</sup>
28 mm	0.1 µm	PES	MBS*	Male Luer Lock	100-150 µl***	6.2 cm <sup>2</sup>
26 mm	0.2   1 µm	PTFE	MBS*	Male Luer Lock   Male Luer Slip	100-150 µl***	5.3 cm <sup>2</sup>
26 mm	0.45 µm	Dome+PTFE	MBS*	Male Luer Slip (inlet and outlet)	100-150 µl***	5.3 cm <sup>2</sup>
25 mm	0.2 µm	PTFE	PP**	Male Luer Slip   Tube Connection	100-200 µl***	4.8 cm <sup>2</sup>
15 mm	0.2 µm	PTFE	MBS*	Male Luer Slip   Slip + Needle	100 µl***	1.7 cm <sup>2</sup>

\* MBS = Methacrylate butadiene styrene polymerizate

\*\* PP = Polypropylene

\*\*\* Hold-up volume after air purge. Volumes can vary depending on membrane and liquid used.



Biocompatibility  
Certificate Minisart® HY



Biocompatibility  
Certificate Minisart® NML



Declaration of Conformity  
Minisart®

## Minisart® Syringe Filters – Specials

 Ordering Information

Type	Ø mm	Membrane	Housing	Pore Size	Connector Outlet	Color	Sterile*	Qty   Pk	Order No.	Marked?
Minisart® NML (SFCA - Cellulose Acetate) Aqueous Filtration										
28 mm	SFCA	MBS	0.2 µm	Male Luer Lock	Blue	Yes	50	16534-----K	CE-marked	
28 mm	SFCA	MBS	0.2 µm	Male Luer Lock	Blue	Yes#	50	16534-----GUK	CE-marked	
28 mm	SFCA	MBS	0.2 µm	Male Luer Lock	Blue	No	500	16534-----Q	CE-marked	
28 mm	SFCA	MBS	0.2 µm	Male Luer Slip	Blue	Yes	50	17597-----K	CE-marked	
28 mm	SFCA	MBS	0.2 µm	Male Luer Slip	Blue	No	500	17597-----Q	CE-marked	
28 mm	SFCA	MBS	0.45 µm	Male Luer Lock	Yellow	Yes	50	16555-----K	CE-marked	
28 mm	SFCA	MBS	0.45 µm	Male Luer Lock	Yellow	Yes#	50	16555-----GUK	CE-marked	
28 mm	SFCA	MBS	0.45 µm	Male Luer Lock	Yellow	No	500	16555-----Q	CE-marked	
28 mm	SFCA	MBS	0.45 µm	Male Luer Slip	Yellow	Yes	50	17598-----K	CE-marked	
28 mm	SFCA	MBS	0.45 µm	Male Luer Slip	Yellow	No	500	17598-----Q	CE-marked	
28 mm	SFCA	MBS	5 µm	Male Luer Lock	Brown	Yes	50	17594-----K	CE-marked	
Minisart® Ophthalsart (SFCA - Cellulose Acetate) Aqueous Filtration										
28 mm	Ophthalsart	MBS	0.2 µm	Male Luer Slip	Pink	Yes	50	17528-----K	CE-marked	
Minisart® High Flow (PES - Polyethersulfone) Aqueous Filtration										
28 mm	PES	MBS	0.1 µm	Male Luer Lock	Dark Red	Yes	50	16553-----K		
Minisart® (PES - Polyethersulfone) Aqueous Filtration										
15 mm	PES	PP	0.2 µm	Male Luer Slip	White	Yes	50	1776D-----ACK		
Minisart® Air (Hydrophobic PTFE) Venting										
15 mm	PTFE	MBS	0.2 µm	Male Luer Slip	Yellow	No	500	1751A-----Q		
15 mm	PTFE	MBS	0.2 µm	Male Luer Slip + Needle	Yellow	Yes#	50	16596-----HNK		
Minisart® HY (hydrophobic PTFE) CE-marked Venting & Gas Filtration										
26 mm	PTFE	MBS	0.2 µm	Male Luer Lock	Clear	Yes	50	16596-----HYK	CE-marked	
26 mm	PTFE	MBS	0.2 µm	Male Luer Lock	Clear	No	500	16596-----HYQ	CE-marked	
26 mm	PTFE	MBS	0.2 µm	Male Luer Lock <sup>a</sup>	Clear	No	500	16599-----HYQ	CE-marked	
26 mm	PTFE	MBS	0.2 µm	Hose Barbs <sup>b</sup>	Clear	No	500	40078-----Q	CE-marked	
26 mm	PTFE	MBS	1 µm	Male Luer Lock	Clear	No	500	1659A-----HYQ		
26 mm	PTFE	MBS	1 µm	Hose Barbs <sup>b</sup>	Clear	No	500	1659B-----HYQ		

Type	Ø mm	Membrane	Housing	Pore Size	Connector	Outlet	Color	Sterile*	Qty   Pk	Order No.	Marked?
<b>Minisart® Acticosart with Dome Reservoir + Hydrophobic PTFE Venting &amp; Ultracleaning of Gases</b>											
	26 mm	active carbon	MBS	0.45 µm	Male Luer Slip <sup>a</sup>		Blue	No	500	17840-----Q	
<b>Minisart® SRP (Hydrophobic PTFE) CE-marked Venting &amp; Gas Filtration</b>											
	25 mm	PTFE	PP	0.2 µm	Male Luer Slip		White, Printed	Yes	50	17575-----ACK	CE-marked
	25 mm	PTFE	PP	0.2 µm	Hose Barb		White	No	500	1757A-----Q	

\* Sterilized Minisarts are individually packaged. If not stated otherwise, Minisarts are sterilized by ethylene oxide.

#-mark indicates sterilization by gamma irradiation

Minisarts NOT presterilized: SFCA can be sterilized by ethylene oxide or gamma irradiation. PTFE can be sterilized by ethylene oxide.

<sup>a</sup> Connector inlet: Male Luer slip (all other Minisarts have female luer lock inlets)

<sup>b</sup> Hose barbs, inlet and outlet, 5 mm diameter

Do you need additional pore sizes?

Would you like to use other membranes?

– Please refer to Minisart® NML

available in 0.65, 0.8 and 1.2 µm pore sizes; also available with PES membrane

Would you like to filter solvents, acids or bases?

Do you want to filter small volumes?

Would you like to use PP housings and other membranes?

– Please refer to Minisart® RC, Minisart® NY or Minisart® SRP for the highest chemical compatibility; also available in 4 mm and 15 mm

## Technical Specifications

### Minisart® Acticosart, Minisart® Air, Minisart® HY, Minisart® Ophthalmart

(Specifications for Minisarts® High Flow, NML, NML Plus, PES and SRP please see on pages 60 and 63)

Application limits	Max. recommended operating pressure: Acticosart 1/14.5; Air 1.5/21.8; HY+Ophthalmart 4.5 bar   65 psi
Housing burst pressure	6 bar   87 psi and higher
Max. temperature	50°C
Bubble point	Min value with water (Ophthalmart): ≥ 3.2 bar   ≥ 46 psi (0.2 µm) Min value with IPA (Air; HY): ≥ 1.4 bar   ≥ 20.3 psi (0.2 µm), ≥ 0.9 bar   ≥ 13.1 psi (0.45 µm) Acticosart: BP not determined
Flow rate	Ophthalmart: 60 ml/min for water at 1 bar   14.5 psi Air: 1.1 l/min (0.2 µm) for air at 0.1 bar   1.5 psi HY: 1.4 l/min (0.2 µm), 2.3 l/min (0.45 µm) for air at 0.1 bar   1.5 psi Acticosart: 2.3 l/min (0.45 µm) for air at 0.1 bar   1.5 psi
Cytotoxicity	Detectably no inhibition with MRC-5 and L-929 (Ophthalmart)

## ■ Sartolab® P20 and Sartolab® P20 Plus

Sterile Filtration of Sample Volumes Up to 5 Liter



### Using Sartolab® P20 or Sartolab® P20 Plus in Available Systems

#### Systems with luer lock connectors

Sartolab® units with a luer lock inlet fit directly onto the corresponding connectors of tubing from peristaltic pumps. This configuration can also be used with luer lock syringes.



### □ Application and Product Description

#### Application

The Sartolab® P20 is a ready-to-use pressure filtration unit for sterile filtration of media and aqueous solutions in batches ranging from 100 ml to 5 l. For media that contain sera and difficult-to-filter solutions, a Sartolab® P20 Plus unit with an incorporated prefilter is also available.

#### Security


The effectiveness of every batch of 0.2 µm cellulose acetate membranes for sterile filtration is confirmed by bacteria challenge tests (HIMA) using *Brevundimonas diminuta*. Only biosafe material is used in the filtration units; they have been proven non-toxic by passing the USP plastics test for toxicity. Tests with MRC-5 human lung cells on cellulose acetate membranes and glass fiber prefilters showed no cytotoxic effects. Finished units are also tested for their sterile filtration capability and for housing and membrane integrity.

#### Systems with Tubing

If the liquid to be filtered will be pumped through open tubing from a pressure tank or a peristaltic pump, then a Sartolab® unit with a hose nipple inlet is required. The stepwise increase in diameter of the nipple, from 6 mm to 12 mm, makes it suitable for a wide range of tubing.

#### Fast Filtration


The combination of a large filtration area (20 cm<sup>2</sup>) and the optimal design of the filter support guarantees high flow rates with high total throughputs. Automatic venting of any trapped air through the PTFE membrane-protected vent ports ensures that the entire filter surface is used for effective filtration.

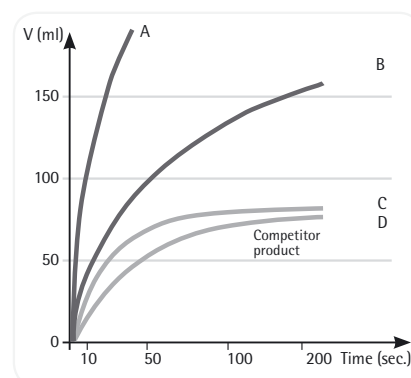

Choice of inlet connector,  
hose nipple or luer lock

With or without  
binder-free glass fiber  
prefilter

Polycarbonate housing  
and surfactant-free  
cellulose acetate  
membrane filter

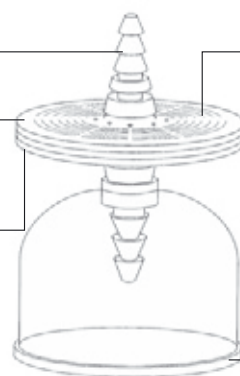

20 cm<sup>2</sup> filter area,  
optimal flow path and  
automatic venting

Outlet side with protective  
filling bell

## Technical Specifications

### Technical Specifications for Sartolab® P20 and Sartolab® P20 Plus Units with SFCA Membrane

Properties	Description	
	Sartolab® P20 18052 18053	Sartolab® P20 Plus 18056 18058
Filter material	SFCA, type 12587, 0.2 µm pore size and PTFE	SFCA, type 12587 0.2 µm pore size plus GF, 100 % free of binding agents and PTFE
Housing material	Polycarbonate	Polycarbonate
Color code	Transparent	Transparent
Filter diameter	64 mm	64 mm
Connector inlet	Female luer lock or stepped hose nipple with 6–12 mm outer diameter	Female luer lock or stepped hose nipple with 6–12 mm outer diameter
Connector outlet	Hose nipple	Hose nipple
Filling bell	Yes	Yes
Filtration area	20 cm <sup>2</sup>	20 cm <sup>2</sup>
Hold-up volume before bubble point	1 ml	Approx 1.5 ml
Housing burst pressure	> 5 bar   72.5 psi	> 5 bar   72.5 psi
Bubble point	≥ 3.2 bar   46.4 psi	≥ 3.2 bar   46.4 psi
Max. recommended inlet pressure	3 bar   43.5 psi	3 bar   43.5 psi
Flow rate for water	≥ 250 ml/min at Δp = 1 bar   14.5 psi	≥ 250 ml/min at Δp = 1 bar   14.5 psi
Filtration range	100 ml–max. 5 l	100 ml–max. 10 l
pH-range	4–8	4–8
Non-specific protein adsorption	No loss of protein detectable (filtration of γ globulin, method acc. to Bradford)	< 80 µg/cm <sup>2</sup> (filtration of γ globulin, method acc. to Bradford)
Sterilization	EO sterilization	EO sterilization
Biosafety	Class VI Plastics Test	Class VI Plastics Test
Operating instructions	Directions for use included in each box	Directions for use included in each box

### Technical Specifications for Sartolab® P20 and Sartolab® P20 Plus Units with PES Membrane

	Sartolab® P20 Plus 18068	Sartolab® P20 18075
Filter material	PES, type 15407 MI 0.2 µm pore size plus GF, 100 % free of binding agents and PTFE	PES, type 15407 MI 0.2 µm pore size
Housing material	Polycarbonate	Polycarbonate
Color code	Transparent	Transparent
Filter diameter	64 mm	64 mm
Connector inlet	Female luer lock or stepped hose nipple with 6–12 mm outer diameter	Female luer lock or stepped hose nipple with 6–12 mm outer diameter
Connector outlet	Hose nipple	Hose nipple
Filling bell	Yes	No
Filtration area	20 cm <sup>2</sup>	20 cm <sup>2</sup>
Hold-up volume before bubble point	Approx 1.5 ml	1 ml
Housing burst pressure	> 5 bar   72.5 psi	> 5 bar   72.5 psi

Properties	Description	
	Sartolab® P20 Plus 18068	Sartolab® P20 18075
Bubble point	≥ 3.2 bar   46.4 psi	≥ 3.2 bar   46.4 psi
Max. recommended inlet pressure	3 bar   43.5 psi	3 bar   43.5 psi
Flow rate for water	400 ml/min at Δp = 1 bar   14.5 psi	400 ml/min at Δp = 1 bar   14.5 psi
Filtration range	100 ml–max. 5 l	100 ml–max. 10 l
pH range	1–8	1–8
Non-specific protein adsorption	< 80 µg/cm <sup>2</sup> (filtration of γ globulin, method acc. to Bradford)	No loss of protein detectable (filtration of γ globulin, method acc. to Bradford)
Sterilization	EO sterilization	EO sterilization
Biosafety	Class VI Plastics Test	Class VI Plastics Test
Operating instructions	Directions for use included in each box	Directions for use included in each box

## Ordering Information

Type	Membrane	Housing	Pore Size	Inlet	Outlet	Sterile	Qty   Pk	Order No.
<b>Sartolab® P20</b>								
	SFCA	PC	0.2 µm	Hose Nipple	Hose Nipple	Yes	10	18052-----D
	SFCA	PC	0.2 µm	Luer Lock	Hose Nipple	Yes	10	18053-----D
	PES	PC	0.2 µm	Luer Lock	Hose Nipple	Yes	10	18075-----D
<b>Sartolab® P20 Plus</b>								
	SFCA + GF	PC	0.2 µm	Hose Nipple	Hose Nipple	Yes	10	18056-----D
	SFCA + GF	PC	0.2 µm	Luer Lock	Hose Nipple	Yes	10	18058-----D
	PES + GF	PC	0.2 µm	Luer Lock	Hose Nipple	Yes	10	18068-----D

SFCA - Cellulose Acetate, PES - Polyethersulfone, GF - Glass Fiber Prefilter, PC - Polycarbonate

## ■ Sartolab® 150v

Disposable PES Vacuum Filtration Unit for Volumes Up to 15 Liter



**Description**  
Sartolab® 150v is a disposable, sterile, ready-to-use membrane filter capsule for highest convenience. Sartolab® 150v capsules are made of a unique hydrophilic Polyethersulfone (PES) membrane providing outstanding total throughput, flow rate performance, low extractables and broadest chemical compatibility.

**Applications**  
Typical applications include sterilizing grade filtration of:

- Biological Fluids
- Purified water
- Media
- Buffers

**Compatibility**  
The PES membrane is compatible with a pH range from pH 1 to pH 14 making Sartolab® 150v ideal for filtration of solutions with high | low pH.

**Microbiological Retention**  
Sartolab® 150v capsules 0.2 µm rated are fully validated as sterilizing grade filters according to HIMA and ASTM F-838-05 guidelines.

**Quality Control**  
Each individual element is tested for integrity by B.-P. and Diffusion-Test before final assembly. Sartolab® 150v capsules are designed, developed and manufactured in accordance with a ISO 9001:2000 certified Quality Management System.

**Performance**  
The unique pleated filter construction combined with the highly asymmetric pore structure of the PES membrane, offers excellent flow rates and superior total throughput performance, especially in comparison to conventional disc filter systems.


### □ Materials and Data

Filter	Polyethersulfone (PES)
Support Fleece	PP-fleece
Housing	Polypropylene
Vent	PTFE
Pore Size	0.2 µm
Filtration Area	150 cm²
Flow Rate	1.5 L/min
Total Throughput	15 L


### □ Ordering Information

Units   Pk	Order No.
3	18080-M


#### Easy to Use

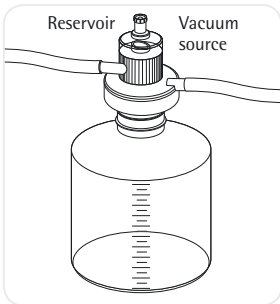


Place Sartolab® 150v on the bottle neck. A vacuum is automatically created when used with a water jet pump or other vacuum source.

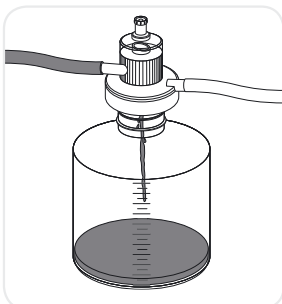

Create a vacuum to start filtration.

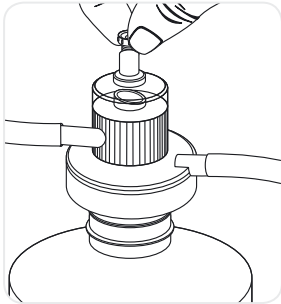

To fill multiple receiving bottles, break the vacuum with the valve on top of the Sartolab® 150v housing.

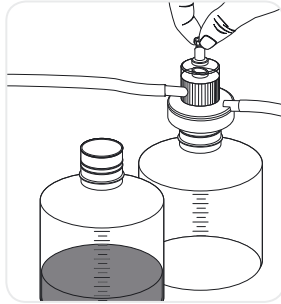

To continue filtration with the next receiving bottle, close the valve.

## ■ Sartolab® RF|BT

### Disposable Filter Systems and Bottle Top Filters



Sartolab® Disposable Sterile Filter Systems and Bottle Top Filters are designed for the vacuum filtration of tissue culture media and components, biological fluids, and other aqueous solutions.

The Sartolab® disposable 150 mL, 250 mL, 500 mL and 1000 mL bottles are designed as storage containers for sterile media, buffers, or other aqueous solutions.



These products are for laboratory use only. Not for human parenteral applications.

#### Materials

The filter funnels, dust covers and receiver bottles are manufactured from virgin, heavy metal-free polystyrene. The tubing adapters, filter adapters, and the plug seal caps are made from heavy metal-free polyethylene. Sartolab® filter systems are available with polyethersulfone membranes. All units are sterilized by gamma irradiation.

#### Performance

The filter units contain membranes integrally sealed to a support grid designed to maximize flow and reduce foaming and protein denaturation.

The membrane is compatible with most aqueous solutions and tested for use in cell culture applications.

#### Filter Systems

The filter adapter utilizes a gasket design to ensure a vacuum-tight seal on the receiver|storage bottle. Each filter unit also contains a convenient tubing adapter that will fit most vacuum hoses.

The bottles are single-use containers.

They cannot withstand autoclaving or use at temperatures greater than 70°C. The suitability of the bottles for storage of solutions below 0°C depends both on the solution and the storage conditions. Many aqueous solutions, including culture media, have been successfully frozen and stored at temperatures down to -20°C. However, a trial run under actual conditions is strongly recommended to test the suitability of the bottles for frozen storage.

Pore Size	Membrane Material	Characteristics
0.22 µm	Polyethersulfone	Very low protein binding and low extractables, fast flow rate
0.1 µm	Polyethersulfone	Very low protein binding and low extractables, fast flow rate

**Bottle Top Filters**

The filter adapter is available in 45 mm thread finish, and is designed to ensure a vacuum tight seal on customer supplied bottles with the appropriate thread finish. Each filter unit also contains a tubing adapter that will fit most vacuum hoses.

**Chemical Compatibility**

The mechanical strength, color, appearance, and dimensional stability of filter systems, bottle top filters, and plastic bottles are affected to varying degrees by the chemicals with which they come in contact. Specific operating conditions, especially temperature, will also affect their chemical resistance. A table is provided to serve as a general guideline for the chemical resistance of Sartolab® disposable sterile filters and bottles.

 Chemical Resistance of Sartolab® Filters


Chemical Class	Membrane (PES)	Housing (PS)
Weak Acids	3	1
Strong Acids	3	2
Alcohols	1	2
Aldehydes	3	3
Aliphatic Amines	1	3
Aromatic Amines	3	3
Bases	3	1
Esters	3	3
Hydrocarbons	3	3
Ketones	3	3

Key: 1, recommended;  
 2, may be suitable for some applications; a trial run is recommended;  
 3, not recommended. PS, polystyrene; PES, polyethersulfone.

## Ordering Information

Volume	Membrane	Filter Area	Qty   Pk	Order No.
<b>Sartolab® RF Filtration System Including Collection Bottle</b>				
150 ml	0.22 µm PES	18 cm <sup>2</sup>	12	180C1-----E
250 ml	0.22 µm PES	24 cm <sup>2</sup>	12	180C7-----E
500 ml	0.22 µm PES	39 cm <sup>2</sup>	12	180C2-----E
1,000 ml	0.22 µm PES	62 cm <sup>2</sup>	12	180C3-----E
1,000 ml	0.1 µm PES	62 cm <sup>2</sup>	12	180C8-----E
<b>Sartolab® BT Bottle Top Filters Without Integrated Collection Bottle</b>				
150 ml	0.22 µm PES	18 cm <sup>2</sup>	48	180C4-----K
500 ml	0.22 µm PES	39 cm <sup>2</sup>	12	180C5-----E
1,000 ml	0.22 µm PES	62 cm <sup>2</sup>	12	180C6-----E

## Water Throughput

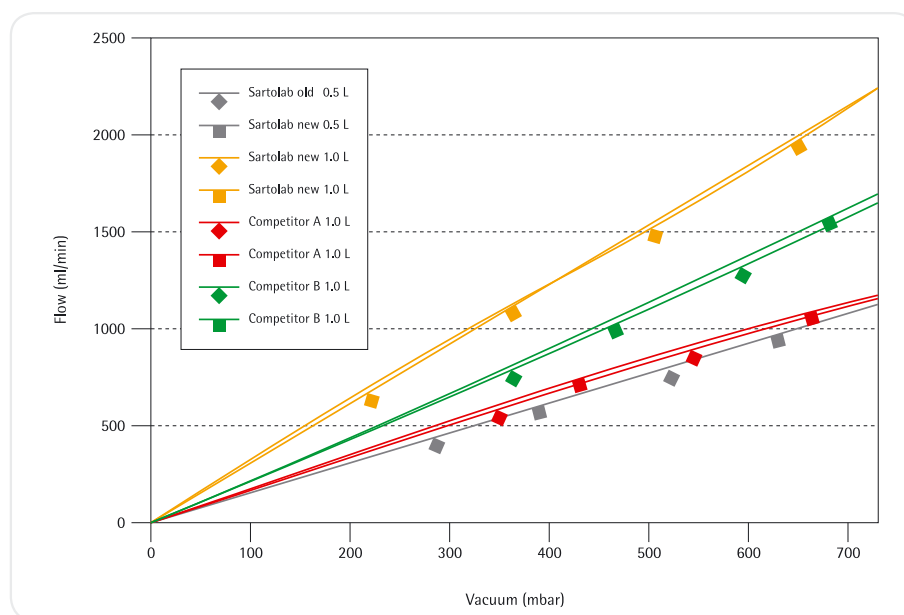

## Chemical Compatibility

	PES membrane	SFCA membrane	PTFE membrane	RC membrane	Nylon XXX	GF depth filter	Housing MBS	Housing PP	Minisart® HighFlow	Minisart® NML   Ophthalsart	Minisart® plus	Minisart® GF	Minisart® HY	Minisart® RC	Minisart® NY	Minisart® SRP
Filter	PES	SFCA	PTFE	RC	PA	GF			PES	SFCA	SFCA + GF	GF	PTFE	RC	PA	PTFE
Housing							MBS	PP	MBS	MBS	MBS	MBS	MBS	PP	PP	PP
Sterilization																
Ethylene oxide	++	++	++	++	++	++	++	++	++	++	++	++	++	++	++	++
Gamma irradiation	++	++	- <sup>1</sup>	++	-	++	++	-	++	++	++	++	-	-	-	-
Autoclaving 121 °C, 30 min	++	++	++	++	++	++	-	++	-	-	-	-	-	++	++	++
Solvents																
Acetone	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++
Acetonitrile	-	-	++	++	n.a.	++	-	++	-	-	-	-	-	++	n.a.	++
Gasoline	+	+	++	++	++	++	++	++	+	+	+	++	++	++	++	++
Benzene	-	+	++	++	++	++	-	++	-	-	-	-	-	++	++	++
Benzyl alcohol	-	-	++	++	++	++	-	+	-	-	-	-	-	+	+	+
n-Butyl acetate	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++
n-Butanol	+	+	++	++	++	++	++	++	+	+	+	++	++	++	++	++
Cellosolve	++	-	++	++	++	++	-	-	-	-	-	-	-	-	-	-
Chloroform	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++
Cyclohexane	-	+	++	++	++	++	+	+	-	+	+	+	+	+	+	+
Cyclohexanone	-	+	++	++	++	++	-	+	-	-	-	-	-	+	+	+
Diethylacetamide	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++
Diethyl ether	-	+	++	++	++	++	-	++	-	-	-	-	-	++	++	++
Dimethyl formamide	-	-	++	+	+	++	-	+	-	-	-	-	-	+	+	+
Dimethylsulfoxide	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++
Dioxane	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++
Ethanol, 98%	++	+	++	++	++	++	-	+	-	-	-	-	-	+	+	+
Ethyl acetate	-	-	++	++	++	++	-	+	-	-	-	-	-	+	+	+
Ethylene glycol	++	+	++	++	++	++	++	++	++	+	+	++	++	++	++	++
Formamide	++	-	++	+	++	++	++	++	++	-	-	++	++	+	++	++
Glycerin	++	+	++	++	++	++	+	+	+	+	+	+	+	+	+	+
n-Heptane	++	+	++	++	++	++	+	++	+	+	+	+	+	++	++	++
n-Hexane	++	+	++	++	++	++	+	+	+	+	+	+	+	+	+	+
Isobutanol	++	+	++	++	++	++	++	-	++	+	+	++	++	-	-	-
Isopropanol	++	+	++	++	++	++	-	++	-	-	+	-	-	++	++	++
Isopropyl acetate	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++
Methanol, 98%	+	-	++	++	++	++	++	+	+	-	+	++	++	+	+	+
Methyl acetate	-	-	++	++	++	++	-	+	-	-	-	-	-	+	+	+
Methylene chloride	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++
Methyl ethyl ketone	-	-	++	++	++	++	-	+	-	-	-	-	-	+	+	+
Methyl isobutyl ketone	-	-	++	++	++	++	-	+	-	-	-	-	-	+	+	+
Monochlorobenzene	-	-	++	++	++	++	-	+	-	-	-	-	-	+	+	+
Nitrobenzene	-	-	++	++	+	++	-	+	-	-	+	-	-	+	+	+
n-Pentane	++	+	++	++	++	++	+	++	+	+	+	+	+	++	++	++
Perchloroethylene	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++


	PES membrane	SFCA membrane	PTFE membrane	RC membrane	Nylon XXX	GF depth filter	Housing MBS	Housing PP	Minisart® HighFlow	Minisart® NML   Ophthalsart	Minisart® plus	Minisart® GF	Minisart® HY	Minisart® RC	Minisart® NY	Minisart® SRP
Filter	PES	SFCA	PTFE	RC	PA	GF			PES	SFCA	SFCA + GF	GF	PTFE	RC	PA	PTFE
Housing							MBS	PP	MBS	MBS	MBS	MBS	MBS	PP	PP	PP
<b>Sterilization</b>																
Pyridine	-	-	++	++	++	++	-	++	-	-	-	-	-	++	++	++
Carbon tetrachloride	-	-	++	++	++	++	-	-	-	-	-	-	-	-	-	-
Tetrahydrofuran	-	-	++	++	++	++	-	++	-	-	++	-	-	++	++	++
Toluene	-	+	++	++	++	++	-	++	-	-	++	-	-	++	++	++
Trichloroethane	-	-	++	++	++	++	-	n.a.	-	-	-	-	-	n.a.	n.a.	n.a.
Trichloroethylene	-	+	++	++	++	++	-	++	-	-	++	-	-	++	++	++
Xylene	-	+	++	++	++	++	-	+	-	-	-	-	-	+	+	+
<b>Acids</b>																
Acetic acid, 25%	+	+	++	++	-	++	-	+	-	-	-	-	-	+	-	+
Acetic acid, 80%	n.a.	-	++	++	-	++	-	+	-	-	-	-	-	+	-	+
Hydrofluoric acid, 25%	+	-	++	+	-	++	+	+	+	-	-	+	+	+	-	+
Hydrofluoric acid, 50%	+	-	++	+	-	++	+	+	+	-	-	+	+	+	-	+
Perchloric acid, 25%	-	-	++	-	-	++	n.a.	+	-	-	-	n.a.	n.a.	-	-	+
Phosphoric acid, 1%	++	+	++	-	-	++	+	+	+	+	+	+	+	-	-	+
Phosphoric acid, 86%	++	+	++	-	-	++	+	+	+	+	+	+	+	-	-	+
Nitric acid, 30%	+	-	++	-	-	++	+	+	+	-	-	+	+	-	-	+
Nitric acid, conc.	-	-	++	-	-	++	-	-	-	-	-	-	-	-	-	-
Hydrochloric acid, 15%	++	+	++	-	-	++	+	+	+	+	+	+	+	-	-	+
Hydrochloric acid, 20%	++	-	++	-	-	++	+	+	+	-	-	+	+	-	-	+
Sulfuric acid, 25%	+	-	++	+	-	++	+	++	+	-	-	+	+	+	-	++
Sulfuric acid, 98%	-	-	++	-	-	++	-	+	-	-	-	-	-	-	-	+
Trichloroacetic acid, 25%	n.a.	-	++	++	-	++	-	+	-	-	-	-	-	+	-	+
<b>Bases</b>																
Ammonia, 1N	++	-	++	+	++	++	-	++	-	-	-	-	-	+	++	++
Ammonium hydroxide, 25%	+	+	++	+	++	+	-	+	-	-	-	-	-	+	+	+
Potassium hydroxide, 32%	++	-	++	-	+	+	-	++	-	-	-	-	-	-	+	++
Sodium hydroxide, 32%	++	-	++	-	+	+	-	+	-	-	-	-	-	-	+	+
Sodium hydroxide, 1N	++	-	++	+	++	+	-	++	-	-	-	-	-	+	++	++
<b>Aqueous solutions</b>																
Formalin, 30%	+	++	++	+	++	++	+	+	+	+	+	+	+	+	+	+
Sodium hypochlorite, 5%	++	-	++	-	-	++	+	+	+	-	-	+	+	-	-	+
Hydrogen peroxide, 35%	++	-	++	-	-	++	+	++	+	-	-	+	+	-	-	++
<b>pH range</b>																
pH 1-14	-	-	++	-	-	++	-	++								
pH 1-13	++	-	++	-	-	++	-	++								
pH 3-14	+	-	++	+	++	++	-	++								
pH 3-12	++	-	++	++	++	++	+	++								
pH 4-8	++	++	++	++	++	++	++	++								

**Legend**

Compatible	++
Limited compatibility	+
Not compatible	-
not analyzed	n.a.

<sup>1)</sup> gamma irradiation feasible for Minisart® Air

Contact time: 24 hours at 20 °C. Chemical compatibilities can be influenced by various factors. Therefore, we recommend that you confirm compatibility with the liquid you want to filter by performing a trial filtration run before you start your actual filtration.

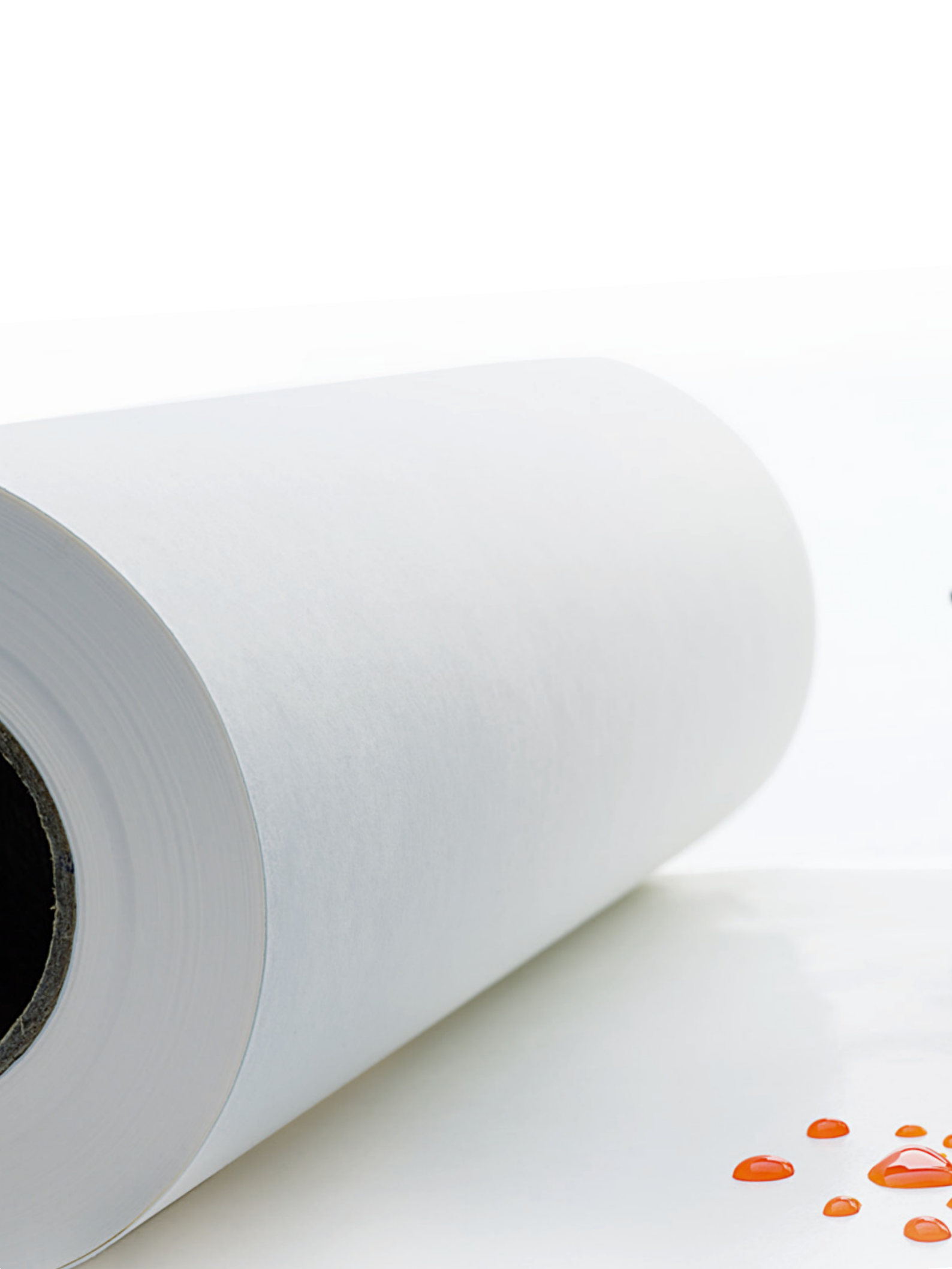

## ■ Basic Filtration

### Table of Contents

76	Filter Papers
80	Glass & Quartz Microfiber Filters
84	Extraction Thimbles
85	Membrane Filters
91	Blotting   Chromatography Papers & Blotting Membranes
92	Filtration Equipment
102	Combisart® manifolds










## Ash-free Filter Papers

For Quantitative and Gravimetric Analyses

These filter papers are used for quantitative and gravimetric analyses as well as for pressure or vacuum filtration. They are made out of 100% cotton linters with an  $\alpha$ -cellulose content of > 98% and are acid-washed to make the papers ashless and achieve high purity.

### Technical Specifications

Grade	Weight (g/m <sup>2</sup> )	Thickness (mm)	Particle Retention (μm)	Filtration (s)	Precipitates	Properties
 388	84	0.21	12-15	10	Coarse crystalline	Wide-pore, loose structure, fast filtering
 389	84	0.19	8-12	20	Medium-fine crystalline	Medium- to wide-pore, medium fast filtering
 389 F	84	0.19	8-12	20	Medium-fine crystalline	Medium- to wide-pore, medium fast filtering
 392	84	0.17	5-8	50	Fine crystalline	Medium dense, medium fast filtering
 390	84	0.16	3-5	100	Fine crystalline	Narrow-pore, dense, slow filtering
 391	84	0.15	2-3	180	Very fine crystalline	Fine-pore, dense, very slow filtering
 393	100	0.18	1-2	300	Very fine crystalline	Very fine-pore, very dense, very slow filtering

### Ordering Information



#### Filter Discs, 100 Pieces

Ø in mm	Grade 388	Grade 389	Grade 389 F	Grade 390	Grade 391	Grade 392	Grade 393
55	FT-3-101-055	FT-3-102-055	FT-3-112-055	FT-3-103-055	FT-3-104-055	FT-3-105-055	FT-3-127-055
70	FT-3-101-070	FT-3-102-070	FT-3-112-070	FT-3-103-070	FT-3-104-070	FT-3-105-070	FT-3-127-070
90	FT-3-101-090	FT-3-102-090	FT-3-112-090	FT-3-103-090	FT-3-104-090	FT-3-105-090	FT-3-127-090
110	FT-3-101-110	FT-3-102-110	FT-3-112-110	FT-3-103-110	FT-3-104-110	FT-3-105-110	FT-3-127-110
125	FT-3-101-125	FT-3-102-125	FT-3-112-125	FT-3-103-125	FT-3-104-125	FT-3-105-125	FT-3-127-125
150	FT-3-101-150	FT-3-102-150	FT-3-112-150	FT-3-103-150	FT-3-104-150	FT-3-105-150	FT-3-127-150
185	FT-3-101-185	FT-3-102-185	FT-3-112-185	FT-3-103-185	FT-3-104-185	FT-3-105-185	FT-3-127-185
240	FT-3-101-240	FT-3-102-240	FT-3-112-240	FT-3-103-240	FT-3-104-240	FT-3-105-240	FT-3-127-240



#### Folded Filters, 100 Pieces

Ø in mm	Grade 388	Grade 389	Grade 389 F	Grade 390	Grade 391	Grade 392
110	FT-4-101-110	FT-4-102-110	FT-4-112-110	FT-4-103-110	FT-4-104-110	FT-4-105-110
125	FT-4-101-125	FT-4-102-125	FT-4-112-125	FT-4-103-125	FT-4-104-125	FT-4-105-125
150	FT-4-101-150	FT-4-102-150	FT-4-112-150	FT-4-103-150	FT-4-104-150	FT-4-105-150
185	FT-4-101-185	FT-4-102-185	FT-4-112-185	FT-4-103-185	FT-4-104-185	FT-4-105-185
240	FT-4-101-240	FT-4-102-240	FT-4-112-240	FT-4-103-240	FT-4-104-240	FT-4-105-240



#### Sheets in 580×580 mm, 100 Pieces

Grade 388	Grade 389	Grade 390	Grade 391	Grade 392	Grade 393
FT-2-101-580580	FT-2-102-580580	FT-2-103-580580	FT-2-104-580580	FT-2-105-580580	FT-2-127-580580

Other dimensions are available on request

## Wet-strengthened Filter Papers

For Qualitative Analyses

These qualitative filter papers are essentially used for analytical purposes and routine analyses, whenever no gravimetric analyses are required. They are wet-strengthened and can be used for pressure and vacuum filtration. They are made of refined pulp and linters with an >95%  $\alpha$ -cellulose content, are very pure with an ash content  $\leq 0.1\%$ .

### Technical Specifications

Grade	Weight (g/m <sup>2</sup> )	Thickness (mm)	Particle Retention (μm)	Filtration (s)	Precipitates	Properties
1288	84	0.21	12-15	10	Coarse crystalline	Wide-pore, loose structure, fast filtering
1289	84	0.21	8-12	20	Medium-fine crystalline	Medium- to wide-pore, medium fast filtering
1292	84	0.17	5-8	50	Fine crystalline	Medium dense, medium fast filtering
1290	84	0.21	3-5	100	Fine crystalline	Narrow-pore, dense, slow filtering
1291	84	0.16	2-3	180	Very fine crystalline	Fine-pore, dense, very slow filtering
293	80	0.15	1-2	300	Very fine crystalline	Very fine-pore, very dense, very slow filtering

### Ordering Information



#### Filter Discs, 100 Pieces

Ø in mm	Grade 1288	Grade 1289	Grade 1290	Grade 1291	Grade 1292	Grade 293
55	FT-3-206-055	FT-3-207-055	FT-3-208-055	FT-3-209-055	FT-3-210-055	FT-3-211-055
70	FT-3-206-070	FT-3-207-070	FT-3-208-070	FT-3-209-070	FT-3-210-070	FT-3-211-070
90	FT-3-206-090	FT-3-207-090	FT-3-208-090	FT-3-209-090	FT-3-210-090	FT-3-211-090
110	FT-3-206-110	FT-3-207-110	FT-3-208-110	FT-3-209-110	FT-3-210-110	FT-3-211-110
125	FT-3-206-125	FT-3-207-125	FT-3-208-125	FT-3-209-125	FT-3-210-125	FT-3-211-125
150	FT-3-206-150	FT-3-207-150	FT-3-208-150	FT-3-209-150	FT-3-210-150	FT-3-211-150
185	FT-3-206-185	FT-3-207-185	FT-3-208-185	FT-3-209-185	FT-3-210-185	FT-3-211-185
240	FT-3-206-240	FT-3-207-240	FT-3-208-240	FT-3-209-240	FT-3-210-240	FT-3-211-240



#### Folded Filters, 100 Pieces

Ø in mm	Grade 1288	Grade 1289	Grade 1290	Grade 1291	Grade 1292	Grade 293
110	FT-4-206-110	FT-4-207-110	FT-4-208-110	FT-4-209-110	FT-4-210-110	FT-4-211-110
125	FT-4-206-125	FT-4-207-125	FT-4-208-125	FT-4-209-125	FT-4-210-125	FT-4-211-125
150	FT-4-206-150	FT-4-207-150	FT-4-208-150	FT-4-209-150	FT-4-210-150	FT-4-211-150
185	FT-4-206-185	FT-4-207-185	FT-4-208-185	FT-4-209-185	FT-4-210-185	FT-4-211-185
240	FT-4-206-240	FT-4-207-240	FT-4-208-240	FT-4-209-240	FT-4-210-240	FT-4-211-240



#### Sheets in 580 × 580 mm, 100 Pieces

Grade 1288	Grade 1289	Grade 1290	Grade 1291	Grade 1292	Grade 293
FT-2-206-580580	FT-2-207-580580	FT-2-208-580580	FT-2-209-580580	FT-2-210-580580	FT-2-211-580580

Other dimensions are available on request

## High-Purity Filter Papers

### For Qualitative Analyses

These paper grades are used for analytical purposes that require a low ash content. Grades 292 and 292a are especially suitable for soil analyses because they are low in nitrogen. For phosphate or sodium determination, we recommend grades 131 and 132.

### Technical Specifications

Grade	Weight (g/m <sup>2</sup> )	Thickness (mm)	Particle Retention (µm)	Filtration (s)	Material
292	87	0.18	5-8	45	Cotton linters, low-nitrogen and nitrates, ash content ≤ 0.06% according to DIN 54370
292a	97	0.19	4-7	60	Cotton linters, low-nitrogen and nitrates, ash content ≤ 0.06% according to DIN 54370
131	80	0.16	3-5	100	Cotton linters and refined pulp, low-phosphate and low sodium, ash content ≤ 0.02% according to DIN 54370
132	80	0.17	5-7	55	Cotton linters and refined pulp, low-phosphate and low sodium, ash content ≤ 0.02% according to DIN 54370

### Ordering Information



#### Filter Discs, 100 Pieces

Ø in mm	Grade 131	Grade 132	Grade 292	Grade 292a
55	FT-3-351-055	FT-3-329-055	FT-3-205-055	FT-3-215-055
70	FT-3-351-070	FT-3-329-070	FT-3-205-070	FT-3-215-070
90	FT-3-351-090	FT-3-329-090	FT-3-205-090	FT-3-215-090
110	FT-3-351-110	FT-3-329-110	FT-3-205-110	FT-3-215-110
125	FT-3-351-125	FT-3-329-125	FT-3-205-125	FT-3-215-125
150	FT-3-351-150	FT-3-329-150	FT-3-205-150	FT-3-215-150
185	FT-3-351-185	FT-3-329-185	FT-3-205-185	FT-3-215-185
240	FT-3-351-240	FT-3-329-240	FT-3-205-240	FT-3-215-240



#### Folded Filters, 100 Pieces

Ø in mm	Grade 131	Grade 132	Grade 292	Grade 292a
110	FT-4-351-110	FT-4-329-110	FT-4-205-110	FT-4-215-110
125	FT-4-351-125	FT-4-329-125	FT-4-205-125	FT-4-215-125
150	FT-4-351-150	FT-4-329-150	FT-4-205-150	FT-4-215-150
185	FT-4-351-185	FT-4-329-185	FT-4-205-185	FT-4-215-185
240	FT-4-351-240	FT-4-329-240	FT-4-205-240	FT-4-215-240



#### Sheets in 580 × 580 mm, 100 Pieces

Grade 131	Grade 132	Grade 292	Grade 292a
FT-2-351-580580	FT-2-329-580580	FT-2-205-580580	FT-2-215-580580

Other dimensions are available on request

## Filter Papers

### For Qualitative-Technical Analyses

These filter papers are used for routine analyses like clarification, determination of substances, but also as discs with a center hole for technical applications. Grades with a wet burst resistance > 30 kPa are referred to as wet-strengthened and are therefore suitable for pressure or vacuum filtration. They are made of refined pulp and linters with an > 95%  $\alpha$ -cellulose content, are very pure with an ash content between <0.1–0.15%. Below you will find an overview of the most commonly used grades.

#### Technical Specifications

Grade	Surface	Weight (g/m <sup>2</sup> )	Thickness (mm)	Filtration (s)	Wet Burst Resistance (kPa)	Properties
3 hw	Smooth	65	0.14	20	> 40	Medium fast filtering, filter paper for routine work in the lab
4 b	Smooth	75	0.15	22	> 15	Medium fast filtering, filtration of coarse precipitates, wick paper for seed testing
603/N	Crêped	75	0.25	8	> 50	Fast filtering, filtration of sugar solutions
6	Smooth	80	0.17	15	> 30	Fast filtering, degassing beer before analysis, clarification of spirits
100/N	Smooth	85	0.18	30	> 80	Medium fast filtering, ash content <0.1%, low potassium and sodium content, determination of the sugar content
5 H/N	Crêped	85	0.28	3	> 40	Very fast filtering, wide-pore, filtration of essential oils
3 S/h	Smooth	200	0.36	55	> 15	Medium fast to slow filtering, narrow-pore, re-wet test for diapers

#### Ordering Information

##### Filter Discs

Ø in mm	Grade 3 hw (100 Pieces)	Grade 4 b (100 Pieces)	Grade 603/N (100 Pieces)	Grade 6 (100 Pieces)	Grade 100/N (100 Pieces)	Grade 5 H/N (100 Pieces)	Grade 3 S/h (50 Pieces)
55	FT-3-303-055	FT-3-309-055		FT-3-312-055	FT-3-328-055		FT-3-307-055
70	FT-3-303-070	FT-3-309-070		FT-3-312-070	FT-3-328-070		
90	FT-3-303-090	FT-3-309-090	FT-3-335-090	FT-3-312-090	FT-3-328-090	FT-3-423-090	FT-3-307-090
110	FT-3-303-110	FT-3-309-110	FT-3-335-110	FT-3-312-110	FT-3-328-110		FT-3-307-110
125	FT-3-303-125	FT-3-309-125	FT-3-335-125	FT-3-312-125	FT-3-328-125	FT-3-423-125	FT-3-307-125
150	FT-3-303-150	FT-3-309-150	FT-3-335-150	FT-3-312-150	FT-3-328-150	FT-3-423-150	FT-3-307-150
185	FT-3-303-185	FT-3-309-185	FT-3-335-185	FT-3-312-185	FT-3-328-185	FT-3-423-185	FT-3-307-185
240	FT-3-303-240	FT-3-309-240	FT-3-335-240	FT-3-312-240	FT-3-328-240	FT-3-423-240	FT-3-307-240

##### Folded Filters, 100 Pieces

Ø in mm	Grade 3 hw	Grade 4 b	Grade 603/N	Grade 6	Grade 100/N	Grade 5 H/N
125	FT-4-303-125	FT-4-309-125	FT-4-335-125	FT-4-312-125		FT-4-423-125
150	FT-4-303-150	FT-4-309-150	FT-4-335-150	FT-4-312-150	FT-4-328-150	FT-4-423-150
185	FT-4-303-185	FT-4-309-185	FT-4-335-185	FT-4-312-185		FT-4-423-185
240	FT-4-303-240	FT-4-309-240	FT-4-335-240	FT-4-312-240	FT-4-328-240	FT-4-423-240
270	FT-4-303-270	FT-4-309-270	FT-4-335-270	FT-4-312-270	FT-4-328-270	FT-4-423-270
320	FT-4-303-320	FT-4-309-320	FT-4-335-320	FT-4-312-320	FT-4-328-320	FT-4-423-320

##### Sheets in 580 × 580 mm, 100 Pieces

Grade 3 hw	Grade 4 b	Grade 603/N	Grade 6	Grade 100/N	Grade 5 H/N
FT-2-303-580580	FT-2-309-580580	FT-2-335-580580	FT-2-312-580580	FT-2-328-580580	FT-2-423-580580

Other dimensions are available on request

## Glass Microfiber Filters

### Without Binder

Binder-free glass microfiber filters are recommended for analytical and gravimetric analyses and also as prefilters. These filters combine fast flow rates with high load capacity and the retention of very fine particles; they are biologically inert, are resistant to most chemicals and withstand temperatures up to 500°C (grade 550-HA up to 550°C).

### Technical Specifications

Grade	Weight (g/m <sup>2</sup> )	Thickness (mm)	Air Resistance (mbar)	Particle Retention (µm)	Filtration (s)	Main Applications
MGA	52	0.25	36.0	1.6	40	Clarification of buffer and reagent solutions
MG 160	75	0.37		1.2		Air monitoring
MGB	143	0.70	95.0	1.0	100	Membrane prefilter
MGC	52	0.25	51.0	1.2	66	Analysis of suspended solids in wastewater according to DIN EN 872
MGD	120	0.53	18.5	2.7	18	Membrane prefilter
MGF	75	0.40	145.0	0.7	185	Clarification of protein solutions
MGG	64	0.28	32.0	1.5	38	Water filtration
MG 550-HA 65		0.30		1.5		Analysis of suspended solids in wastewater according to 2540D
13440	92	0.44		0.7		Membrane prefilter

### Ordering Information



#### Filter Discs

Ø in mm	MGA (100 Pieces)	MG 160 (50 Pieces)	MGB (50 Pieces)	MGC (100 Pieces)	MGD (50 Pieces)
21			FT-3-1102-021		
25	FT-3-1101-025		FT-3-1102-025	FT-3-1103-025	FT3-1104-025
37	FT-3-1101-037	FT-3-01110-037			
47	FT-3-1101-047		FT-3-1102-047	FT-3-1103-047	FT3-1104-047
50	FT-3-1101-050	FT-3-01110-050	FT-3-1102-050	FT-3-1103-050	FT3-1104-050
55	FT-3-1101-055		FT-3-1102-055	FT-3-1103-055	
70	FT-3-1101-070	FT-3-01110-070	FT-3-1102-070	FT-3-1103-070	FT3-1104-070
80	FT-3-1101-080				
90	FT-3-1101-090	FT-3-01110-090	FT-3-1102-090	FT-3-1103-090	FT3-1104-090
100	FT-3-1101-100	FT-3-01110-100	FT-3-1102-100	FT-3-1103-100	FT3-1104-100
110	FT-3-1101-110	FT-3-01110-110	FT-3-1102-110	FT-3-1103-110	FT3-1104-110
125	FT-3-1101-125		FT-3-1102-125	FT-3-1103-125	FT3-1104-125
150	FT-3-1101-150		FT-3-1102-150	FT-3-1103-150	FT3-1104-150

Ø in mm	MGF (100 Pieces)	MGG (100 Pieces)	MG 550-HA (50 Pieces)	13440*
25	FT-3-1105-025	FT-3-1106-025		
42				13440-042-----Q
44				13440-044-----Q
47	FT-3-1105-047	FT-3-1106-047	FT-3-01147-047	13440-047-----Q
50	FT-3-1105-050	FT-3-1106-050	FT-3-01147-050	13440-050-----Q
55	FT-3-1105-055	FT-3-1106-055		
70	FT-3-1105-070	FT-3-1106-070	FT-3-01147-070	
90	FT-3-1105-090	FT-3-1106-090	FT-3-01147-090	
100				13440-100-----K
110	FT-3-1105-110	FT-3-1106-110	FT-3-01147-110	
125	FT-3-1105-125		FT-3-01147-125	
130				13440-130-----K
150	FT-3-1105-150			

\* Q = 500 pieces | K = 50 pieces

Other dimensions as well as sheets are available on request

## Glass Microfiber Filters

### With Binder

These filters are mostly used either for monitoring air and gas or as prefilter. They are manufactured with synthetic binding agents to ensure that the filter has a defined strength. They are mechanically and chemically stable, have a temperature resistance up to 180°C and – depending on the binding agent used – are either hydrophobic or hydrophilic.

### Technical Specifications

Grade	Weight (g/m <sup>2</sup> )	Thickness (mm)	Air resistance (mbar)	Retention Rating 0.12–0.3 µm Particulate Matter (%)*	Filter Class	Binding Agent
MG 227/1/60	60	0.27	21.0	99.3	E 12	Hydrophobic
MG 1336/2	73	0.40	28.0	99.93	H 13	Hydrophobic
MG 1387/1	90	0.40	31.0			Hydrophilic
13400	74	0.34				Hydrophilic
13430	223	1.30				Hydrophilic

\* Tested and classified according to the standards EN 1822-1:2009

### Ordering Information



#### Filter Discs

Ø in mm	13400**	13430**	MG 221/1/60 (100 Pieces)	MG 1336/2 (100 Pieces)	MG 1387/1 (50 Pieces)
13	13400--13-----S				
20	13400--20-----S				
42	13400--42-----Q				
44	13400--44-----Q				
45					FT-3-01125-045
47	13400--47-----Q			FT-3-01120-047	FT-3-01125-047
50	13400--50-----Q				FT-3-01125-050
55				FT-3-01120-055	FT-3-01125-055
80	13400--80-----N				
100	13400-100-----K	13430-100-----K			
110				FT-3-01120-110	FT-3-01125-110
120	13400-120-----K				
124	13400-124-----K				
125					FT-3-01125-125
127	13400-127-----K	13430-127-----K			
130	13400-130-----K	13430-130-----K			FT-3-01125-130
150				FT-3-01124-150	

\*\* K = 50 pieces, N = 100 pieces, Q = 500 pieces, S = 200 pieces

Other dimensions as well as sheets are available on request

## ■ Quartz Microfiber Filters

These quartz microfiber filters are free of glass fibers and binding agents. They are especially suited for emission monitoring at temperatures of up to 900°C and wherever filters of the highest purity are needed (such as for air monitoring according to US EPA PM 10 or trace element analyses).

### They are available in two grades:

- Grade T 293, quartz microfiber filters unconditioned
- Grade MK 360, quartz microfiber filters conditioned (heat pre-treated); Certificate on trace elements available for every batch.

### □ Technical Specifications

Grade	Material	Weight (g/m <sup>2</sup> )	Thickness (mm)	Retention 0.3 µm Particulate Matter (%)*	Penetration 0.3 µm Particulate Matter (%)*
T 293   MK 360	100% Quartz microfiber silicium dioxide (SiO <sub>2</sub> )	85	0.43	99.998	< 0.002

\* Tested and classified according to the standards EN 1822-1:2009

### □ Ordering Information

#### ○ Filter Discs

Ø in mm	T 293 (50 Pieces)	MK 360 (25 Pieces)
25	FT-3-1109-025	
37	FT-3-1109-037	FT-3-01108-037
45	FT-3-1109-045	FT-3-01108-045
47	FT-3-1109-047	FT-3-01108-047
50	FT-3-1109-050	FT-3-01108-050
70	FT-3-1109-070	
90	FT-3-1109-090	
100	FT-3-1109-100	
110	FT-3-1109-110	
125	FT-3-1109-125	
142	FT-3-1109-142	
150	FT-3-1109-150	
293	FT-3-1109-293	

Other dimensions as well as sheets are available on request

## Extraction Thimbles

Thimbles are supplied in three different thimble designs to cover the majority of application areas. The cellulose or glass microfiber thimbles are primarily used in Soxhlet extraction units to extract defined substances from solids for further analyses. Quartz microfiber thimbles are preferred for emission control due to their high temperature resistance. They are supplied in a large variety of diameters and lengths. Below you will find an overview of the most commonly used grades.

### Technical Specifications

	Grade 30	Grade 40	Grade MK 360
Material	100% Cotton linters	100% Borosilicate glass	100% Quartz-microfiber (SiO <sub>2</sub> )
Inner diameter in mm*	+0 / -3	+1 / -3	+0 / -3
Outer length in mm*	± 1	± 1	± 1
Wall thickness in mm*	1.5 ± 0.5	2 ± 0.5	2 ± 0.5
Penetration % (0.3 µm)**		< 0.002	< 0.002
Max. temperature resistance	120 °C	500 °C	900 °C

\* According to the German DIN standard 12449, the tolerances of the inner diameter, the length and the wall thickness depend on the size of the extraction thimble. The tolerances above are for thimbles with an inner diameter < 48 mm.

\*\* According to the Standard EN 1822-1

### Ordering Information




#### Thimbles, 25 Pieces

Inner Ø × Outer Length (mm)	Grade 30	Grade 40	Grade MK 360
9 × 50	FT-1201-009050		
10 × 50	FT-1201-010050	FT-1204-010050	
22 × 80	FT-1201-022080	FT-1204-022080	
25 × 80	FT-1201-025080		
25 × 100	FT-1201-025100	FT-1204-025100	FT-01208-025100
26 × 60	FT-1201-026060	FT-1204-026060	
28 × 80	FT-1201-028080		
28 × 100	FT-1201-028100		
30 × 80	FT-1201-030080	FT-1204-030080	
30 × 100	FT-1201-030100	FT-1204-030100	FT-01208-030100
33 × 80	FT-1201-033080	FT-1204-033080	
33 × 94	FT-1201-033094		
33 × 100	FT-1201-033100		
33 × 118	FT-1201-033118		
33 × 205	FT-1201-033205		
43 × 123	FT-1201-043123	FT-1204-043123	FT-01208-043123
44 × 230		FT-1204-044230	
54 × 165	FT-1201-054165		

Other dimensions are available on request

## Membrane Filtration – Quick Selection Guide

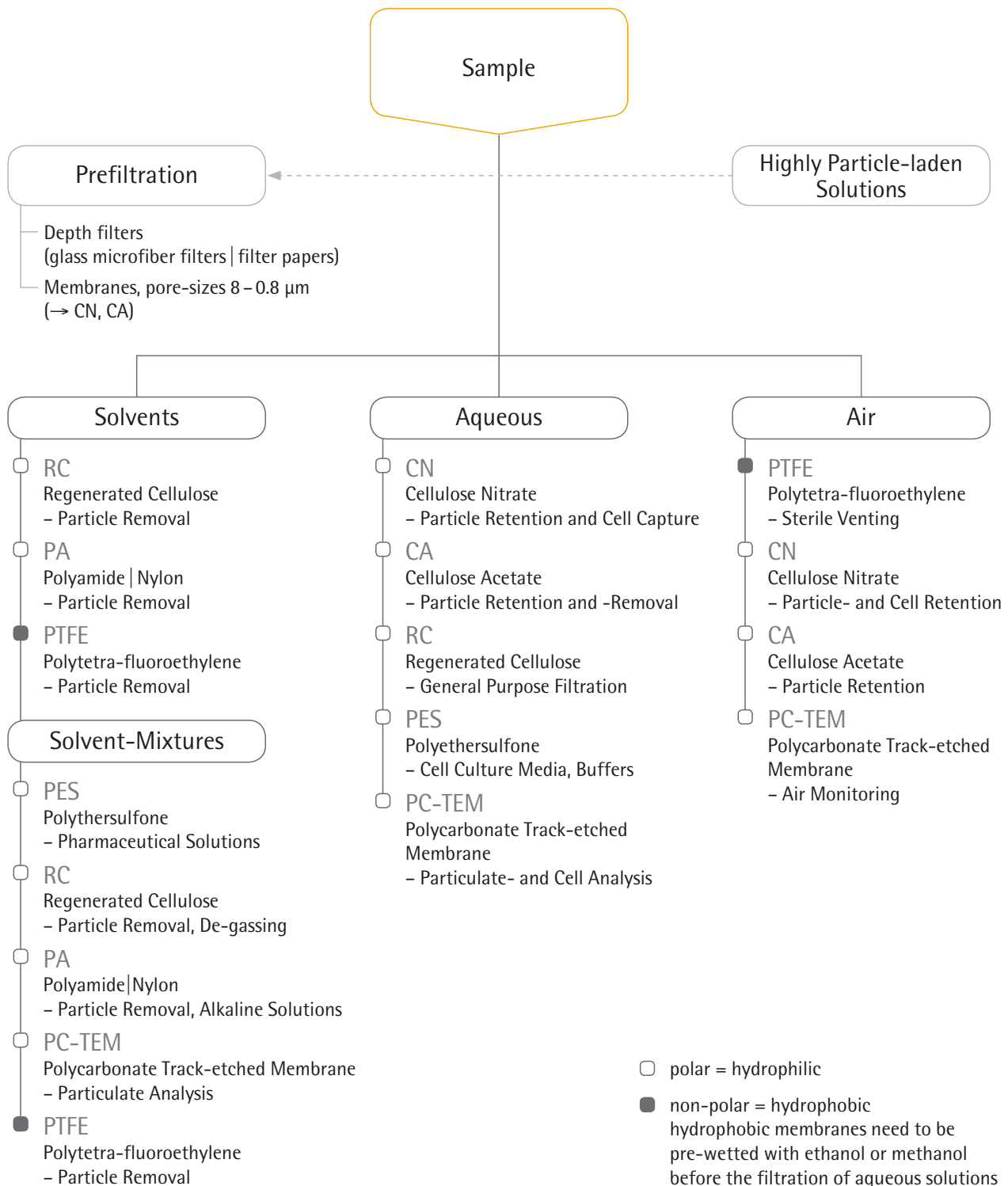

## Cellulose Nitrate

Cellulose nitrate membrane filters are indicated for many general laboratory applications where a membrane with a high non-specific adsorption is suitable. They are hydrophilic, have high flow rates thanks to their symmetric structure and are compatible with aqueous solutions (pH 4-8), hydrocarbons and several other organic solvents. The cellulose nitrate membranes are available in different pore sizes from 0.2  $\mu\text{m}$  to 8  $\mu\text{m}$ .

### Technical Specifications

Type	Pore Size ( $\mu\text{m}$ )	Thickness ( $\mu\text{m}$ )	Bubble Point (bar)	Water Flow Rate (ml/min/cm <sup>2</sup> /bar)	Burst Pressure (bar)
11327	0.2	130	4.2	25	$\geq 0.35$
11306	0.45	130	2.4	70	$\geq 0.3$
11305	0.65	130	2	130	$\geq 0.25$
11304	0.8	130	1.4	200	$\geq 0.2$
11303	1.2	130	1	200	$\geq 0.2$
11302	3	130	0.5	430	$\geq 0.2$
11342	5	130	0.5	570	$\geq 0.15$
11301	8	130	0.3	750	$\geq 0.1$

### Ordering Information

#### Filter Discs

$\varnothing$ in mm	11301*	11302*	11303*	11304*
13	11301--13-----N	11302--13-----N	11303--13-----N	11304--13-----N
20				11304--20-----N
25	11301--25-----N	11302--25-----N	11303--25-----N	11304--25-----N
30				11304--30-----N
37	11301--37-----N			11304--37-----N
47	11301--47-----N	11302--47-----N	11303--47-----N	11304--47-----N
50	11301--50-----N	11302--50-----N	11303--50-----N	11304--50-----N
68		11302--68-----G		
70	11301--70-----G			
90		11302--90-----G	11303--90-----G	11304--90-----G
100	11301-100-----N	11302-100-----G	11303-100-----G	11304-100-----G
$\varnothing$ in mm	11305*	11306*	11327*	11342*
13	11305--13-----N	11306--13-----N		11342--13-----N
20		11306--20-----N		
25	11305--25-----N	11306--25-----N	11327--25-----N	11342--25-----N
30		11306--30-----N		
37		11306--37-----N		
47	11305--47-----N	11306--47-----N		11342--47-----N
50	11305--50-----N	11306--50-----N		11342--50-----N
55		11306--55-----N		
85		11306--85-----N		
90		11306--90-----N		11342--90-----G
100	11305-100-----N	11306-100-----N		11342-100-----G
110		11306-110-----N		

\* G = 25 pieces, N = 100 pieces

Other dimensions and packaging units are available on request

## Cellulose Acetate

Cellulose acetate membranes combine high flow rates and thermal stability with very low adsorption characteristics, and are therefore excellently suited for use in pressure filtration devices. They are hydrophilic, have high flow rates thanks to their symmetric structure and are compatible with aqueous solutions (pH 4–8), oils, alcohols and other organic solvents. The 0.2 µm membrane is the filter of choice for sterile filtration of aqueous solutions, such as nutrient media, buffers and sera. The cellulose acetate membranes are available in different pore sizes from 0.2 to 5 µm.

### Technical Specifications

Type	Pore Size (µm)	Thickness (µm)	Bubble Point (bar)	Water Flow Rate (ml/min/cm <sup>2</sup> /bar)	Burst Pressure (bar)
11107	0.2	120	2.9	24	0.8
11106	0.45	120	1.9	69	0.7
11105	0.65	120	1.5	115	0.7
11104	0.8	120	1	200	0.5
12303	1.2	140	0.8	320	0.4
12342	5	140	0.4	570	0.25

### Ordering Information

#### Filter Discs

Ø in mm	11104*	11105*	11106*	11107*	12303*	12342*
13	11104--13-----N		11106--13-----N	11107--13-----N		
25	11104--25-----N	11105--25-----N	11106--25-----N	11107--25-----N	12303--25-----N	12342--25-----N
30			11106--30-----N	11107--30-----N		
37			11106--37-----N			
45						
47	11104--47-----N	11105--47-----N	11106--47-----N	11107--47-----N	12303--47-----N	12342--47-----N
50	11104--50-----N	11105--50-----N	11106--50-----N	11107--50-----N	12303--50-----N	
70						
85			11106--85-----N			
90			11106--90-----G	11107--90-----G		
100			11106--100-----N	11107--100-----N	12303--100-----G	
110			11106--110-----N			

\* G = 25 pieces, N = 100 pieces

Other dimensions and packaging units are available on request

## Regenerated Cellulose

The very low adsorption membranes are hydrophilic, solvent-resistant (pH 3-12) and therefore suited for the particle removal from solvents. The membrane is asymmetric and reinforced with nonwoven cellulose. They are available in two pore sizes: 0.45 µm and 0.2 µm.

### Technical Specifications

Type	Pore Size (µm)	Thickness (µm)	Bubble Point (bar)	Water Flow Rate (ml/min/cm <sup>2</sup> /bar)
18407	0.2	170	4.4	15
18406	0.45	170	2.9	30

### Ordering Information

#### Filter Discs

Ø in mm	18406*	18407*
13	18406--13-----N	18407--13-----N
25	18406--25-----N	18407--25-----N
47	18406--47-----N	18407--47-----N
50	18406--50-----N	18407--50-----N
100	18406-100-----G	18407-100-----G

\* G = 25 pieces, N = 100 pieces

Other dimensions and packaging units are available on request

## Polyethersulfone

Polyethersulfone (PES) membrane filters are hydrophilic, have high flow rates thanks to their symmetric structure, have a low non-specific protein adsorption and are chemically resistant over a pH range of 1 – 14. They are therefore recommended for the filtration of aqueous solutions as well for the protein filtration. Furthermore, the low level of extractables makes them suitable for environmental analysis.

### Technical Specifications

Type	Pore Size (µm)	Thickness (µm)	Bubble Point (bar)	Water Flow Rate (ml/min/cm <sup>2</sup> /bar)	Burst Pressure (bar)
15458	0.1	150	3.8	10	≥ 0.6
15407MI	0.2	150	3.5	25	≥ 0.5
15406	0.45	150	2.6	46	≥ 0.5

### Ordering Information

#### Filter Discs

Ø in mm	15406*	15407*	15458*
25	15406--25-----N	15407--25----MIN	15458--25-----N
47	15406--47-----N	15407--47----MIN	15458--47-----N
50	15406--50-----N	15407--50----MIN	15458--50-----N
90		15407--90----MIK	

\*K = 50 pieces, N = 100 pieces

Other dimensions are available on request

## ■ Polyamide

Polyamide membrane filters are hydrophilic and chemically resistant to alkaline solutions and organic solvents. They are therefore recommended for the particle removal from aqueous solutions and solvents for analytical determination such as HPLC, for the sterile filtration of these liquids as well as for applications where a membrane with a relatively high non-specific adsorption is suitable.

### □ Technical Specifications

Type	Pore Size (µm)	Thickness (µm)	Bubble Point (bar)	Water Flow Rate (ml/min/cm <sup>2</sup> /bar)	Burst Pressure (bar)
25007	0.2	115	3.2	15	≥ 0.25
25006	0.45	115	2.3	35	≥ 0.23

### □ Ordering Information

#### ○ Filter Discs

Ø in mm	25006*	25007*
13	25006--13-----N	25007--13-----N
25	25006--25-----N	25007--25-----N
47	25006--47-----N	25007--47-----N
50	25006--50-----N	25007--50-----N
90	25006--90-----G	25007--90-----G

\* G = 25 pieces, N = 100 pieces

Other dimensions and packaging units are available on request

## ■ Hydrophobic PTFE

The main application of these membrane filters is the filtration of air, gases or chemicals. They are made of PTFE (polytetrafluorethylene) only and are therefore permanently hydrophobic. Unlike other (hydrophilic) filter types, they are not wetted by air humidity, allowing unhindered passage of air at low differential pressures as well. PTFE membrane filters have an excellent chemical compatibility (pH 1-14), so that they are also used for the filtration of solvents and acids, to which other filter types are not resistant. Due to their hydrophobic characteristics, they must be pre-wetted with ethanol or methanol before the filtration of aqueous media.

### □ Technical Specifications

Type	Pore Size (µm)	Thickness (µm)	Bubble Point (bar)	Isopropanol Flow Rate (ml/min/cm <sup>2</sup> /bar)
11807	0.2	65	1.0	11
11806	0.45	80	0.8	20
11803	1.2	100	0.45	80
11842	5	100	0.10	210

### Ordering Information



#### Filter Discs

Ø in mm	11803*	11806*	11807*	11842*
13	11803--13-----N	11806--13-----N	11807--13-----N	
25	11803--25-----N	11806--25-----N	11807--25-----N	11842--25-----N
37		11806--37-----N		
42				11842--42-----N
47	11803--47-----N	11806--47-----N	11807--47-----N	11842--47-----N
50	11803--50-----N	11806--50-----N	11807--50-----N	11842--50-----N
90	11803--90-----G	11806--90-----G	11807--90-----G	
100	11803-100-----G	11806-100-----G	11807-100-----G	11842-100-----G

\*G= 25 pieces, N= 100 pieces

Other dimensions and packaging units are available on request

## Polycarbonate Track-Etched

Those white and hydrophilic polycarbonate track-etched membranes are manufactured from high grade polycarbonate film using track-etch technology. Their capillary pore structure is uniform and precise, with a narrow pore size distribution to retain particles on their surface. Track-etched membranes are an excellent choice for accurate fractionation of particulates because of their precise pore size. Track-etch technology offers the user distinct performance advantages when excellent surface capture and high sample visibility are required. Their main applications are particulate analysis, epifluorescence microscopy, fluid clarification, cytology, cell biology, bioassays, water microbiology and environmental analysis.

### Technical Specifications

Type	Pore Size (µm)	Thickness (µm)	Bubble Point (bar)	Water Flow Rate (ml/min/cm <sup>2</sup> /0.7 bar)	Burst Pressure (bar)
23007	0.2	25	4.8	≥ 10	≥ 0.7
23006	0.4	25	2.5	≥ 30	≥ 0.7

### Ordering Information



#### Filter Discs, 100 Pieces

Ø in mm	23006	23007
25	23006--25-----N	23007--25-----N
47	23006--47-----N	23007--47-----N
50		23007--50-----N

Other dimensions are available on request

## ■ Blotting | Chromatography Papers

These papers are made of cotton linters only with  $\alpha$ -cellulose content of  $> 98\%$ . These highly pure papers are not only ideal for blotting and chromatography, but also for a wide range of absorption applications like those common in the life sciences and diagnostics. Below you will find an overview of the most commonly used grades.

### □ Technical Specifications

Grade	Weight (g/m <sup>2</sup> )	Thickness (mm)	Capillary Rise (mm/30 min)	Capillary Rise (mm/10 min)	Properties
FN 4	125	0.24	95		Chromatography paper, ash content $< 0.04\%$
FN 7	150	0.32	145		Chromatography paper, ash content $< 0.04\%$
FN 30	320	0.90	240		Chromatography paper, ash content $< 0.04\%$ , paper for antibiotic test strips
FN 100	195	0.35	115	70	The most commonly used chromatography and blotting paper
BF 3	330	0.76	30	130	Blotting paper to increase and maintain the transport of liquids

### □ Ordering Information



**Sheets in 580 × 600 mm**

Grade FN 4 (100 Sheets)	Grade FN 7 (50 Sheets)	Grade FN 30 (25 Sheets)	Grade FN 100 (50 Sheets)	Grade BF 3 (50 Sheets)
FT-2-504-580600N	FT-2-507-580600K	FT-2-526-580600G	FT-2-527-580600K	FT-2-520-580600K

Other dimensions and shapes are available on request

## ■ Nitrocellulose Membrane for Blotting



### □ Technical Specifications

**Typical Parameter Value for NC 0.22  $\mu\text{m}$  and 0.45  $\mu\text{m}$**

	0.22 $\mu\text{m}$	0.45 $\mu\text{m}$
Material	Cellulose nitrate	Cellulose nitrate
Thickness	120 $\mu\text{m}$	130 $\mu\text{m}$
Water flow rate	27 ml/(min.cm <sup>2</sup> bar)	70 ml/(min.cm <sup>2</sup> bar)
Bubble point	4.4 bar	2.4 bar
Wettability in water	$\leq 1$ s	$\leq 1$ s
Extractable content in water	$\leq 1\%$	$\leq 1\%$
Burst pressure	0.8 bar	0.2 bar
Binding capacity for IgG	200 $\mu\text{g}/\text{cm}^2$	200 $\mu\text{g}/\text{cm}^2$

### □ Ordering Information

	Roll Size	Order No.
NC 0.22 $\mu\text{m}$	30 cm × 3 m	11327-----41BL
NC 0.45 $\mu\text{m}$	30 cm × 3 m	11306-----41BL

All indicated data to be understood as typical average values

## ■ Re-usable, 13 mm Syringe Filter Holders

For the Ultracleaning of Small Volumes Up to About 10 ml



### PTFE Holder for Solvents and Chemicals

Made completely of PTFE, this holder is unaffected by chemicals and contains no trace elements which could be released into the liquid being filtered. It is therefore extremely well suited for particle removal from samples and reagents for analytical methods, such as NMR samples. Other

benefits of this application are the low hold-up volume, the easy cleaning and the drying at a temperature of 180°C. The construction of the holder ensures leak proof sealing without a sealing ring, and avoids twisting of the membrane filter when the top is tightened onto the base.

### □ Technical Specifications

Connectors	Female Luer Lock inlet, luer slip outlet
Chemical compatibility	As for PTFE
Filtration area	0.5 cm <sup>2</sup>
Materials	PTFE top and bottom parts
Max. operating pressure	5 bar   500 kPa   72.5 psi
Membrane filter diameter	13 mm
Sterilization	By autoclaving (max. 134 °C) or by dry heat (max. 180 °C)
Hold-up volume	Less than 0.03 ml after overcoming the bubble point (0.3 ml before)

### □ Ordering Information

Description	Order No.
13 mm PTFE Syringe Filter Holder	16574



### Polycarbonate Holder for Aqueous Solutions

This inexpensive filter holder is made of clear, autoclavable polycarbonate. The silicone gasket enables a leak-free filtration

at pressures of up to 7 bar by simply manually screwing it together. Filter supports in the top and bottom parts allow filtration in either direction.

### □ Technical Specifications

Connectors	Female Luer Lock inlet, luer slip outlet
Chemical compatibility	As for polycarbonate and silicone
Filtration area	0.5 cm <sup>2</sup>
Materials	Polycarbonate top and bottom part, silicone gasket
Max. operating pressure	7 bar   700 kPa   101.57 psi
Membrane filter diameter	13 mm
Sterilization	By autoclaving at 121 °C
Hold-up volume	Less than 0.2 ml after overcoming the bubble point (0.3 ml before)

### □ Ordering Information

Description	Order No.
13 mm Polycarbonate Syringe Filter Holder, pack of 12	16514-----E
Silicon gasket, 10 × 14.9 mm, pack of 10	6980569



## ■ Re-usable 25 mm Syringe Filter Holders

For the Ultracleaning and Sterilizing Filtration of Volumes of Up to About 100 ml



### Stainless Steel Holder for Solvents and Chemicals

Made of stainless steel, this holder is heat-resistant, and the chemical compatibility depends only on the inserted filter type.

The top part can easily be mounted on the bottom part using the enclosed tightening tool. Filter supports in the top and bottom parts allow filtration in either direction.

### □ Technical Specifications

Connectors	Female Luer Lock inlet, luer slip outlet
Chemical compatibility	As for stainless steel
Filtration area	3 cm <sup>2</sup>
Materials	Stainless steel (1.4305) top and bottom parts
Max. operating pressure	7 bar   700 kPa   101.5 psi
Membrane filter diameter	25 mm
Sterilization	By autoclaving (max. 134 °C) or by dry heat (max. 180 °C)
Hold-up volume	Less than 0.1 ml after overcoming the bubble point (0.3 ml before)

### □ Ordering Information

Description	Order No.
25 mm Stainless Steel Holder	16214
Tightening tool, Polyman 24/5	6980595



### Polycarbonate Holder for Aqueous Solutions

This inexpensive filter holder is made of clear, autoclavable polycarbonate. The silicone gasket enables a leak-free filtration

at pressures of up to 7 bar by simply manually screwing it together. Filter supports in the top and bottom parts allow filtration in either direction.

### □ Technical Specifications

Connectors	Female Luer Lock inlet, luer slip outlet
Chemical compatibility	As for polycarbonate and silicone
Filtration area	3 cm <sup>2</sup>
Materials	Polycarbonate top and bottom parts, silicone gasket
Max. operating pressure	7 bar   700 kPa   101.5 psi
Membrane filter diameter	25 mm
Sterilization	By autoclaving at 121 °C
Hold-up volume	Less than 0.3 ml after overcoming the bubble point (0.6 ml before)

### □ Ordering Information

Description	Order No.
25 mm Polycarbonate Syringe Filter Holder, pack of 12	16517-----E
Silicone gasket, 20×25 mm, pack of 10	1EDS-D0053

■ 25 mm Glass Vacuum Filter Holder

For Hybridisation Tests, Particle Testing and Clarification



This filter holder is available in two versions differing from each other only in the type of the filter support. The filter with glass frit ensures uniform distribution of retained particles and is therefore recommended, when the residue on the filter surface is of interest. Because it is easy to clean, the device with the PTFE-coated screen support is preferable when the filtrate is required, or when liquids difficult to remove from glass frits must be

examined. The PTFE ring, which holds the glass frit and the screen support, allows for the autoclaving of the devices with a filter in position and protects the edge of the glass frit from breakage and potential leakage. It has a rim around the upper edge to simplify the positioning of the membrane filter when inserted and a silicone O-ring in the underside for a leak-proof seal on the filtrate side. The funnel-shaped top part simplifies filling in the sample.


□ Technical Specifications

Outlet spout	12 mm Ø
Parts and materials	Borosilicate glass funnel and base PTFE   glass filter support (type 16306) or PTFE   stainless steel filter support, coated with PTFE (type 16315) Silicone O-ring 25×3 mm Anodised Aluminium clamp
Chemical compatibility	As for glass, PTFE and silicone. The silicone O-ring can be replaced by a Viton O-ring (order no. 00118)
Funnel capacity	30 ml
Filtration area	3 cm <sup>2</sup>
Max. operating pressure	Only for vacuum
Suitable membrane filter diameter	25 mm (or 24 mm)
Sterilization	By autoclaving (max. 134 °C) or by dry heat (max. 180 °C)

□ Ordering Information

Description	Order No.
Glass vacuum filtration holder for 25 mm (or 24 mm) membrane filter, with glass frit filter support	16306
Glass vacuum filtration holder for 25 mm (or 24 mm) membrane filter, with PTFE-coated screen filter support	16315

Replacement parts are shown in the diagram.

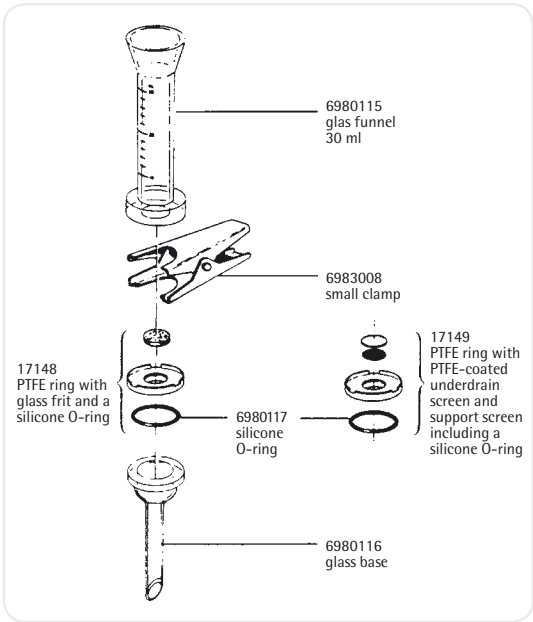

## 50 mm Glass Vacuum Filter Holder

For Particle Testing or Clarification and Sterile Filtration



This filter holder is available in two versions differing from each other only in the type of the filter support. The device with glass frit ensures uniform distribution of retained particles and is therefore recommended, when the residue on the filter surface is of interest. Because it is easy to clean, the device with the PTFE-coated screen support is preferable when the filtrate is required, or when liquids difficult to remove from

glass frits must be examined. The PTFE ring, which holds the glass frit and the screen support, allows the autoclaving of the devices with a filter in position and protects the edge of the glass frit from breakage and potential leakage. It has a rim around the upper edge to simplify the positioning of the membrane filter when inserted and a silicone O-ring in the underside for a leak-proof seal on the filtrate side.




### Technical Specifications

Outlet spouts	15 mm Ø
Parts and materials	Borosilicate glass funnel and base Silicone caoutchouc lid PTFE   glass filter support (type 16307) or PTFE   stainless steel filter support, coated with PTFE (type 16316) Silicone O-ring 45×3 mm Anodised Aluminium clamp
Chemical compatibility	As for glass, PTFE and silicone. The silicone O-ring can be replaced by a Viton O-ring (order no. 00124).
Funnel capacity	250 ml
Filtration area	12.5 cm <sup>2</sup>
Max. operating pressure	Only for vacuum
Suitable membrane filter diameter	50 mm (or 47 mm)
Sterilization	By autoclaving (max. 134 °C) or by dry heat (max. 180 °C)

### Ordering Information

Description	Order No.
Glass vacuum filtration holder for 50 mm (or 47 mm) membrane filter, with glass frit filter support	16307
Glass vacuum filtration holder for 50 mm (or 47 mm) membrane filter, with PTFE-coated screen filter support	16316

Replacement parts are shown in the diagram.

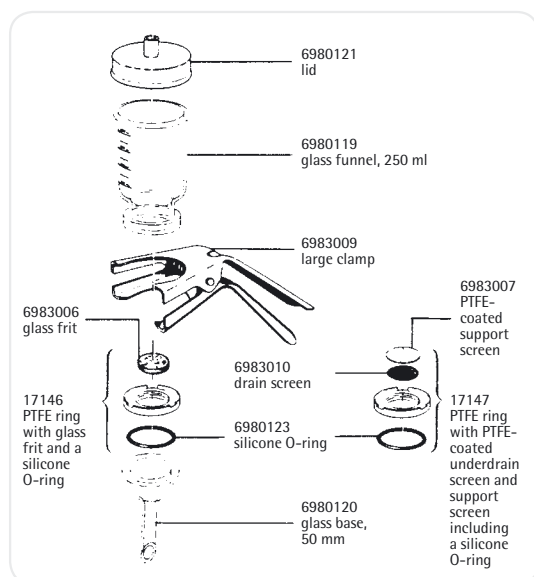

# All-Glass Vacuum Filter Holder

For Analytical Determinations, Particle Removal from Solvents



All areas, where liquid and device can come into direct contact, are made of glass or PTFE. The device, in combination with solvent-resistant, hydrophilic RC-membranes, is therefore ideal for ultracleaning and degassing solvents and solvent mixtures for HPLC, GC and AA. Convenience of handling is ensured by several beneficial features. A 6 mm wide non-ground rim

above the ground glass neck of the suction flask prevents the filtrate from contacting grease on the ground glass surface and so avoids its contamination while being poured out of the flask. The hose nipple connector is made of polypropylene for safe connection of the vacuum hose. The filtrate outlet spout ends well below the entrance to this hose nipple.




## Technical Specifications

Parts and materials	Borosilicate glass funnel, base and flask, sintered glass frit in a PTFE ring and viton O-ring (45 × 3 mm) underneath, anodized aluminium clamp
Chemical compatibility	As for glass and PTFE
Funnel capacity	250 ml
Capacity of the filtrate flask	1 liter
Filtration area	12.5 cm <sup>2</sup>
Max. operating pressure	Only for vacuum
Suitable membrane filter diameter	50 mm (or 47 mm), 40 mm prefilter
Sterilization (without connector)	By autoclaving (max. 134 °C) or by dry heat (max. 180 °C)

## Ordering Information

Description	Order No.
All-glass vacuum filter holder for 50 mm (or 47 mm) membrane filter, with vacuum-resistant flask, capacity 1 liter	16309

Replacement parts are shown in the diagram.

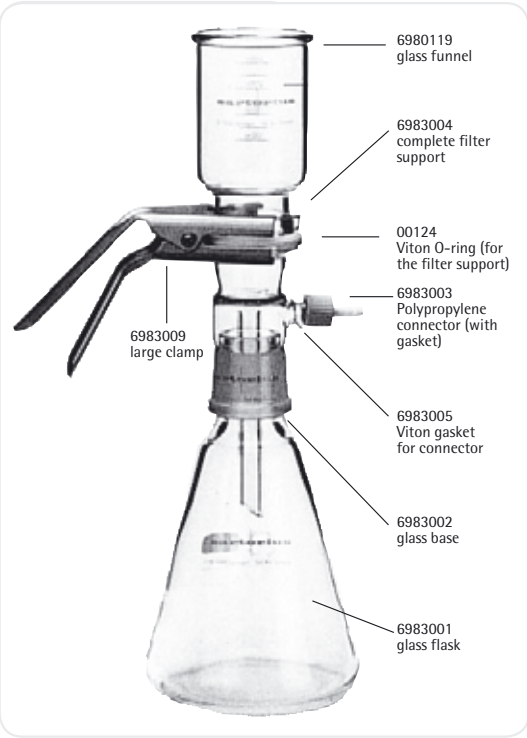

## Polycarbonate In-Line Filter Holder

For the Filtration of Liter Volumes of Aqueous Solutions



This holder is made of stable, autoclavable polycarbonate. This practical holder is suitable for many simple laboratory filtrations. It can be connected to a peristaltic pump or a pressure container. The bell-shaped base protects the filtrate from repeated contamination while flowing in a receiver. The

holder is characterized by an excellent resistance to pressure and density setting by simply hand-tightening. The transparent top part allows the visual control of the correct fit of the O-ring. The hose nipples can be replaced by luer connectors to use it as a large area syringe filter holder.


### Technical Specifications

Chemical compatibility	As for polycarbonate, polypropylene and silicone
Filtration area	12.5 cm <sup>2</sup>
Weight	83 g
Threads for connectors	M 12 × 1 female thread
Materials	Polycarbonate top part, base part and hose nipple, polypropylene filter support, silicone O-ring (40 × 5 mm)
Max. operating pressure	7 bar   700 kPa   101.5 psi
Suitable membrane filter diameter	50 mm (40 mm prefilter)
Sterilization	By autoclaving at 121 °C The material withstands repeated cycles, provided aggressive cleaning agents are completely washed off and that the boiler water does not contain anti-corrosive or anti-scaling additives.

### Ordering Information

Description	Order No.
Polycarbonate in-line filter holder for 50 mm membrane filter, pack of 5.	16508-----B

Replacement parts are shown in the diagram.

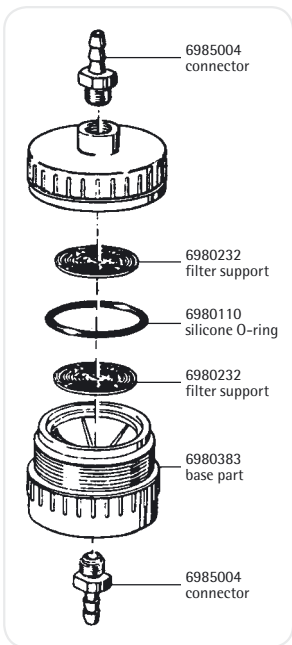

■ 25 mm Stainless Steel Filter Holder

For In-Line Filtration



The G $\frac{1}{4}$  connection threads with density barrel, guarantee leak-proof sealing of the hose nipple and the holder without sealing rings. Other connectors, available as accessories, fit the holder onto reducing

valves or pumps with G $\frac{1}{4}$  female thread (order no. 01030) or G $\frac{3}{8}$  female thread order no. 01029) or onto pressure tanks with G $\frac{3}{8}$  male thread (order no. 00177).


□ Technical Specifications

Connectors	Hose nipples DN10
Filtration area	3 cm <sup>2</sup>
Weight	ca. 170 g
Materials	Stainless steel, except silicone O-ring (21 × 2 mm) and aluminium closing ring
Max. operating pressure	5 bar   500 kPa   72.5 psi
Suitable membrane filter	25 mm (20 mm prefilter for the filtration of liquids only)
Sterilization	By autoclaving (max. 134 °C) or by dry heat (max. 180 °C)

□ Ordering Information

Description	Order No.
Stainless steel pressure filter holder for 25 mm Ø membrane filter.	16251

Replacement parts are shown in the diagram.

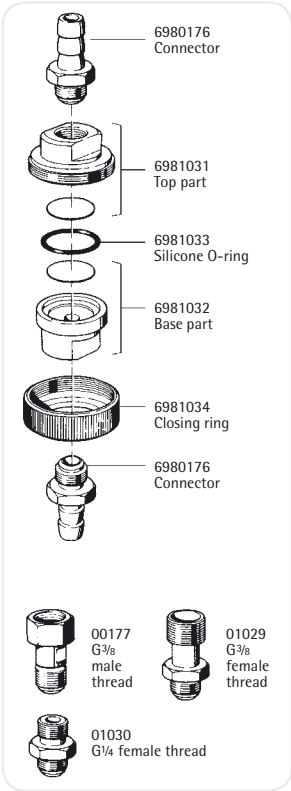

## 47 mm Stainless Steel Filter Holder

For In-Line Filtration



The filter holder is suitable for a pressure of up to 20 bar. The inlet side valve is convenient for the intermittent run-off of waste water. Other connectors, available as accessories, fit the holder onto reducing

valves or pumps with G $\frac{3}{8}$  female thread (order no. 17089) or onto pressure tanks with G $\frac{3}{8}$  male thread (order no. 17069) or on taps with G $\frac{3}{4}$  male thread (order no. 17068).


### Technical Specifications

Connectors	Hose nipples DN10
Connection thread	M12 × 1
Filtration area	13 cm <sup>2</sup>
Weight	ca. 490 g
Materials	Stainless steel, except silicone O-ring (42 × 3 mm), PTFE and Viton valve seals
Max. operating pressure	20 bar   2,000 kPa   290 psi
Suitable membrane filter	47 mm (40 mm prefilter)
Sterilization	By autoclaving (max. 134 °C) or by dry heat (max. 180 °C)

### Ordering Information

Description	Order No.
Stainless steel filter holder for 47 mm Ø membrane filter.	16254

Replacement parts are shown in the diagram.

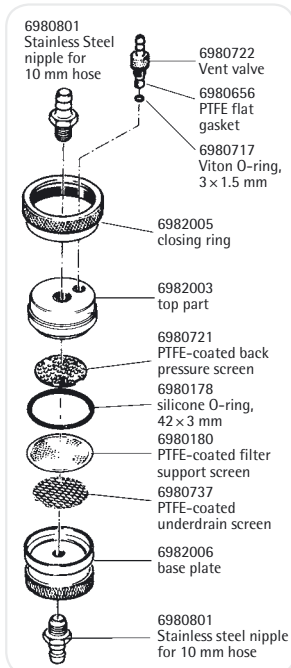

# Chemical-resistant PTFE Filter Holder

For the Filtration of Aggressive Liquids



The holder hinders the release trace elements into the filtrate and is resistant to almost all chemicals. The Viton O-ring in the top part allows easy hand tightening, and can be replaced by a PTFE O-ring,

order no. 17039. The 6 mm outlet nipple is an integral part of the base, the 10 mm inlet hose nipple can be replaced by a G<sup>3</sup>/<sub>8</sub> connector, order no. 17051.

## Technical Specifications

Chemical compatibility	As for PTFE and Viton
Filtration area	12.5 cm <sup>2</sup>
Thread for inlet connector	M 14 × 1.5 male thread
Materials	Top part, barrel, base part: corrugated iron, hose nipples and filter support with 40 × 3.5 mm O-ring: PTFE, locking rings: aluminium 39 × 3.5 mm Viton O-ring (top part)
Max. operating pressure	5 bar   500 kPa   72.5 psi
Suitable membrane filter diameter	47 mm
Sterilization	By autoclaving (max 134 °C) or by dry heat (180 °C)


## Ordering Information

Description	Order No.
PTFE pressure filter holder, 47 mm, with 200 ml capacity.	16579

## Replacement Parts


Description	Order No.
PTFE O-ring, 40 × 3.5 mm	6985000
PTFE Connector	6985002
PTFE filter support screen, 47 mm	6985001
Viton O-ring, 39 × 3.5 mm	6985011
PTFE O-ring, 39 × 3.5 mm	17039
G <sup>3</sup> / <sub>8</sub> connector	17051

Other replacement parts are shown in the diagram.

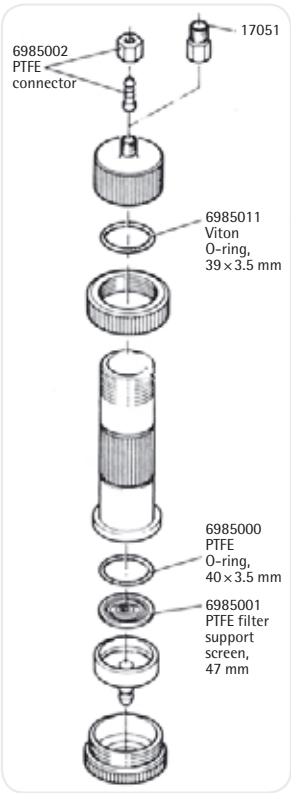

# Stainless Steel Pressure Filter Holder

For the Filtration of Up to 5 Liter Volumes

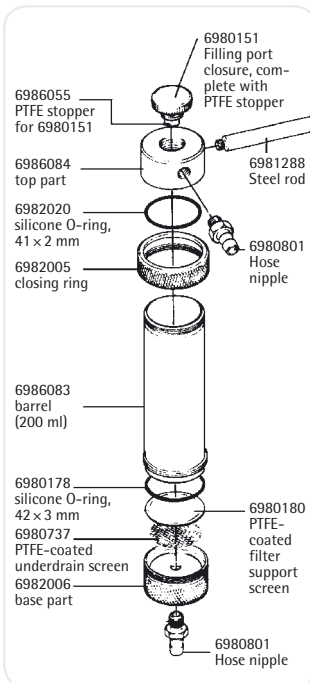

A practical filter holder for many laboratory filtrations. It can be attached to a tripod with the help of a steel rod which can be screwed in. The hose nipple is screwed into the side of the top part, leaving room for a large filling opening. This makes pouring in the sample easier, and the sample can be refilled without removing the tube connection to the pressure source. Leak-proof sealing is achieved by hand-tightening the

closing ring. For the filtration of small volumes (up to about 200 ml of soil samples or viscous liquids, such as oils), the holder is connected directly to a pressure source. For the filtration of up to 5 liter volumes of relatively easily filterable liquids (e.g. buffer solutions, solutions for cell counters and tissue culture solutions), it is used in combination with a pressure tank.

## Technical Specifications

Chemical compatibility	As for stainless steel, PTFE and silicone. If required, the silicone O-ring in the filter support can be replaced by a Viton O-ring 00179 or a PTFE O-ring 17038 (by reducing the max. operating pressure to 4 bar   58 psi); the silicone O-ring in the top part can be replaced by a Viton O-ring 17145.
Filtration area	13 cm <sup>2</sup>
Weight	960 g
Threads for connectors	M 12 x 1 female thread
Materials	Top part, barrel, base part, corrugated iron, closing ring, closure cap, back pressure screen and stainless steel hose nipples 1.4401 (AISI 316), PTFE-coated stainless steel filter support, silicone O-rings, 41 x 2 mm (top part) and 42 x 3 mm (filter support), PTFE-sealing (cap).
Max. operating pressure	10 bar   1,000 kPa   145 psi
Suitable membrane filter diameter	47 mm (prefilter, 42 mm)
Sterilization	By autoclaving (max 134 °C) or by dry heat (180 °C)

## Ordering Information

Description	Order No.
Stainless steel pressure filter holder	16249

## Replacement Parts

Description	Order No.
Viton O-ring, 42 x 3 mm	00179
PTFE O-ring, 42 x 3 mm	17038
Viton O-ring for upper part, 41 x 2 mm	17145

Other replacement parts are shown in the diagram.

## ■ Combisart® Manifolds

1-, 3- and 6-Branch



Made of high-grade stainless steel (B.S. 304S3 | AISI 304); accommodates any type of vacuum funnel. Stainless steel three-way valves (taps) allow the vacuum for each filter holder to be individually

controlled and each holder to be sterilely vented. The low height of the manifold ports is particularly advantageous for working on a clean bench.

### □ Ordering Information

Combisart® Manifolds, without Base Support and Frit	Order No.
Combisart® 1-branch manifold	16844
Combisart® 3-branch manifold	16842
Combisart® 6-branch manifold	16843

Combisart® Sets, Stainless Steel Capacity	Order No.
1-branch 1 × 100 ml	16844-CS
1-branch 1 × 500 ml	16845-CS
3-branch 3 × 100 ml	16824-CS
3-branch 3 × 500 ml	16828-CS
6-branch 6 × 100 ml	16832-CS
6-branch 6 × 500 ml	16831-CS


In each set stainless steel funnels with lids are preassembled.

### □ Accessories and Replacement Parts

Description	Pack Size	Order No.
Plug, conical, to close the venting hole beside the 3-way-valve	10	6980225
Silicone O-ring for manifold female threads	3	6980235

## ■ Glass Filter Holders; 30, 250 ml

For Particle Counting

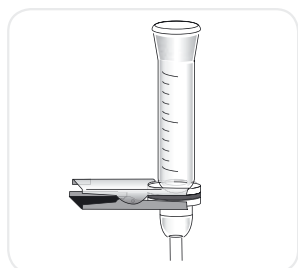


### Glass Filter Holders

Two compact vacuum filter holders for easy particulate analysis. Both the top and bottom part of the filter holders are easily and securely fastened together using the

metal clamp. The center-ing rim on the filter support ensures correct positioning of the membrane filter. The glass frit filter support guarantees uniform distribution of retained particles on the filter surface.

### □ Ordering Information

Order No.	Description	
16306	Glass filter holder	30 ml
	Filter diameter	25 mm (or 24 mm)
		Prefilter, 20 mm
	Filtration area	3 cm <sup>2</sup>
	Capacity	30 ml
	Outlet	12 mm outer diameter
16307	Glass filter holder	250 ml
	Filter diameter	47 mm (or 50 mm)
		Prefilter, 40 mm
	Filtration area	12.5 cm <sup>2</sup>
	Capacity	250 ml
	Outlet	15 mm outer diameter

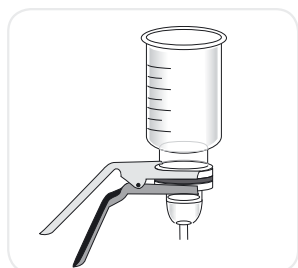


### Adapter, 16836 | Adapter, 16837

For use of a glass filter holder, 16306 or 16307, on a Combisart® stainless steel manifold.


### □ Ordering Information

Description	Order No.
Adapter with 11 mm opening in stopper; for using filter holder 16306 on a Combisart® manifold	16836
Replacement stopper for 16836	00280
Adapter with 14 mm opening in stopper; for using filter holder 16307 on a Combisart® manifold	16837
Replacement stopper for 16837	00281

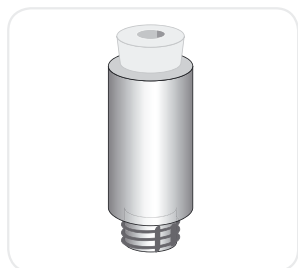

# Polycarbonate Filter Holders

For Particle Counting

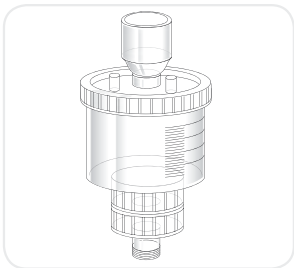

## Polycarbonate Filter Holder, 250 ml

This reusable, practical filter holder made of autoclavable plastic is ideal for analytical testing outside the laboratory. For use of 47 mm membrane filters.

Outlet: TR 20×2 mm male thread

### Ordering Information

Description	Order No.
Polycarbonate filter holder, 250 ml	16511

# Ready-to-Use Biosart® 250 Funnels

For Particle Counting



## Biosart® 250 Funnel

The Biosart® 250 Funnel has been specially designed for analytical quality assurance. The sterile 250 ml plastic funnel guarantees fast filtration and high sample throughputs

during routine testing. Its large inner diameter allows high flow rates, and the tapered inner walls permit thorough flushing of the funnel, after filtration.

### Ordering Information

Description	Order No.
Biosart® 250 Funnel, 50 units, sterile-packaged	16407-25-ALK



## Single Support, 16840

For adapting a Biosart® 250 Funnel for use on a Combisart® stainless steel manifold.


### Ordering Information

Description	Order No.
Stainless steel filter support for stainless steel manifold.	16840

## CONFIDENCE® – Validation Services

Our Services at a Glance

Multi-site laboratories offer globally integrated Validation Services.

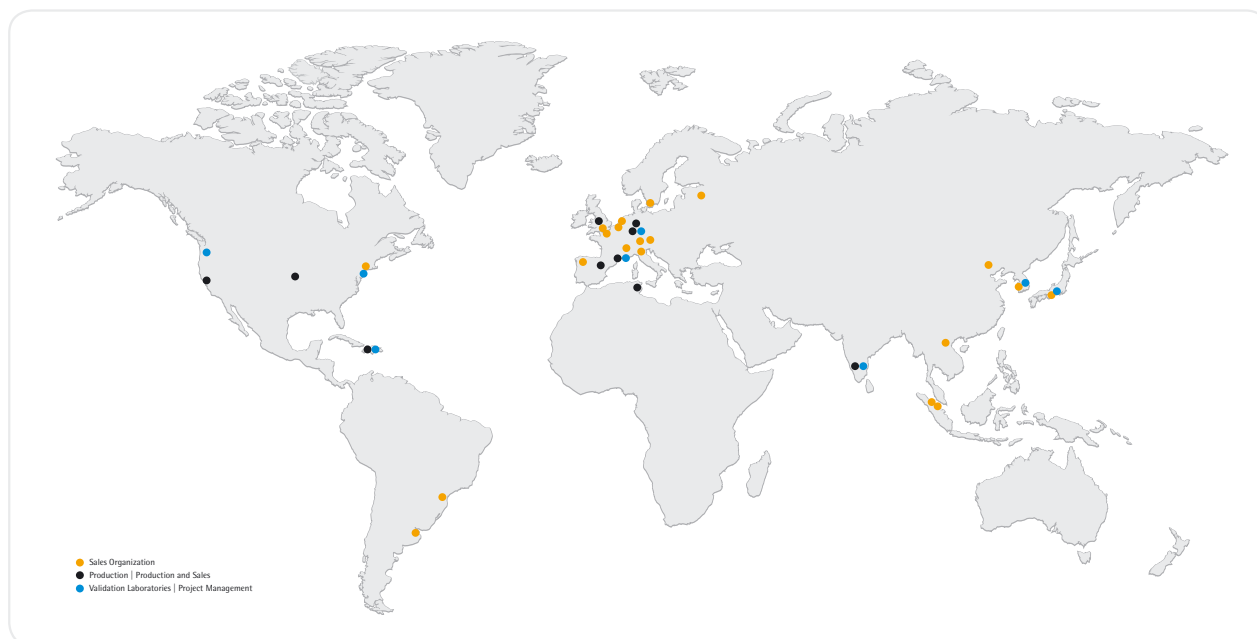

Every Sartorius Stedim Biotech product used in critical steps in the pharmaceutical industry is manufactured with complete traceability of materials. The extensive data obtained in our validation guides and extractables guides form a comprehensive basis for the safety assessment of your respective process and drug product. Given the latest developments in regulatory requirements, our in-depth knowledge of actual drug product testing is a key element of our CONFIDENCE® Validation Services.

Our product array includes, but is not limited to, the following:

Filter Elements	Fluid Management Containers	Polymer-based Process Components
Extractables   Leachables Studies	Extractables   Leachables Studies	Extractables   Leachables Studies
Physico-Chemical Studies	Physico-Chemical Studies	
Microbiological Studies	Microbiological Studies	

### Time To Market

Make the choice that helps bring your products to market rapidly – on time, every time. Ask us about our day commitment program.

Explore your possibilities – we continue to lead and innovate in the field of validation studies to increase the added value of our services offered to the biopharmaceutical and pharmaceutical industries.

Visit our microsite at [www.sartorius-stedim.com/confidence](http://www.sartorius-stedim.com/confidence)





# Sales and Service Contacts

For further contacts, visit [www.sartorius.com](http://www.sartorius.com)

## Europe

### Germany

Sartorius Stedim Biotech GmbH  
August-Spindler-Strasse 11  
37079 Goettingen

Phone +49.551.308.0  
Fax +49.551.308.3289

Sartorius Stedim Systems GmbH  
Robert-Bosch-Strasse 5 – 7  
34302 Guxhagen

Phone +49.5665.407.0  
Fax +49.5665.407.2200

### France

Sartorius Stedim FMT S.A.S.  
ZI Les Paluds  
Avenue de Jouques – CS 91051  
13781 Aubagne Cedex

Phone +33.442.845600  
Fax +33.442.845619

Sartorius Stedim France SAS  
ZI Les Paluds  
Avenue de Jouques – CS 71058  
13781 Aubagne Cedex

Phone +33.442.845600  
Fax +33.442.846545

### Austria

Sartorius Stedim Austria GmbH  
Franzosengraben 12  
1030 Vienna

Phone +43.1.7965763.18  
Fax +43.1.796576344

### Belgium

Sartorius Stedim Belgium N.V.  
Leuvensesteenweg, 248/B  
1800 Vilvoorde

Phone +32.2.756.06.80  
Fax +32.2.756.06.81

### Hungary

Sartorius Stedim Hungária Kft.  
Kagyló u. 5  
2092 Budakeszi

Phone +36.23.457.227  
Fax +36.23.457.147

### Italy

Sartorius Stedim Italy S.p.A.  
Via dell'Antella, 76/A  
50012 Antella-Bagno a Ripoli (FI)

Phone +39.055.63.40.41  
Fax +39.055.63.40.526

### Poland

Sartorius Stedim Poland Sp. z o.o.  
ul. Wrzesinska 70  
62-025 Kostrzyn

Phone +48.61.647.38.40  
Fax +48.61.879.25.04

### Russian Federation

LLC "Sartorius ICR"  
Uralskaya str. 4, Lit. B  
199155, Saint-Petersburg

Phone +7.812.327.5.327  
Fax +7.812.327.5.323

### Spain

Sartorius Stedim Spain SA  
C/Isabel Colbrand 10,  
Oficina 70  
Polígono Industrial de Fuencarral  
28050 Madrid

Phone +34.90.2110935  
Fax +34.91.3589623

### Switzerland

Sartorius Stedim Switzerland AG  
Ringstrasse 24 a  
8317 Tagelswangen

Phone +41.52.354.36.36  
Fax +41.52.354.36.46

### U.K.

Sartorius Stedim UK Ltd.  
Longmead Business Centre  
Blenheim Road, Epsom  
Surrey KT19 9 QQ

Phone +44.1372.737159  
Fax +44.1372.726171

## America

### USA

Sartorius Stedim North America Inc.  
5 Orville Drive, Suite 200  
Bohemia, NY 11716

Toll-Free +1.800.368.7178  
Fax +1.631.254.4253

### Argentina

Sartorius Argentina S.A.  
Int. A. Ávalos 4251  
B1605ECS Munro  
Buenos Aires

Phone +54.11.4721.0505  
Fax +54.11.4762.2333

### Brazil

Sartorius do Brasil Ltda  
Av. Dom Pedro I, 241  
Bairro Vila Pires  
Santo André  
São Paulo  
Cep 09110-001

Phone +55.11.4451.6226  
Fax +55.11.4451.4369

### Mexico

Sartorius de México S.A. de C.V.  
Circuito Circunvalación Poniente  
No. 149  
Ciudad Satélite  
53100, Estado de México  
México

Phone +52.5555.62.1102  
Fax +52.5555.62.2942

## Asia | Pacific

### Australia

Sartorius Stedim Australia Pty. Ltd.  
Unit 5, 7-11 Rodeo Drive  
Dandenong South Vic 3175

Phone +61.3.8762.1800  
Fax +61.3.8762.1828

### China

Sartorius Stedim Biotech (Beijing) Co. Ltd.  
No. 33 Yu'an Road  
Airport Industrial Park Zone B  
Shunyi District, Beijing 101300

Phone +86.10.80426516  
Fax +86.10.80426580

Sartorius Stedim Biotech (Beijing) Co. Ltd.  
Shanghai Branch Office  
3rd Floor, North Wing, Tower 1  
No. 4560 Jin Ke Road  
Pudong District, Shanghai 201210

Phone +86.21.68782300  
Fax +86.21.68782332 | 68782882

Sartorius Stedim Biotech (Beijing) Co. Ltd.  
Guangzhou Representative Office  
Unit K, Building 23  
Huihua Commerce & Trade Building  
No. 80 Xianlie Middle Road  
Guangzhou 510070

Phone +86.20.37618687 | 37618651  
Fax +86.20.37619051

### India

Sartorius Stedim India Pvt. Ltd.  
#69/2-69/3, NH 48, Jakkasandra  
Nelamangala Tq  
562 123 Bangalore, India

Phone +91.80.4350.5250  
Fax +91.80.4350.5253

### Japan

Sartorius Stedim Japan K.K.  
4th Fl., Daiwa Shinagawa North Bldg.  
8-11, Kita-Shinagawa 1-chome  
Shinagawa-ku, Tokyo, 140-0001 Japan

Phone +81.3.4331.4300  
Fax +81.3.4331.4301

### Malaysia

Sartorius Stedim Malaysia Sdn. Bhd.  
Lot L3-E-3B, Enterprise 4  
Technology Park Malaysia  
Bukit Jalil  
57000 Kuala Lumpur, Malaysia

Phone +60.3.8996.0622  
Fax +60.3.8996.0755

### Singapore

Sartorius Stedim Singapore Pte. Ltd.  
1 Science Park Road,  
The Capricorn, #05-08A,  
Singapore Science Park II  
Singapore 117528

Phone +65.6872.3966  
Fax +65.6778.2494

### South Korea

Sartorius Korea Biotech Co., Ltd.  
8th Floor, Solid Space B/D,  
PanGyoYeok-Ro 220, Bundang-Gu  
SeongNam-Si, GyeongGi-Do, 463-400

Phone +82.31.622.5700  
Fax +82.31.622.5799



◀ [www.sartorius.com](http://www.sartorius.com)